


Scherzo SS-C18
Scherzo SM-C18
Scherzo SW-C18

100 x 4.6 mm

Application

Oral hypoglycemic drug (metformin)

ピグアナイド系経口糖尿病治療薬 (メトホルミン)

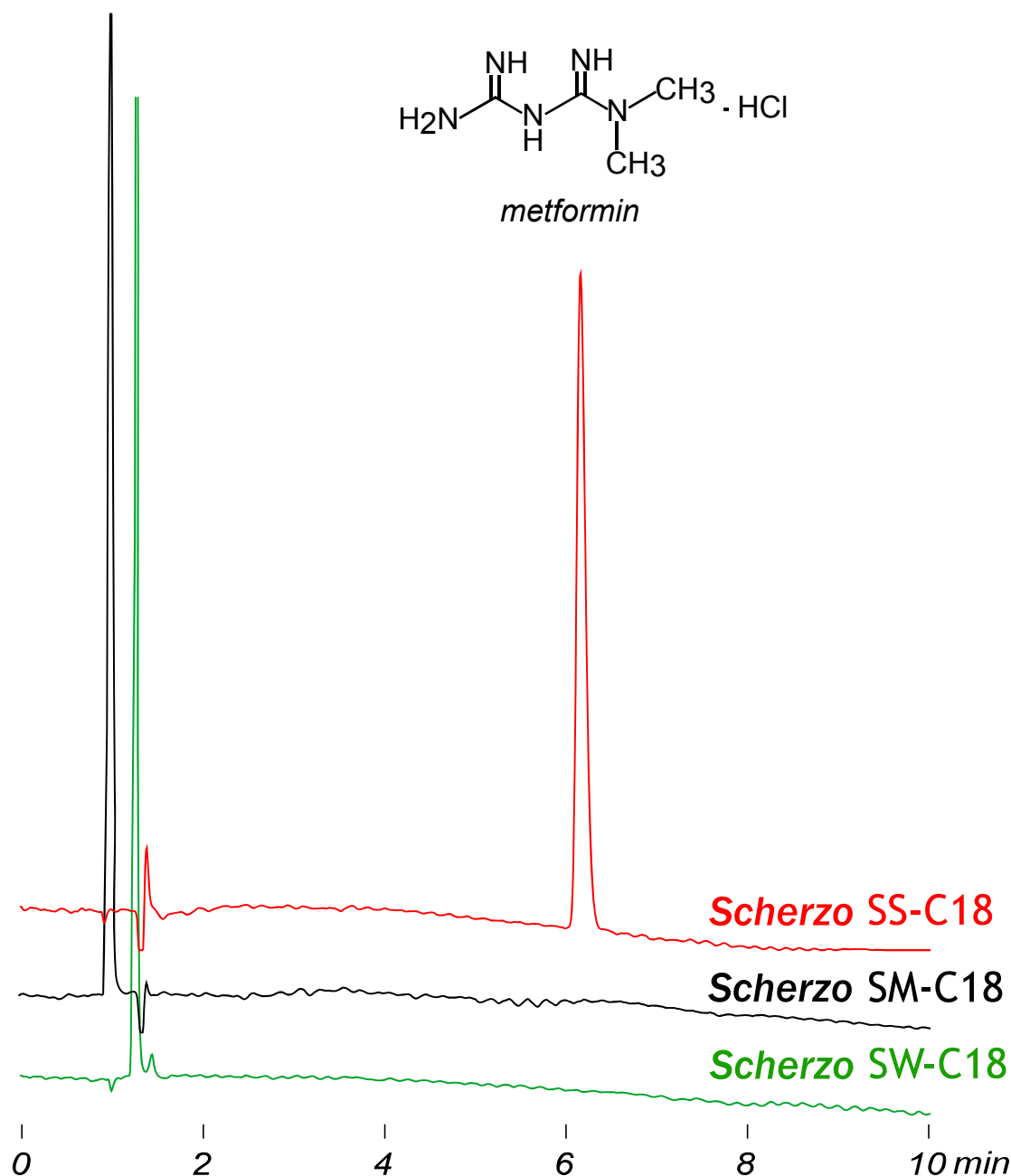

150 x 4.6 mm

A: acetonitrile / water / formic acid = 5 / 95 / 0.1

B: acetonitrile / 100mM HCOONH₄ = 50 / 50

40%B (0-1min), 40-60%B (1-10min)

1 mL/min (8MPa), 37°C

5 μL (0.1mg/mL, water), 230nm