





Intrada Amino Acid

100 x 3 mm

Application

LC-MS: HbA1c and AGEs related amadori compounds

アマドリ化合物(ヘモグロビンA1cおよびAGEs関連)のLC-MS分析

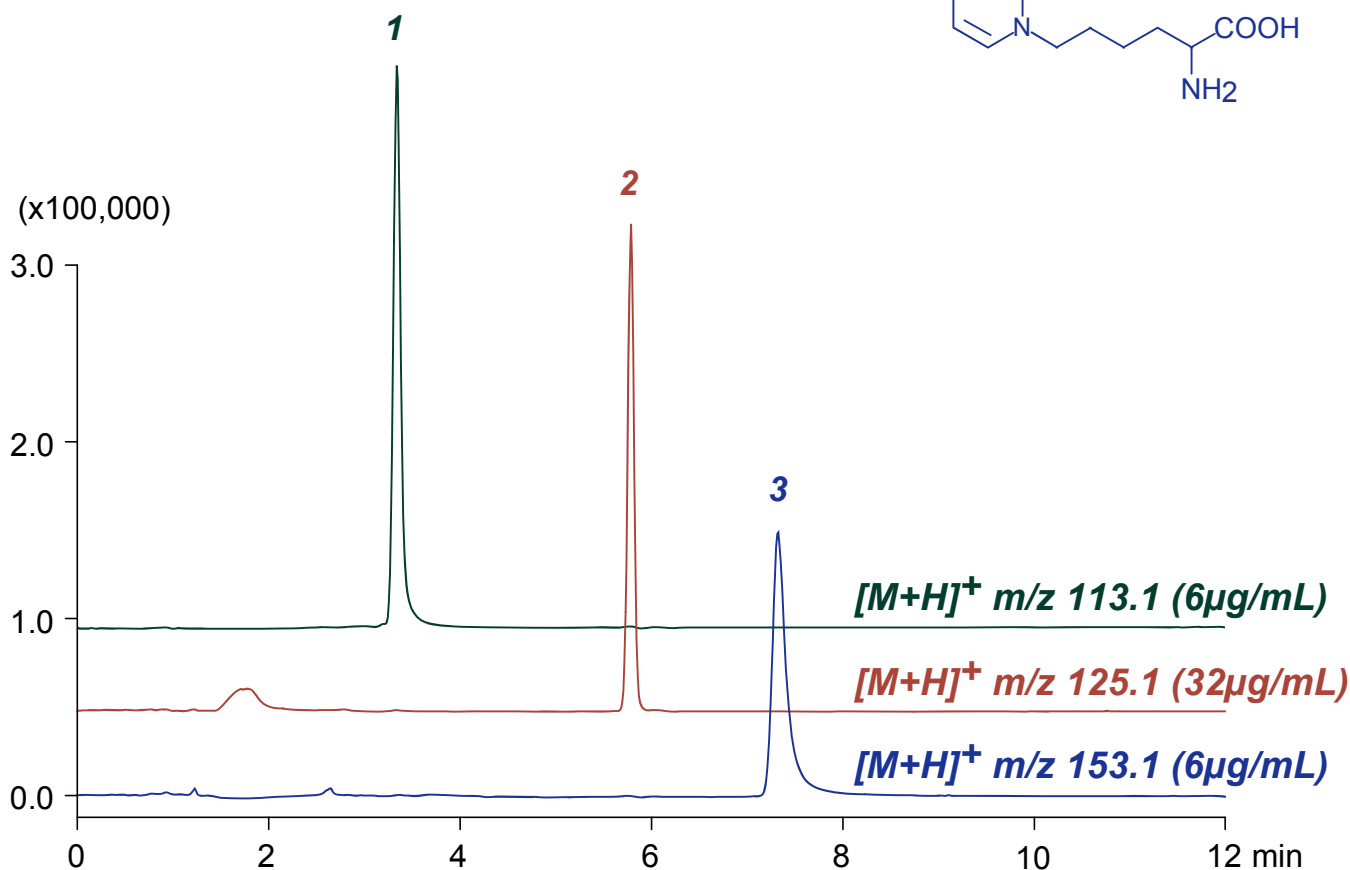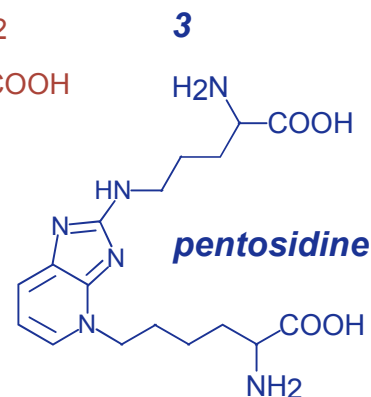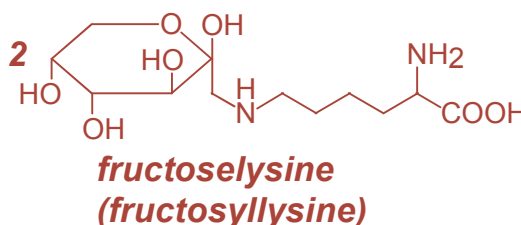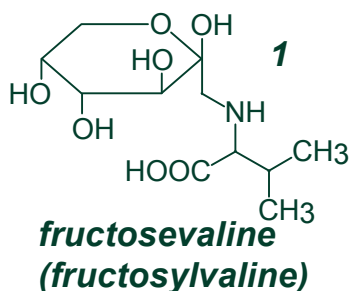

Intrada Amino Acid, 100 x 3 mm

A: acetonitrile / formic acid = 100 / 0.1

B: acetonitrile / 100mM ammonium formate = 20 / 80

10%B (0-1min), 10-100%B (1-4min), 100%B (4-12min)

0.6mL/min (4MPa), 37°C, MS (ESI, SIM)

5μL (6-32μg/mL, 1% HCOOH)