



Cadenza **CX-C18**
 Unison **UK-C8**
 Cadenza **CD-C18**
 Cadenza **CW-C18**
 Unison **UK-C18**

75 x 4.6 mm

Application

Enrofloxacin (fluoroquinolone)

エンロフロキサシン(キノロン系抗菌薬)

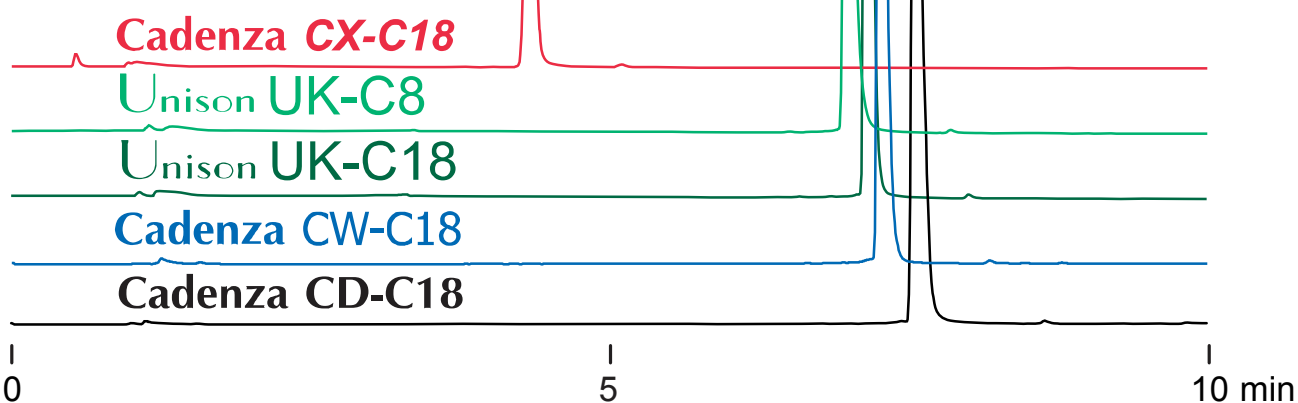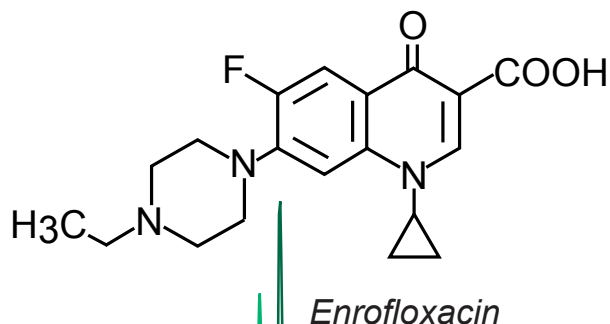

75 x 4.6 mm
 A: water / formic acid = 100 / 0.1
 B: acetonitrile / formic acid = 100 / 0.1
 0-30%B (0-10min)
 1.0 mL/min (7MPa), 37°C, 290nm
 10 µL (0.1mg/mL, 1%HCCOH)