



Cadenza CX-C18  
 Unison UK-C8  
 Cadenza CW-C18  
 Unison UK-C18  
 Cadenza CD-C18

75 x 4.6 mm

Application

**Norfloxacin (fluoroquinolone)**

ノルフロキサシン (キノロン系抗菌薬)

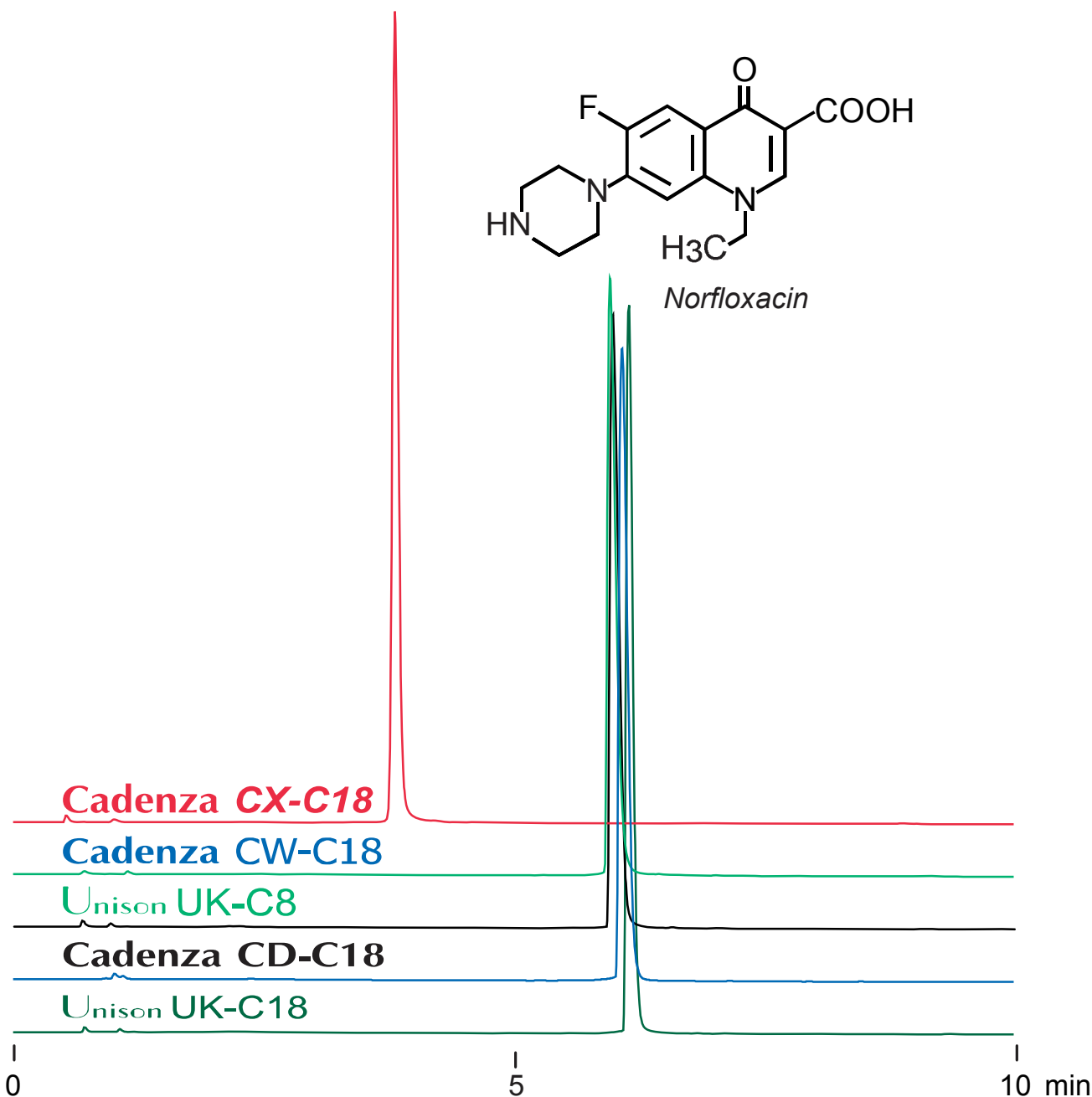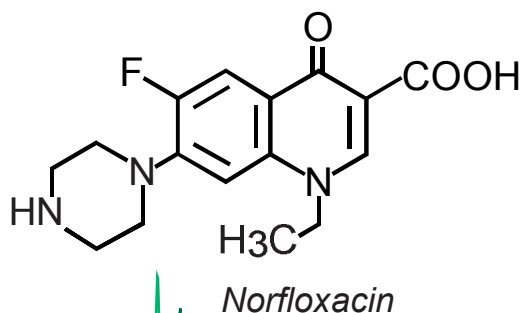

75 x 4.6mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-30%B (0-10min)

1.0 mL/min (7MPa), 37°C , 290nm

10 µL (0.1mg/mL, 1%HCCOH)