


Cadenza **CX-C18**
 Unison **UK-C8**
 Cadenza **CD-C18**
 Cadenza **CW-C18**
 Unison **UK-C18**

75 x 4.6 mm


Application

Ciprofloxacin (fluoroquinolone)

シプロフロキサシン(キノロン系抗菌薬)

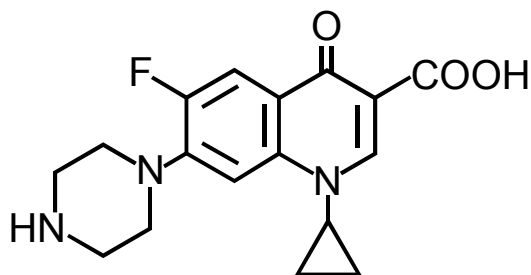

Ciprofloxacin

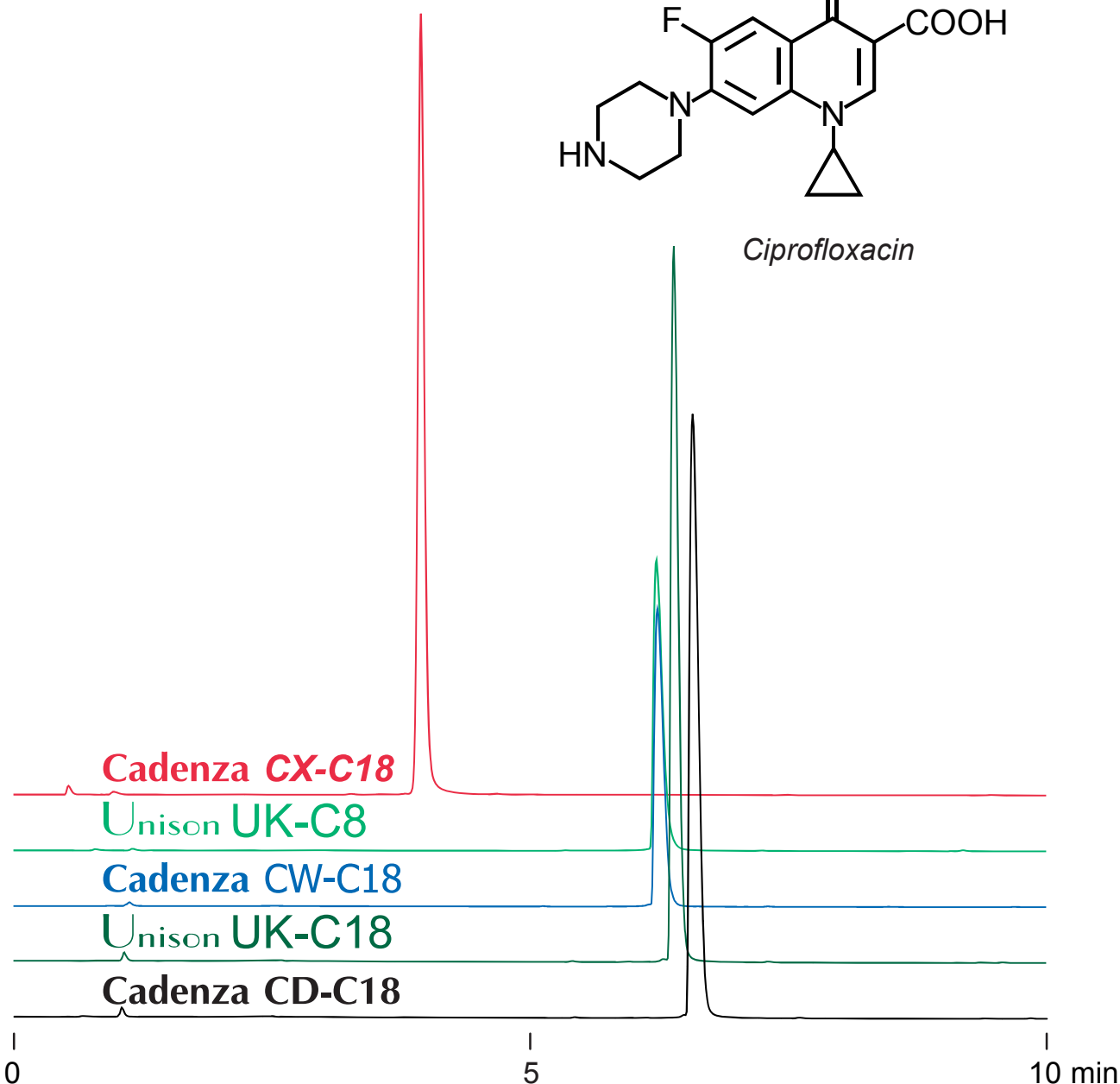

75 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-30%B (0-10min)

1.0 mL/min (7MPa), 37°C, 290nm

10 µL (0.1mg/mL, 1%HCCOH)