



Nardis ND-NX
Unison UK-Amino

250 x 4.6 mm

HPLC COLUMN APPLICATION

Bronopol (Antiseptic)

ブロナポール (殺菌消毒剤)

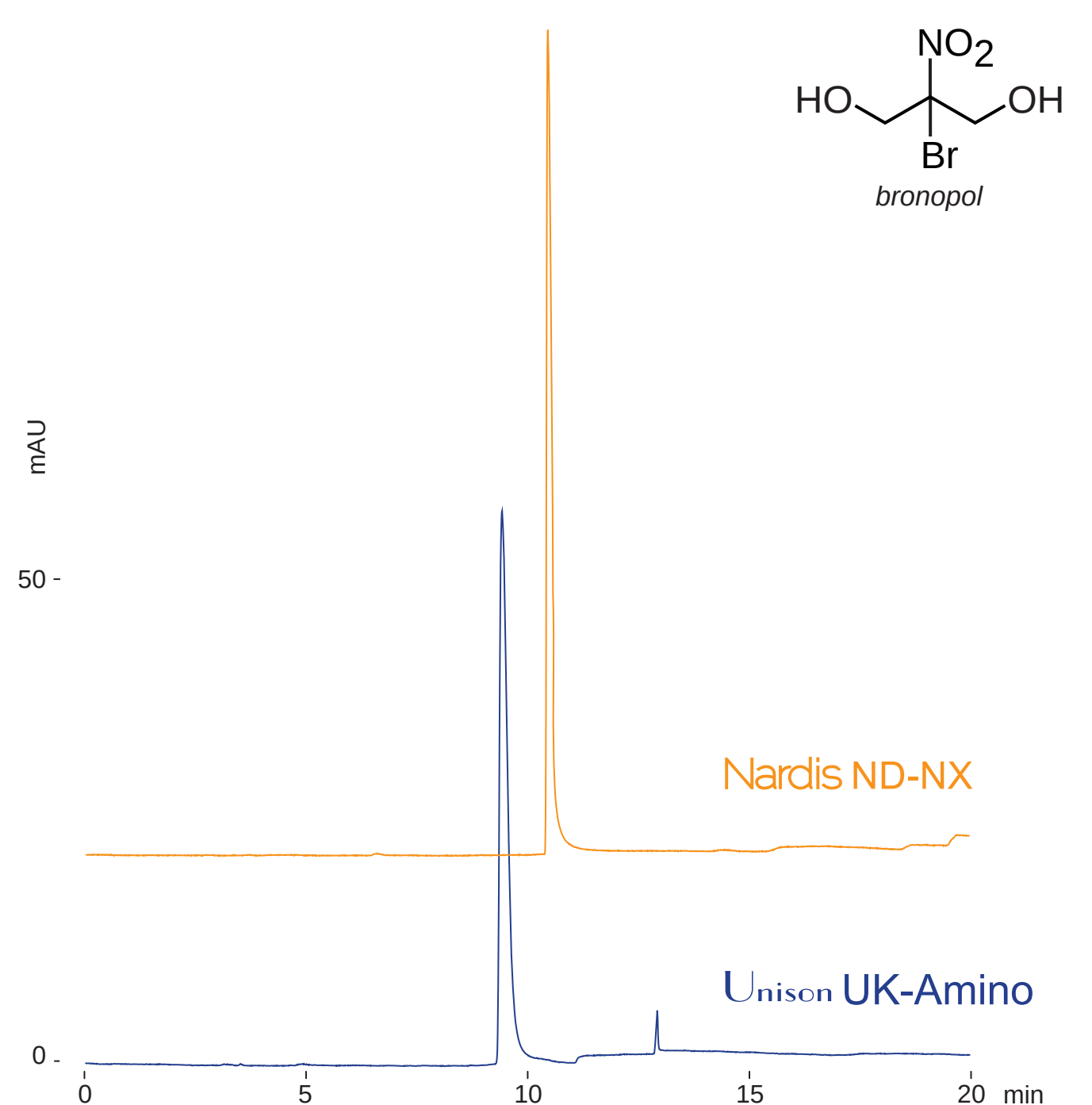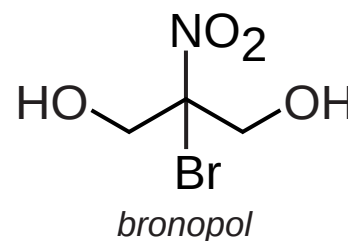

250 x 4.6 mm

A: acetonitrile

B: acetonitrile / 100mM ammonium acetate = 50 / 50

0%B (0-1min), 0-10%B (1-20min), 0%B (20-30min)

1.0mL/min (6-7MPa), 37°C, 250nm

10μL (2.0mg/mL, acetonitrile)