


Nardis ND-NX
Unison UK-Amino

150 x 4.6 mm

HPLC COLUMN APPLICATION

Cidofovir (antiviral drug) シドフォビル (抗ウイルス薬)

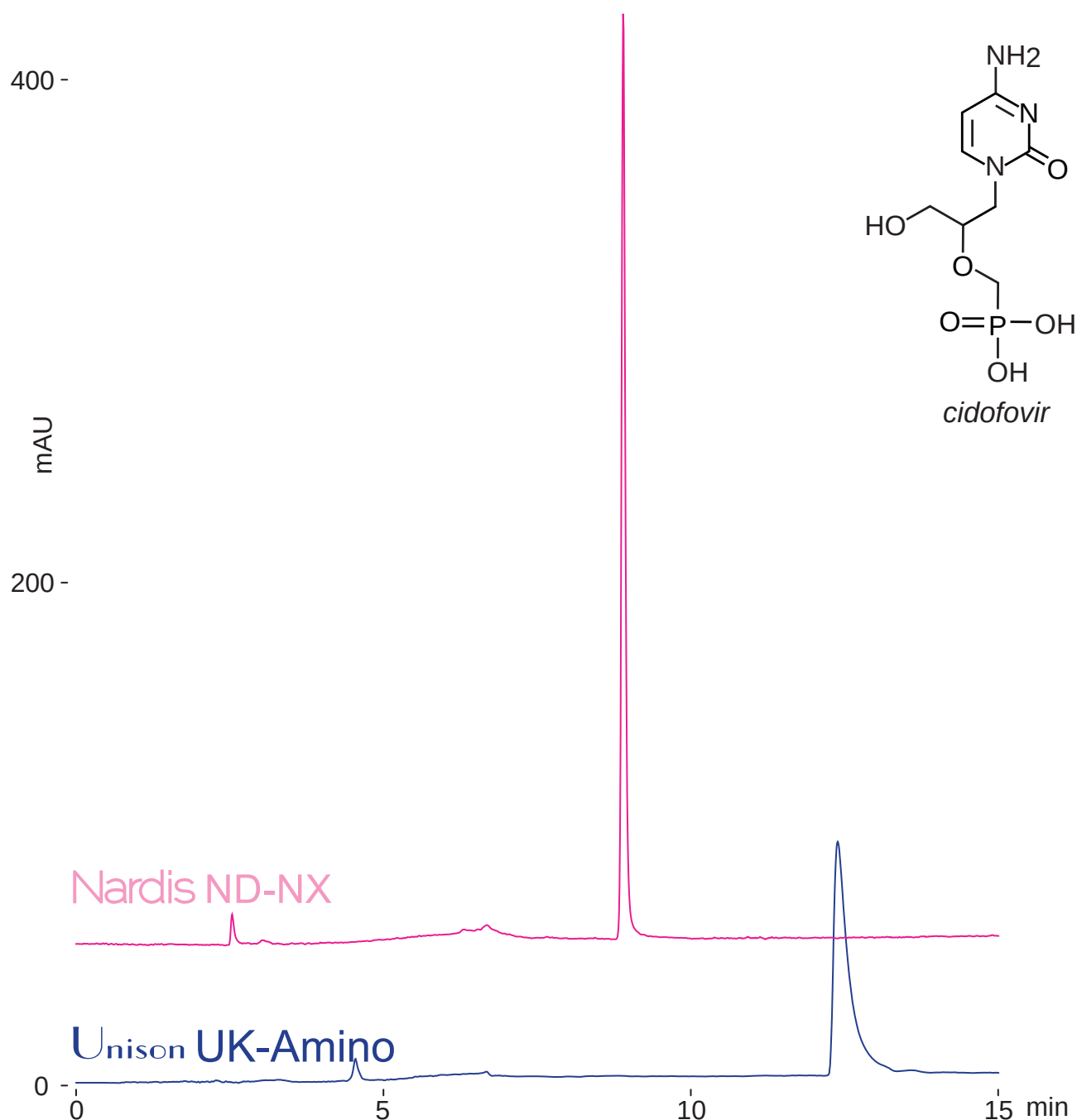

150 x 4.6 mm

A: 100mM ammonium formate / acetonitrile = 5 / 95

B: 100mM ammonium formate / acetonitrile = 90 / 10

0-100%B (0-15min)

1.0mL/min , 37°C, 260nm

10μL (0.3mg/mL, 100mM ammonium formate / acetonitrile = 50/ 50)