

**Unison UK-C18**



**Cadenza CD-C18**

100 x 4.6 mm

**HPLC COLUMN APPLICATION**

**Reactive Black 5 (Reactive dye)**

リアクティブブラック5(反応性染料)

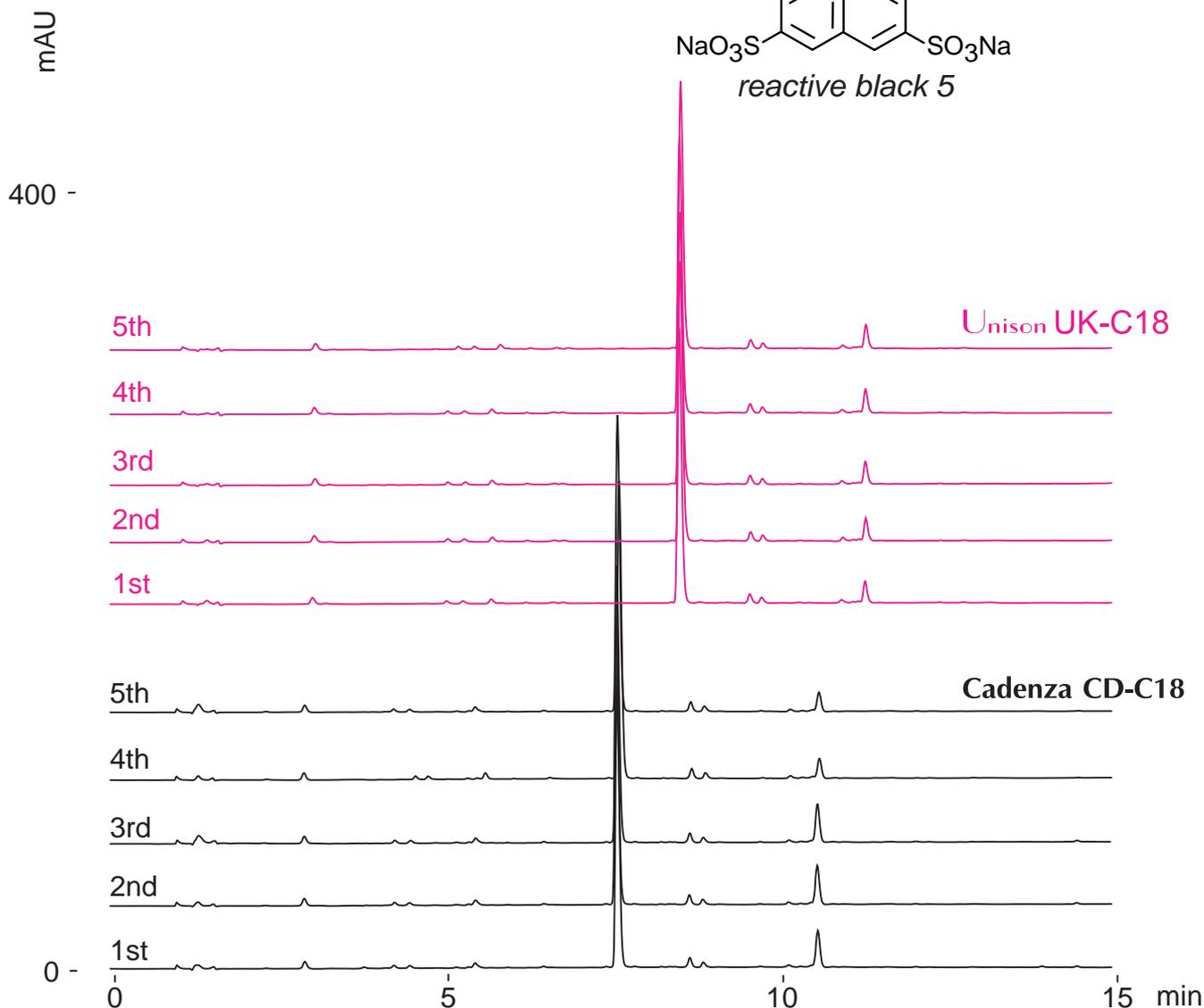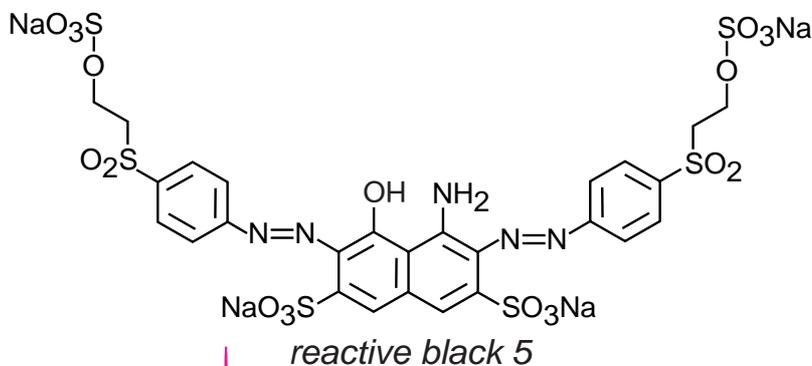

100 x 4.6 mm

A: 10mM ammonium formate / acetonitrile = 95 / 5

B: 100mM ammonium formate / acetonitrile = 50 / 50

0%B (0-1min), 0-75%B (1-15min)

1.0mL/min (13MPa), 37°C, 260nm

10µL (0.25mg/mL, 100mM ammonium formate)