


Scherzo SW-C18

100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Trialkylamine

トリアルキルアミン類の LC-MS 分析

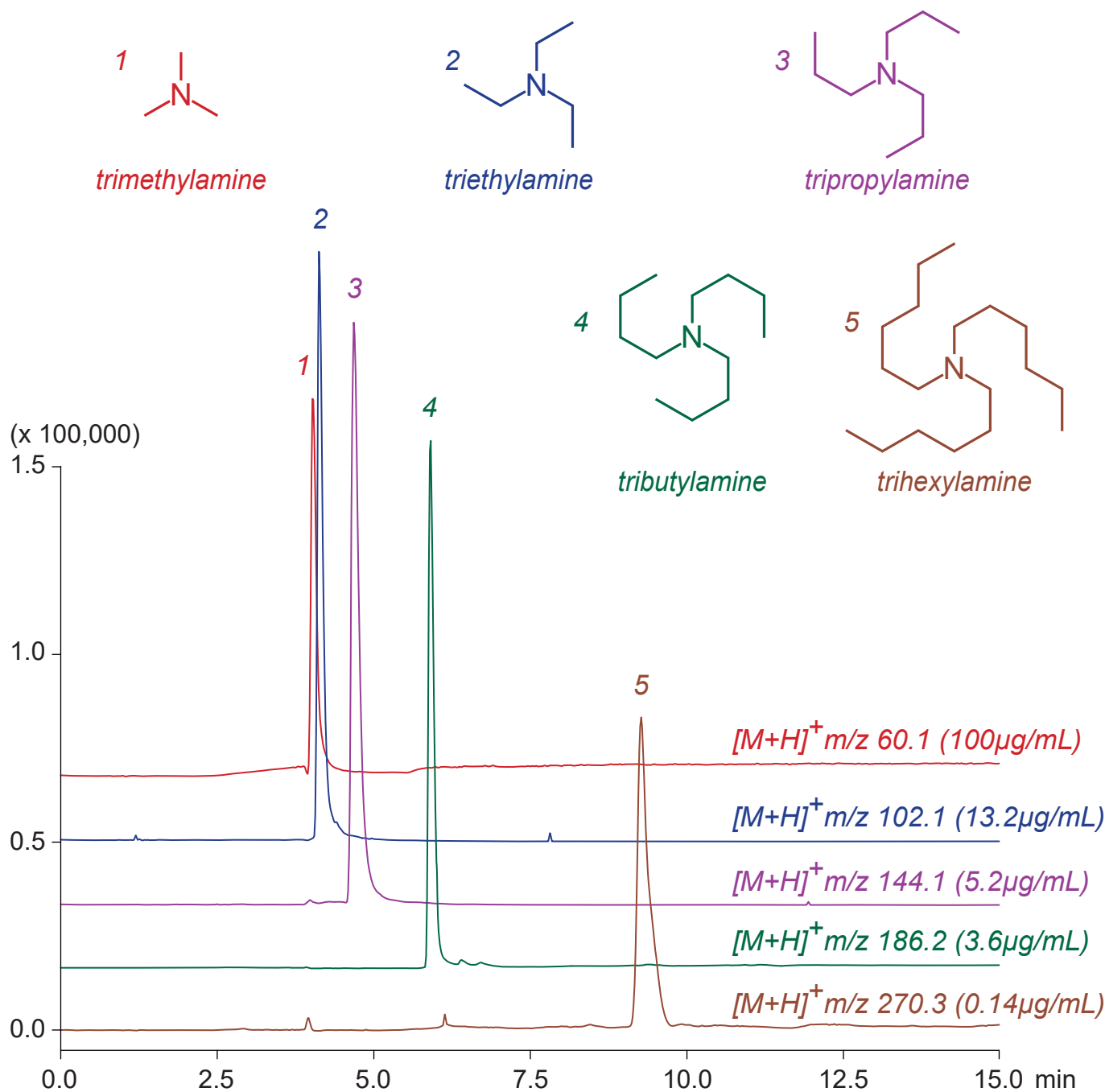

Scherzo SW-C18, 100 x 3 mm

A: methanol / water / acetic acid = 50 / 50 / 0.1

B: methanol / 100mM ammonium acetate = 80 / 20

0-20%B (0-4min), 100%B (4-15min)

0.4mL/min (13MPa), 37°C, Single Quad. MS (ESI)

5µL (0.14-100µg/mL, EtOH / AcOH = 100 / 1)