


Cadenza CD-C18


Application

30 x 4.6 mm 30 x 3 mm 30 x 2 mm

Agricultural Chemicals (4)

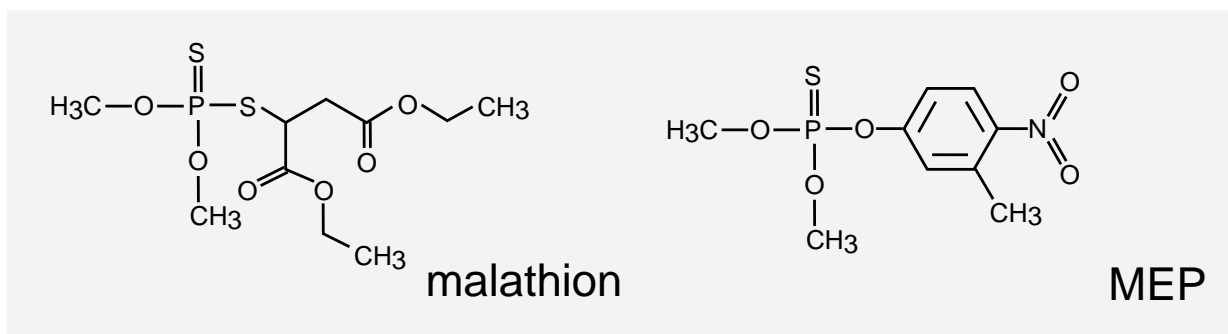


30 x 4.6 mm

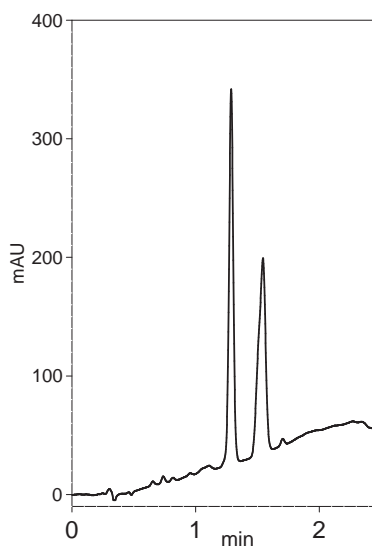


1.0 mL/min, 2.0 MPa
2.0 uL

30 x 3 mm

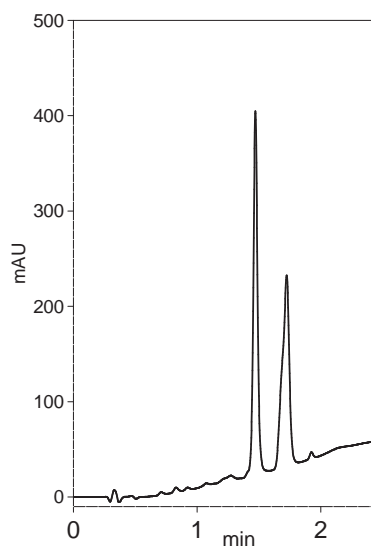

0.43 mL/min, 2.3 MPa
1.0 uL

30 x 2 mm

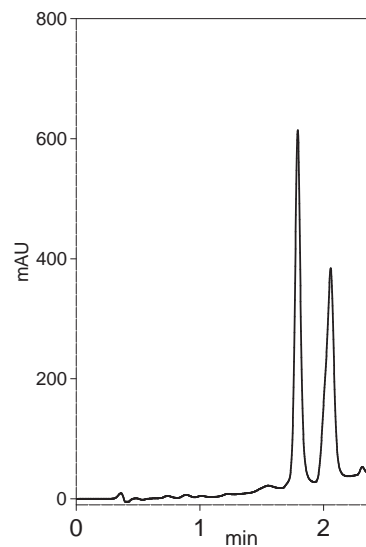

0.2 mL/min, 2.2 MPa
1.0 uL

Cadenza CD-C18
A: water / TFA = 100 / 0.1
B: acetonitrile / TFA = 100 / 0.1
50 - 100 %B (0-2min)
UV at 220 nm, 37 °C