


Cadenza CD-C18

150 x 3 mm


Application

Turmelic (Ukon)

Curcuma aromatica Salisb. (Haru-Ukon)

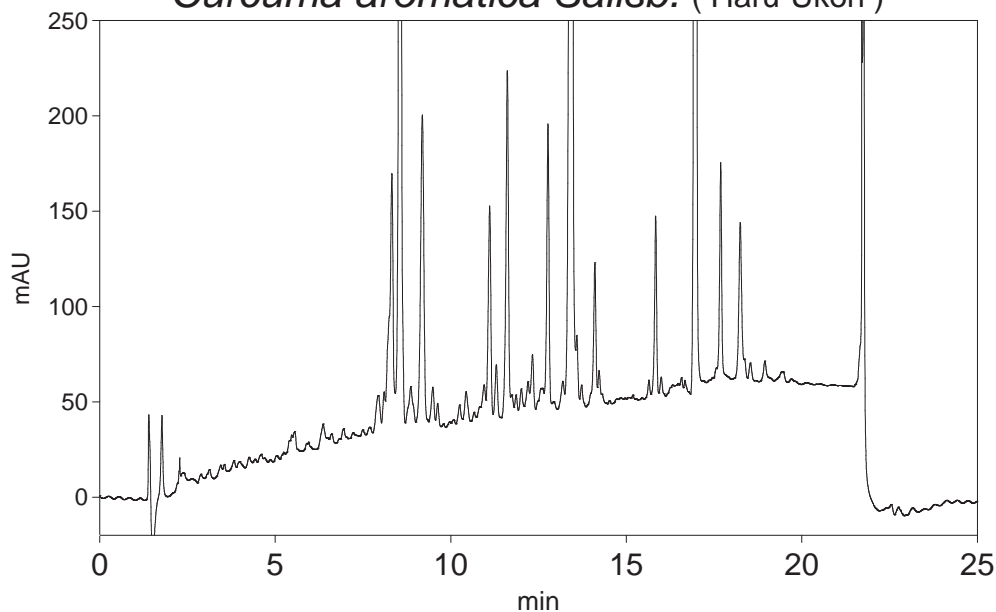

Curcuma zedoaria Roscoe. (Murasaki-Ukon)

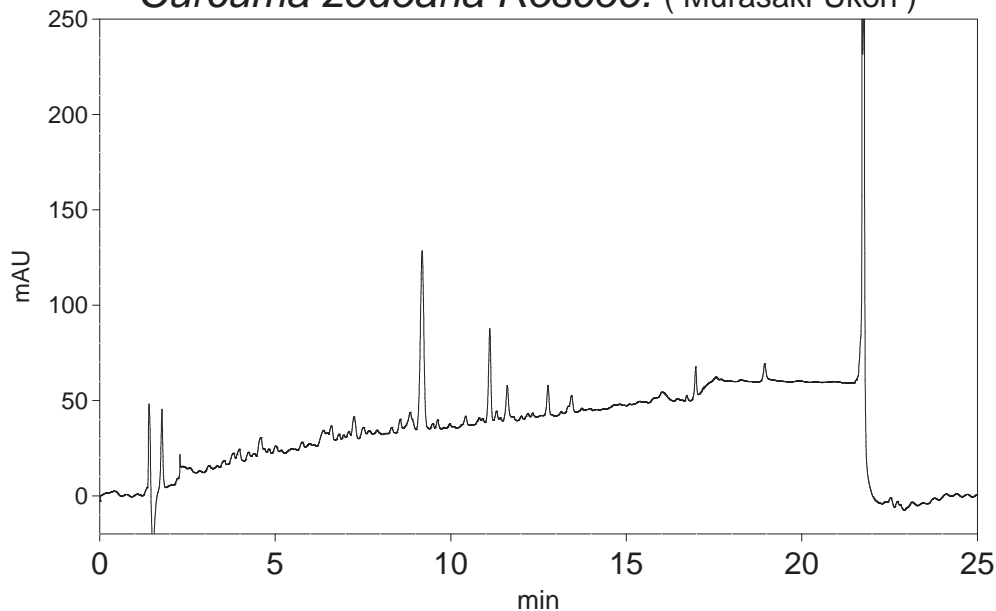

Cadenza CD-C18, 150 x 3.0 mm
 A : water / formic acid = 100 / 0.05
 B : acetonitrile / formic acid = 100 / 0.05
 30 - 100%B (0 - 15 min) , 100%B (15 - 20 min)
 0.5 mL/min, 37 °C, 12.9 MPa, UV at 220 nm