


Unison US-C18

150 x 4.6 mm

Application

Organic Acids

有機酸

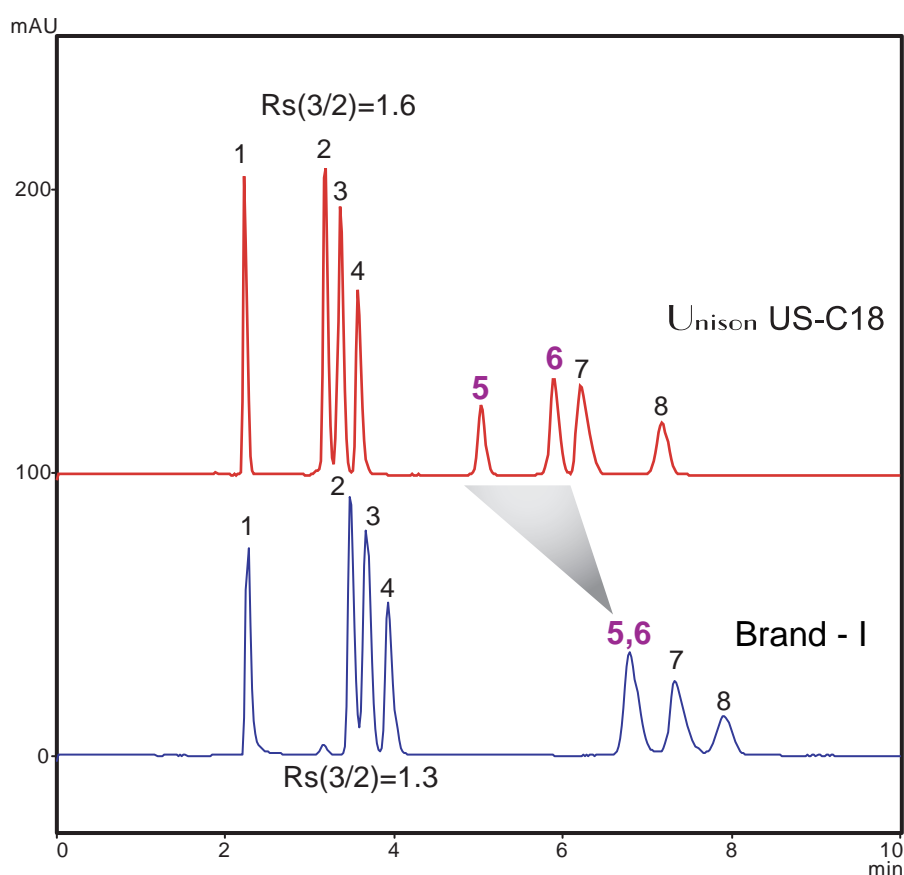

- 1 HCOOH
formic acid
- 2 HOOCCH2COOH
malonic acid
- 3 CC(O)C(=O)O
lactic acid
- 4 CC(=O)O
acetic acid
- 5 OC(=O)C=CC(=O)O
maleic acid
- 6 OC(=O)C(C(=O)O)CC(=O)O
citric acid
- 7 OC(=O)CC(=O)O
succinic acid
- 8 OC(=O)C=CC(=O)O
fumaric acid

150 x 4.6 mm
20 mM H₃PO₄
1.0 mL/min, 37 °C, UV at 210 nm