



Unison UK-Silica
Cadenza CD-C18

150 x 3 mm


Application

Tocopherol Nicotinate Capsule

ニコチン酸トコフェロールカプセル

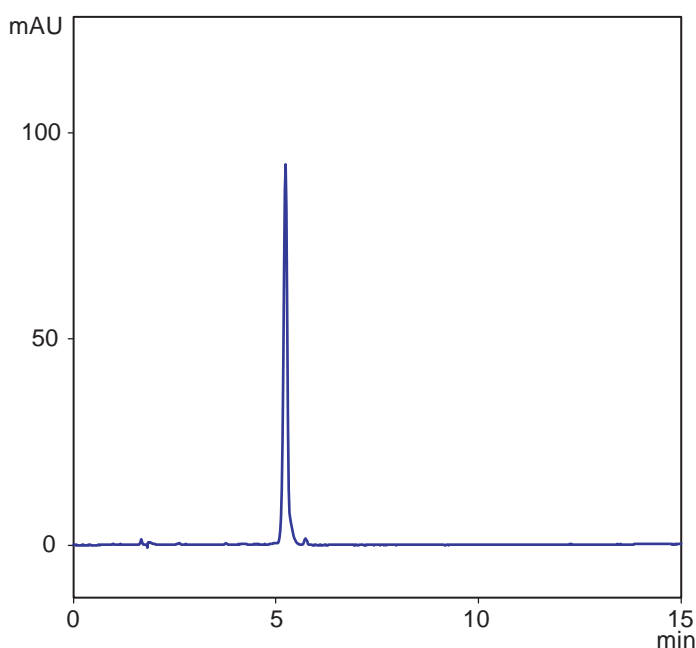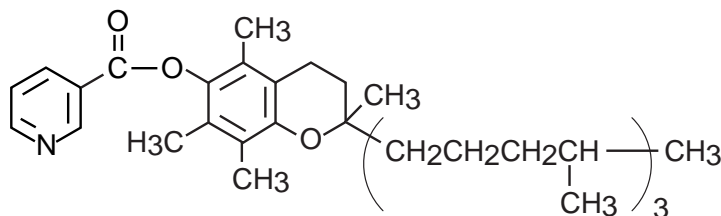

Unison UK-Silica, 150 x 3 mm
hexane / THF = 95 / 5
0.5 mL/min (4MPa), 37 °C, 295 nm
tocopherol nicotinate 5 μ g

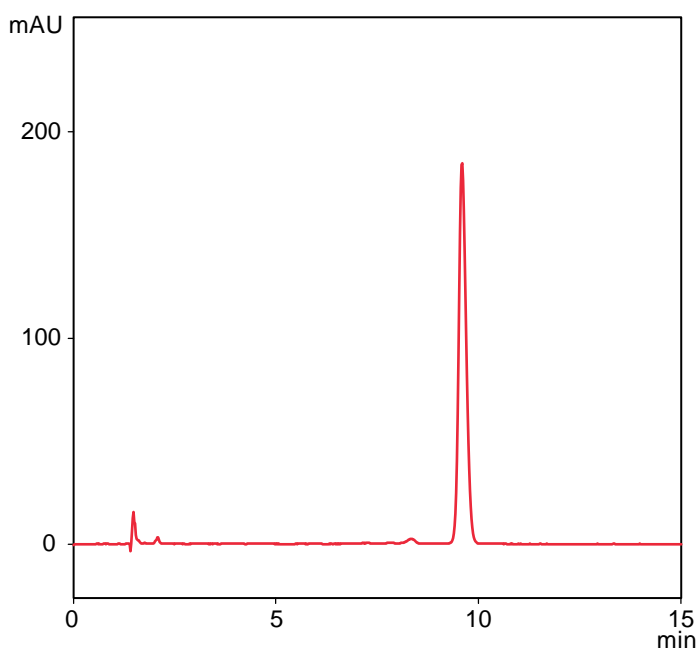

Cadenza CD-C18, 150 x 3 mm
methanol, 0.5 mL/min (6MPa)
37 °C, 295 nm
tocopherol nicotinate 5 μ g