




Unison UK-Silica  
 Unison UK-Phenyl  
 Unison UK-C8  
 Unison UK-C18  
 Cadenza CL-C18

150 x 4.6 mm  
 75 x 4.6 mm


Application

## Calcium Pantothenate and Ascorbic Acid

パントテン酸カルシウムとアスコルビン酸

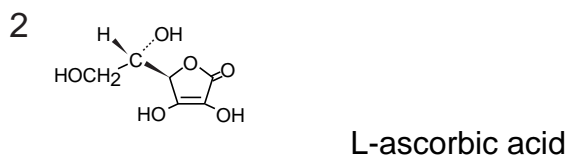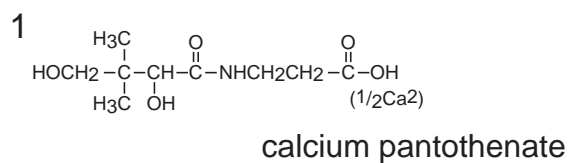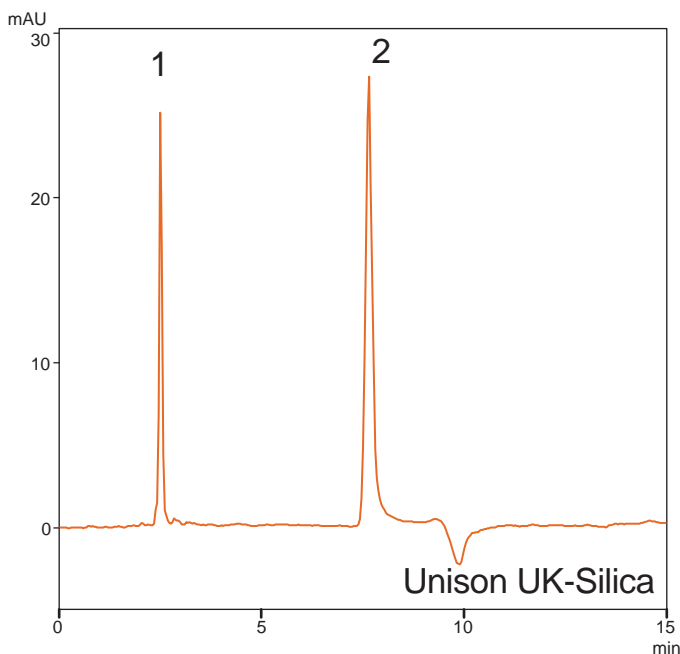

Unison UK-Silica, 150 x 4.6 mm  
 50mM ammonium acetate / acetonitrile  
 = 10 / 90  
 1 mL/min, 37 °C, 210 nm, 2µL(1-2µg)

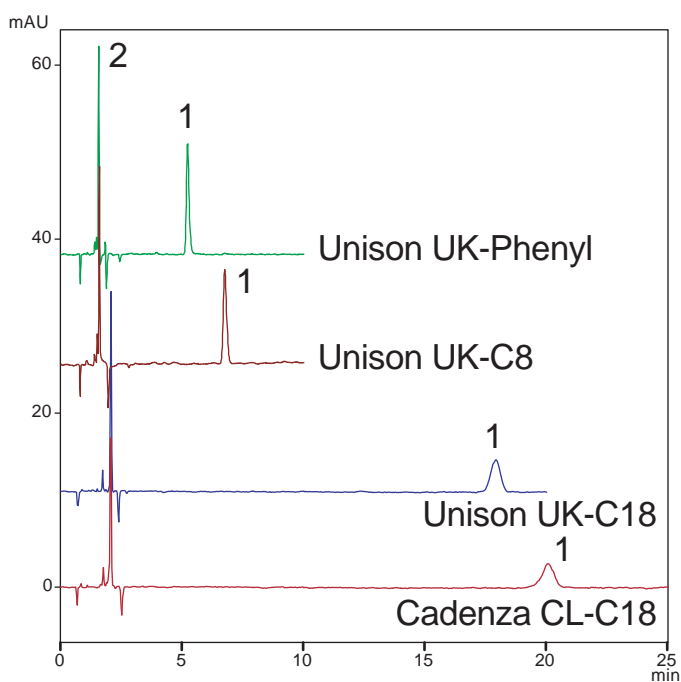

75 x 4.6 mm  
 20mM NaH<sub>2</sub>PO<sub>4</sub> +  
 10mM dibutylammonium acetate  
 1 mL/min, 37 °C, 210 nm, 2µL(0.5-1µg)