

Unison UK-C8



50 x 2 mm

Application


LC-MS Application for Calcium Antagonist

LC-MSアプリケーション (カルシウム拮抗薬, 塩酸ジルチアゼム)

API4000 (TAKARA BIO INC.)

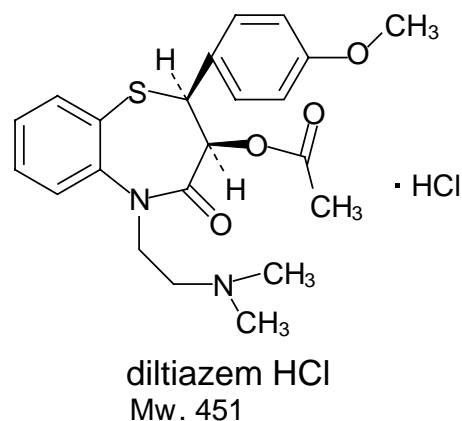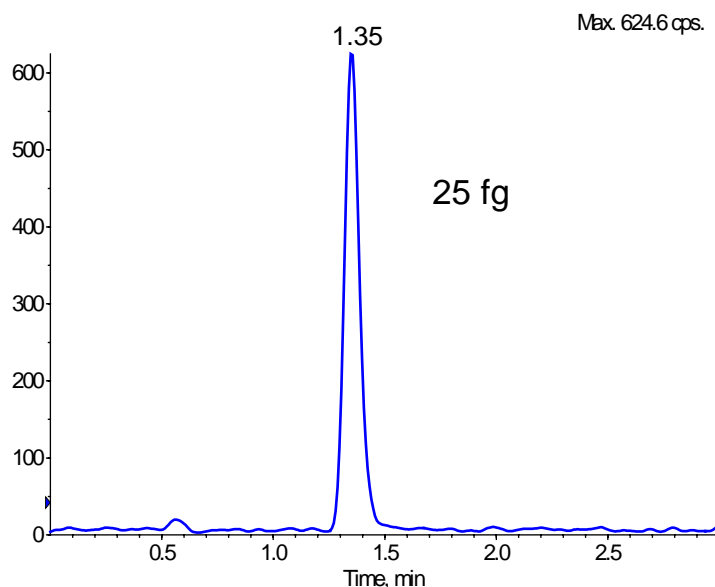

Results.rdb (Diltiazem): "Linear" Regression ("1 / x" weighting): $y = 5.98e+005 x + -22.1$ ($r = 0.9999$)

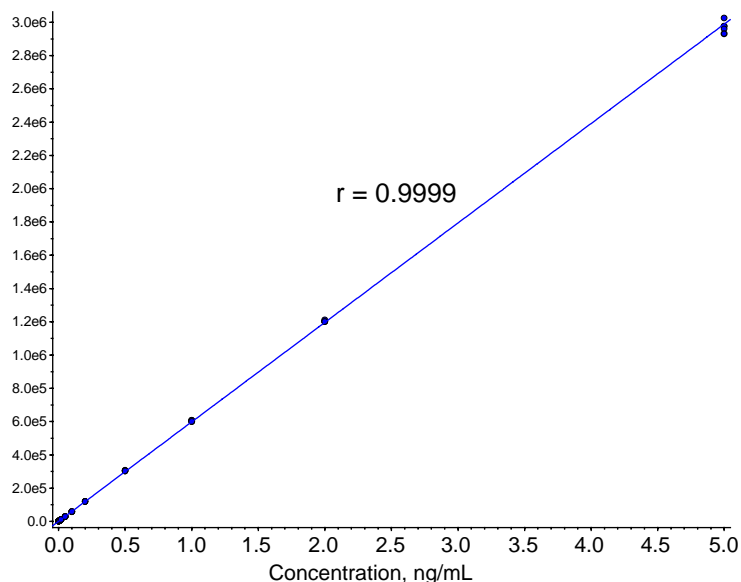

Unison UK-C8, 50 x 2 mm

0.1% formic acid / acetonitrile = 2 / 1, 0.25 mL/min, 5 uL(5 fg/uL)

API4000: ESI, MRM Positive, Q1/Q3 : 415 / 178

Courtesy of J.Watanabe, TAKARA BIO INC.