






Unison UK-Phenyl
 Unison UK-C8
 Cadenza CL-C18
 Cadenza CD-C18
 Unison UK-C18

75 x 4.6 mm

Application

Veterinary Medicine (Synthetic Antimicrobial)

動物用医薬品(合成抗菌剤)

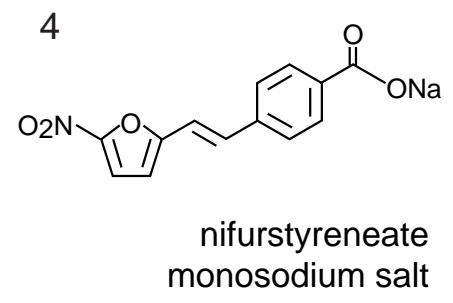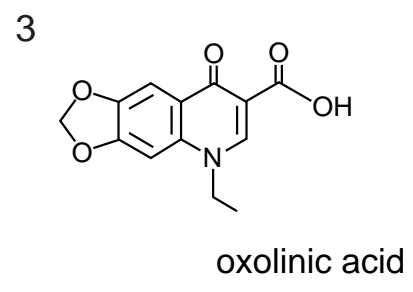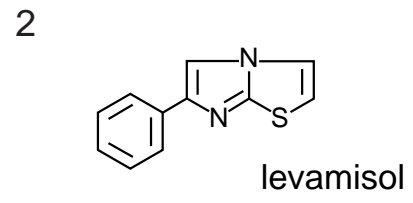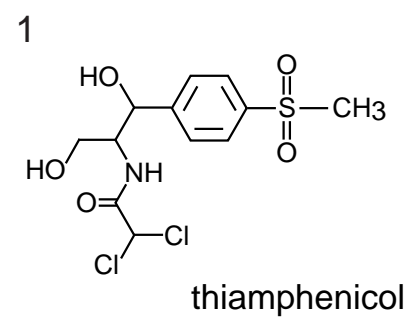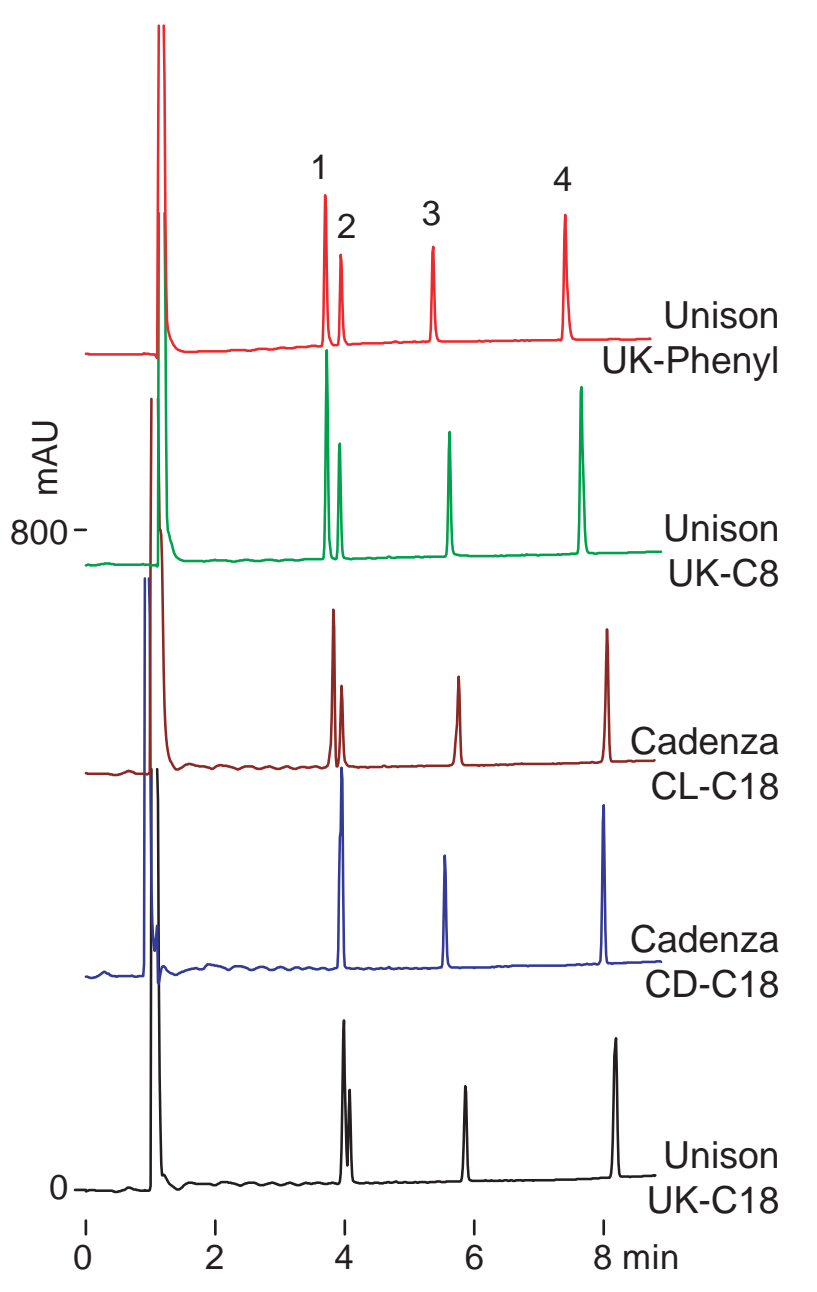

75 x 4.6 mm,
 A: water / trifluoroacetic acid = 100 / 0.1
 B: acetonitrile / trifluoroacetic acid = 100 / 0.1
 5-60 %B (0-8 min), 1 mL/min (7-8 MPa),
 37 deg.C, 220 nm, 2 uL (0.2-0.6 ug)