



Cadenza CD-C18

75 x 4.6 mm




Application

Natural Medicine Properties (Shoseiryuto)

生薬成分 (小青竜湯)

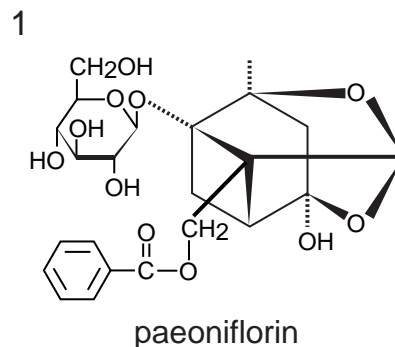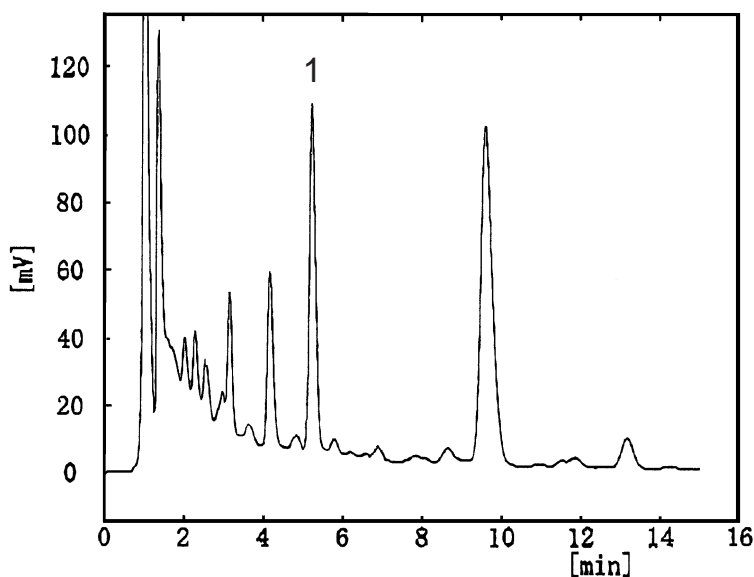

Cadenza CD-C18,  
75 x 4.6 mm  
water / acetonitrile / H<sub>3</sub>PO<sub>4</sub>  
= 700 / 100 / 1  
1 mL/min, 25 deg.C, 230 nm

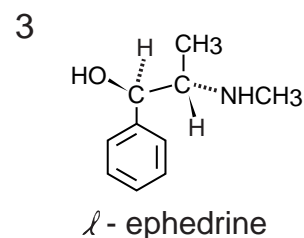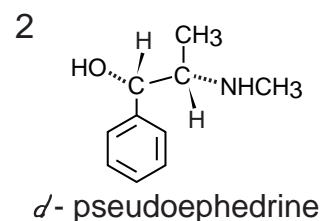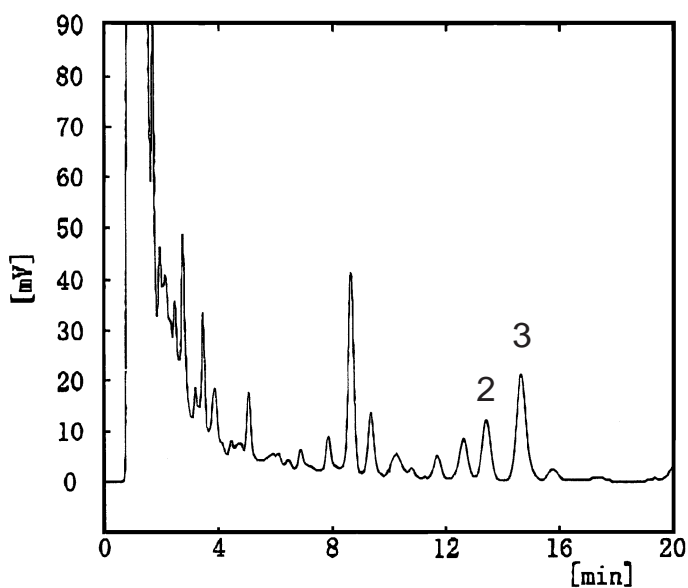

Cadenza CD-C18, 75 x 4.6 mm  
water / acetonitrile / H<sub>3</sub>PO<sub>4</sub> / SDS  
= 650 / 350 / 2 / 5  
1 mL/min, 40 deg.C, 210 nm