


**Unison UK-Amino**  
**Unison UK-C18**

250 x 4.6 mm



Application

### Steroids

ステロイド

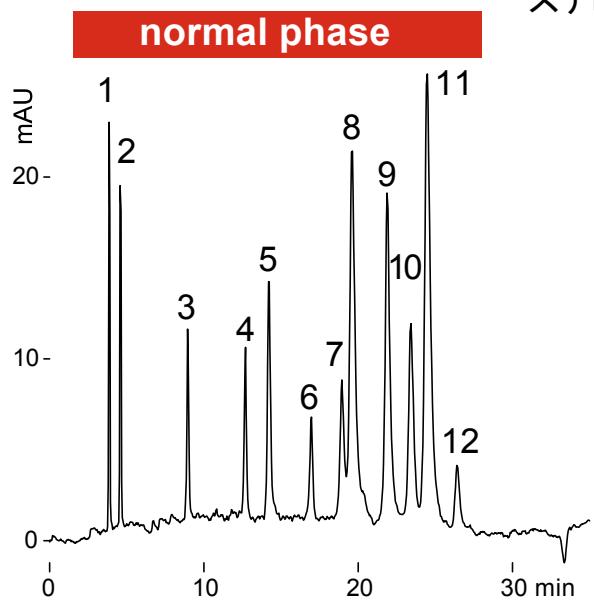

**Unison UK-Amino, 250 x 4.6 mm**  
 A: hexane / acetic acid = 100 / 0.1  
 B: ethanol / acetic acid = 100 / 1  
 10-30% B (0-30 min), 1 mL/min  
 37 deg.C, 260 nm, 2 uL (1-2ug)

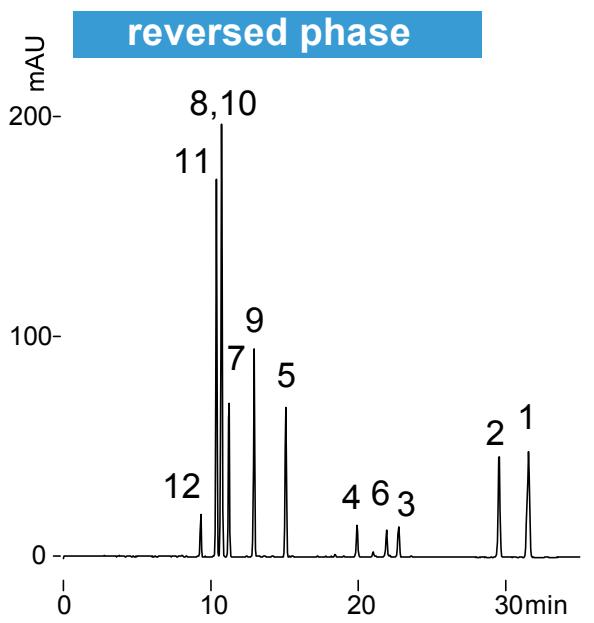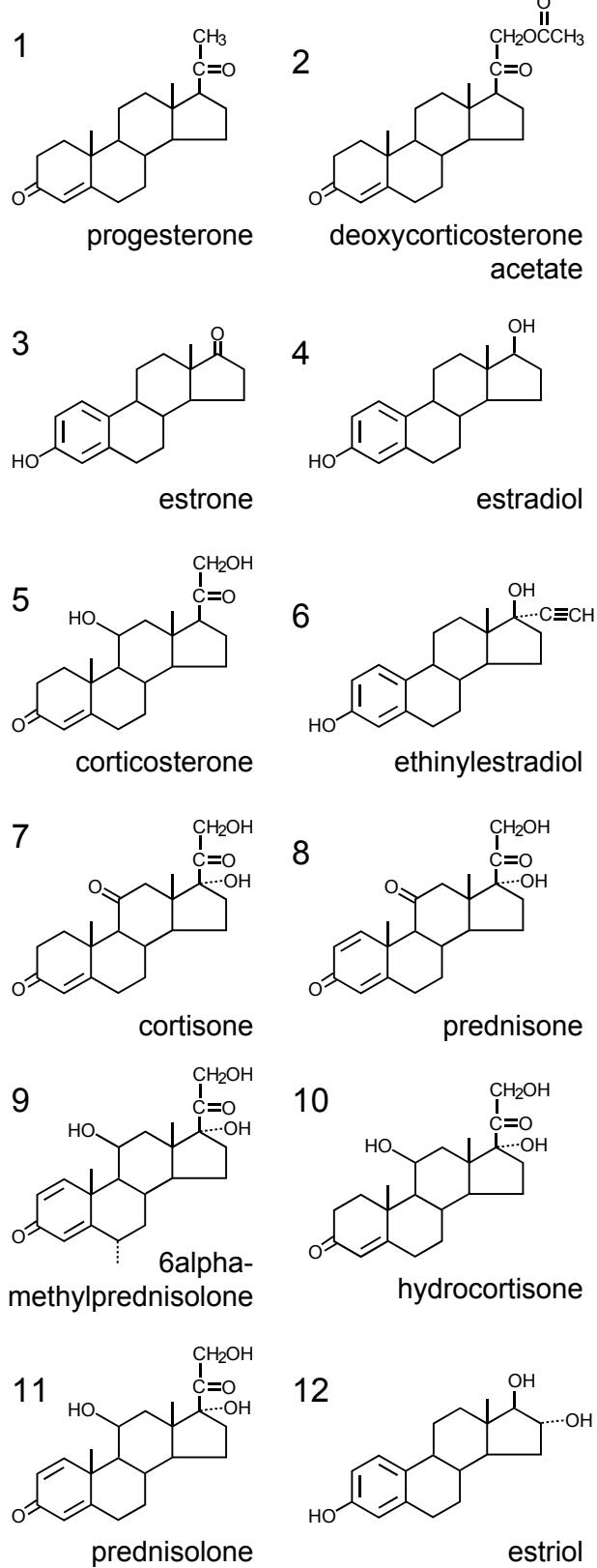

**Unison UK-C18, 250 x 4.6 mm**  
 A: water, B: acetonitrile  
 30-70% B (0-30 min), 0.8 mL/min  
 37 deg.C, 260 nm, 2 uL (1-2ug)