


Cadenza HS-C18

30 x 3 mm

Application

Cefalexin in CMC-Na Aqueous Solution

カルボキシメチルセルロース水溶液中のセファレキシ

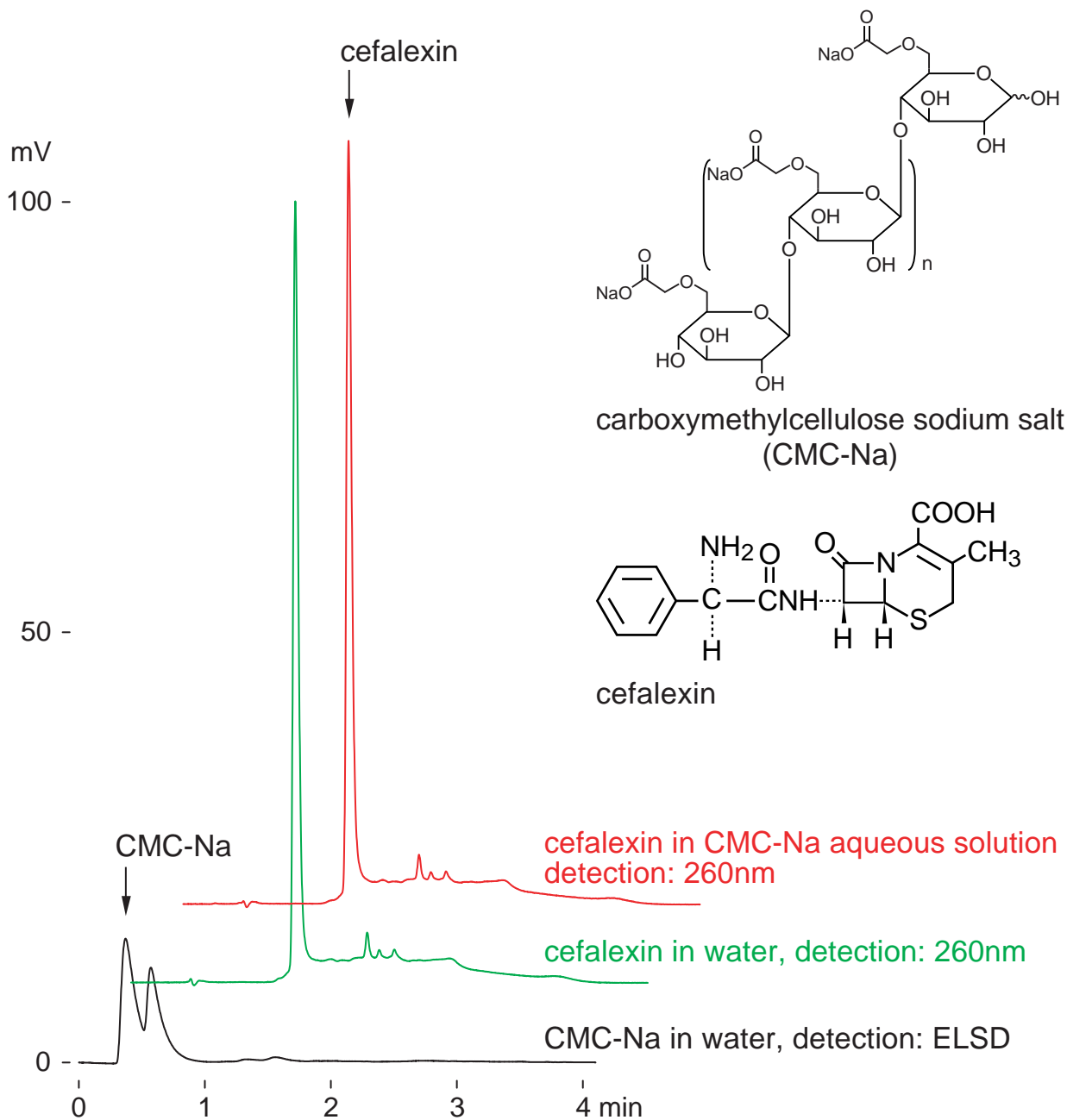

Cadenza HS-C18, 30 x 3mm

A: water /heptafluorobutyric acid = 100 /0.1

B: acetonitrile /heptafluorobutyric acid = 100 /0.1

0-70 %B (0-2min), 0.4 mL/min (3MPa), 37 deg.C

1 uL (cefalexin: 0.2ug, CMC-Na: 1ug)