



Cadenza CW-C18

250 x 3 mm

Application

Standard PA-Sugar Chains

ピリジルアミノ化糖標準品

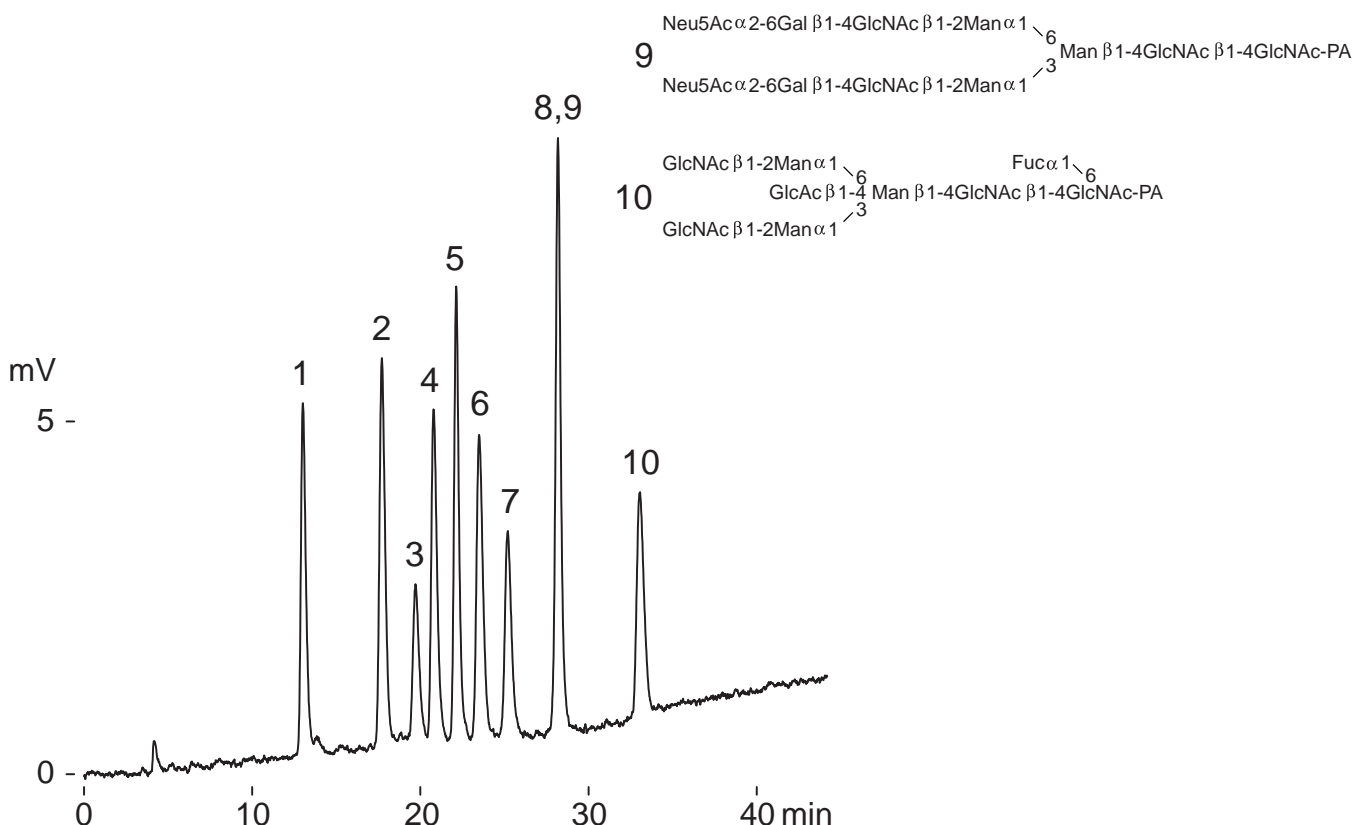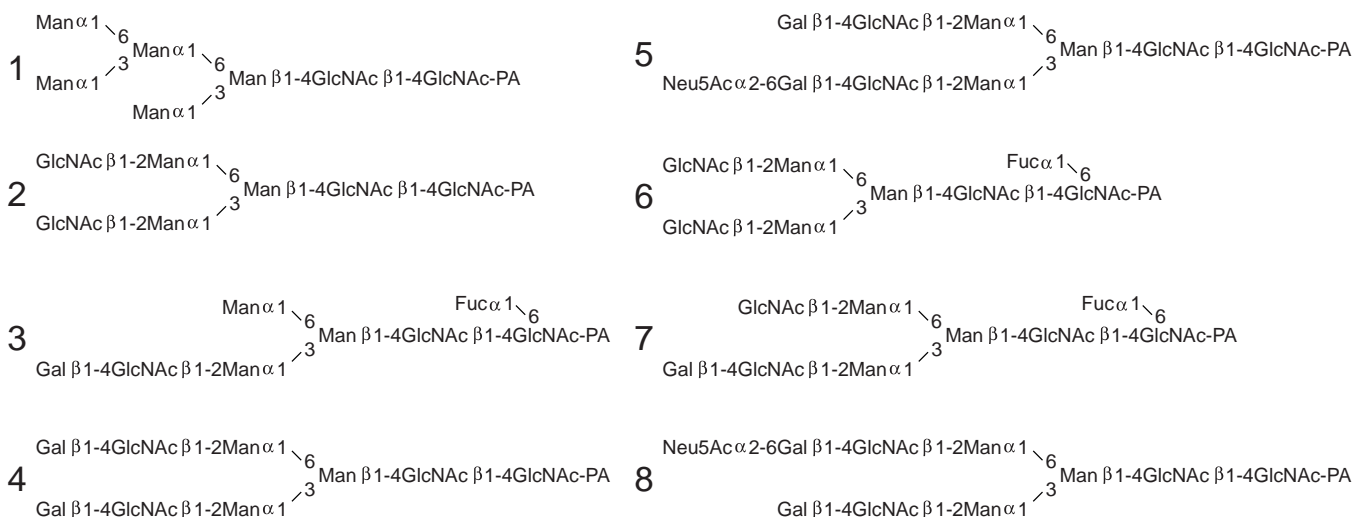

Cadenza CW-C18, 250 x 3 mm

A: water / acetic acid = 100 / 0.1, B: acetonitrile /acetic acid = 100 / 0.1

5-60%B (0-40min), 0.4 mL/min (17MPa), 37 deg.C, FLD (Ex: 300nm, Em: 400nm)

1 uL (0.3pmol)

Courtesy of Masuda Chemical Industries Co., LTD, Japan