



Unison UK-Amino

30 x 3 mm


Application

Adenosine phosphates (AMP, ADP, ATP)


アデノシンリン酸 (AMP, ADP, ATP)

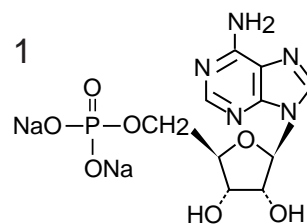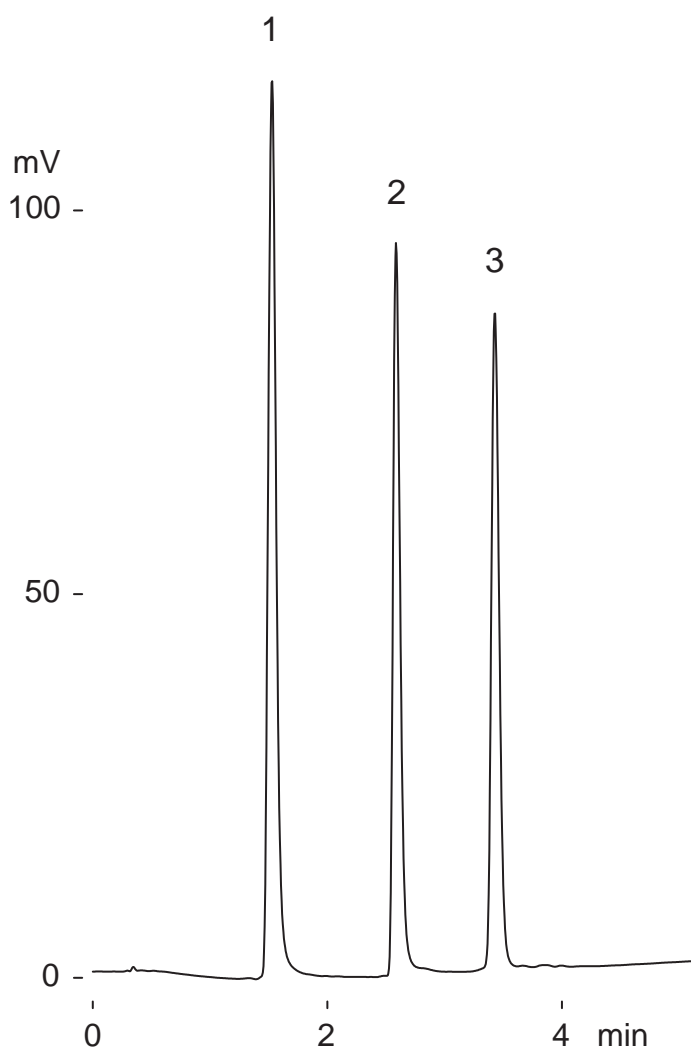

adenosine 5'- monophosphate disodium (AMP)

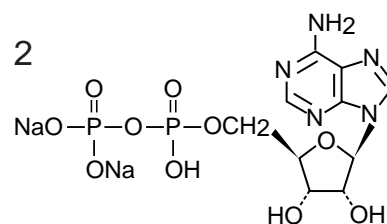

adenosine 5'- diphosphate disodium (ADP)

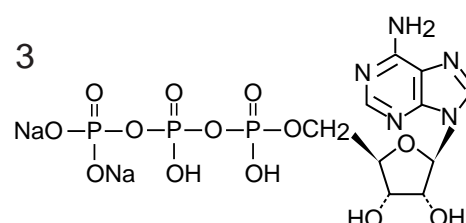

adenosine 5'- triphosphate disodium (ATP)

Unison UK-Amino, 30 x 3 mm

A: acetonitrile

B: water /triethylamine /acetic acid = 100 /4.1 /1.6
30-100%B (0-5min), 0.6 mL/min (2MPa), 50 deg.C
260 nm, 0.6 uL (1ug)