



Unison UK-Phenyl  
 Unison UK-C8  
 Unison UK-C18  
 Cadenza CL-C18  
 Cadenza CD-C18  
 Cadenza CW-C18

75 x 4.6 mm

Application

**Sudan Dyes**

スダン色素

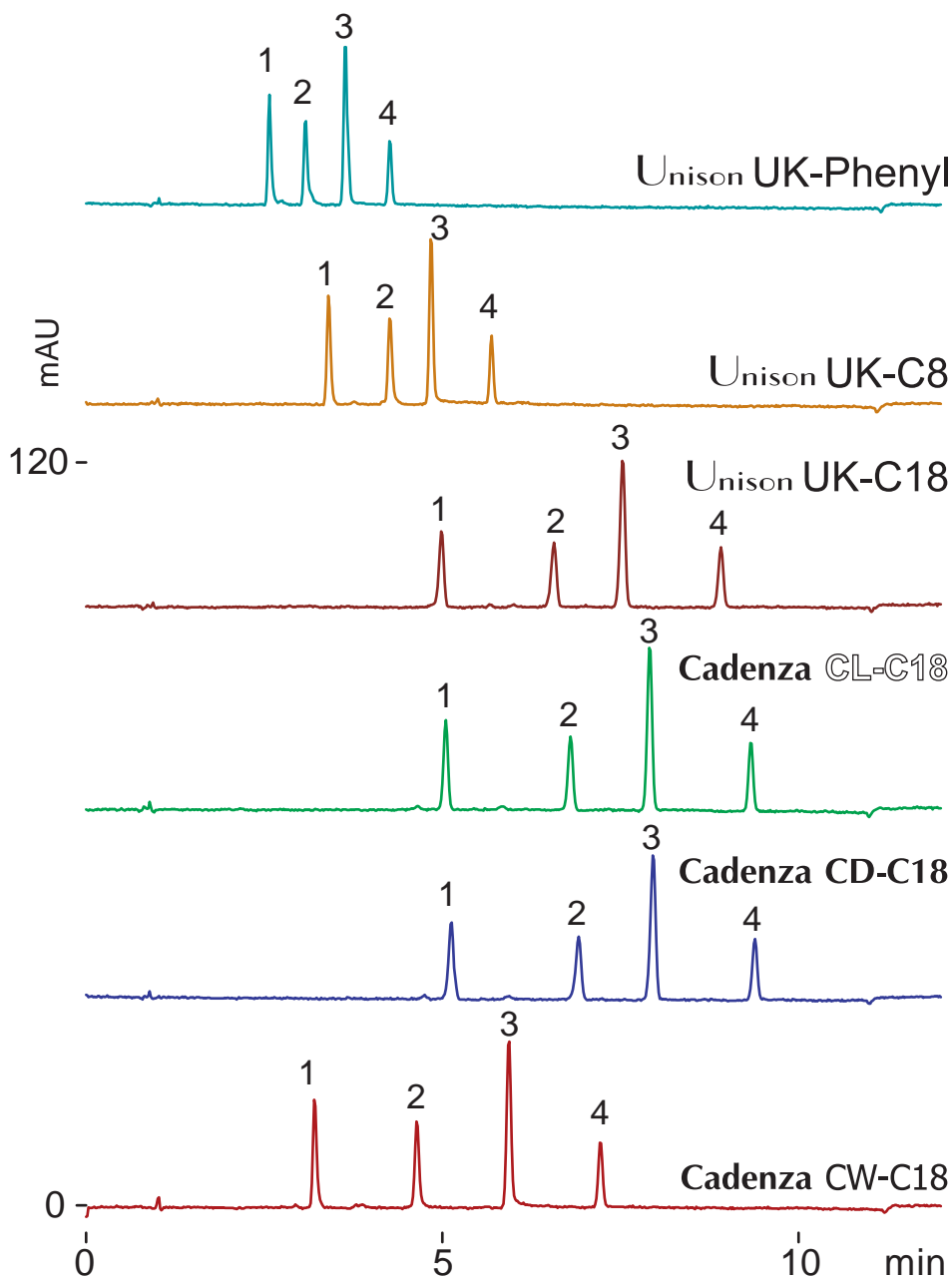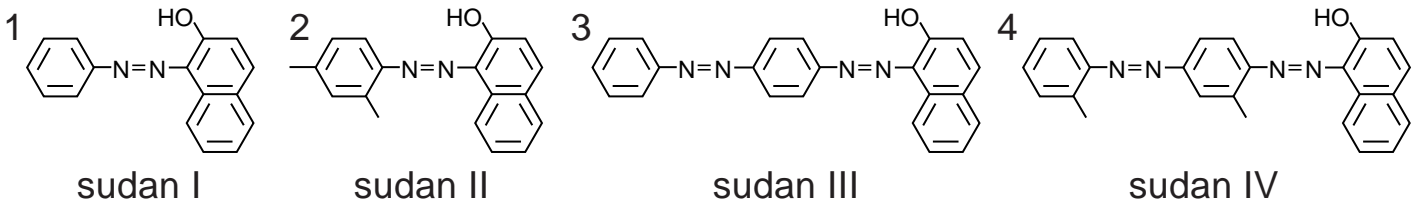

75 x 4.6 mm  
 A: water / formic acid  
 = 100 / 0.1  
 B: acetone / formic acid  
 = 100 / 0.1  
 60 - 90 %B (0-10 min)  
 1 mL/min (7-8 MPa)  
 37 deg.C, 500 nm  
 2 uL (40 ng)