


Cadenza CD-C18

150 x 4.6 mm

Application

LC-MS/MS Application for 25-Hydroxy Vitamin D3

25-ヒドロキシビタミンD3 のLC-MS/MS分析

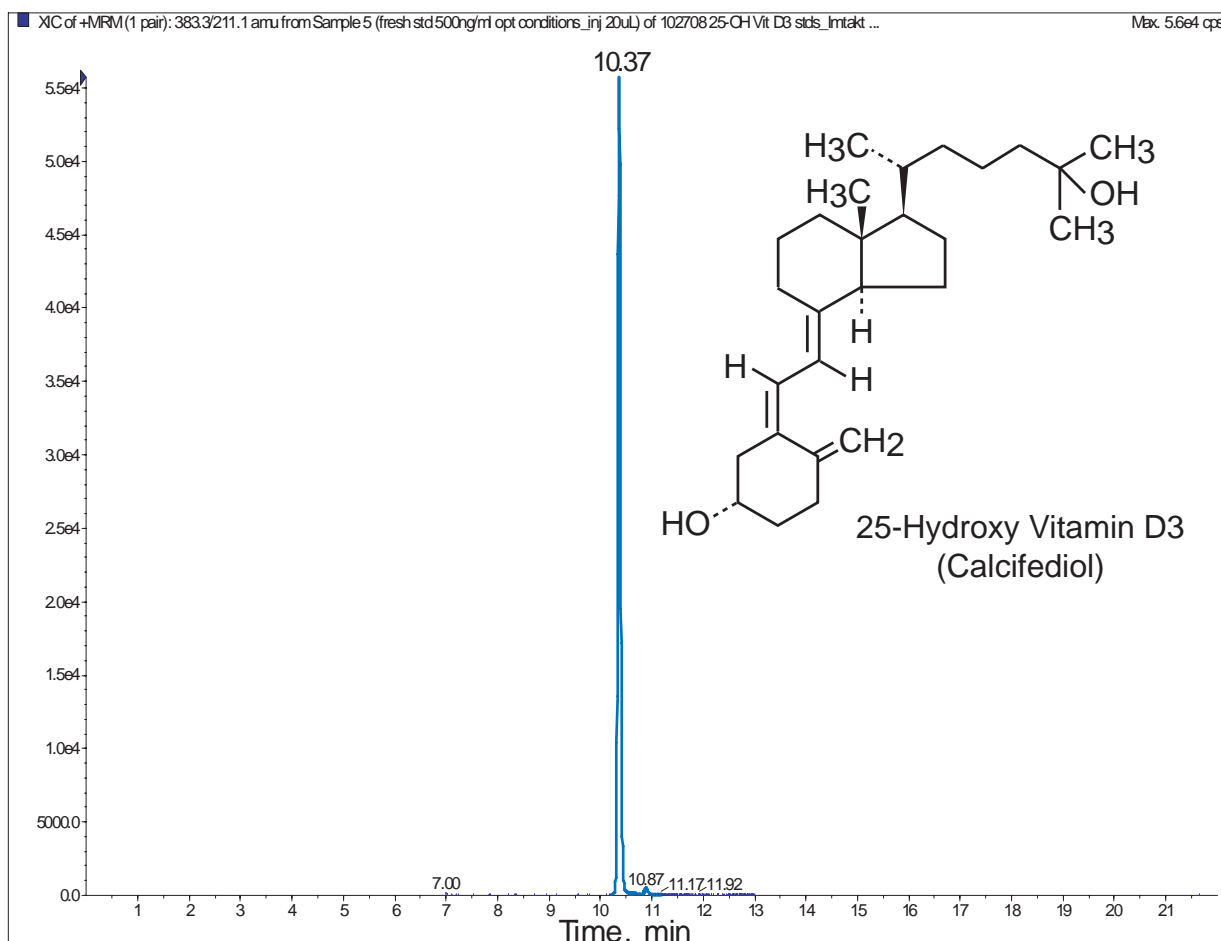

Cadenza CD-C18, 150 x 4.6 mm

A: water / methanol = 75 / 25

B: acetonitrile / ethyl acetate = 70 / 30

40%B (0-1min), 40-95%B (1-8min), 95%B (8-15min), 40%B (15-22min)

0.7mL/min, ambient

Applied Biosystems 4000 Q trap

Ionization: APCI Positive mode

LOQ: 1ng/mL

Courtesy of Penn State Hershey Medical Center, USA