


Scherzo SM-C18

75 x 4.6 mm

Application

Batch to batch reproducibility of Scherzo SM-C18

製造バッチ間の再現性 (Scherzo SM-C18)

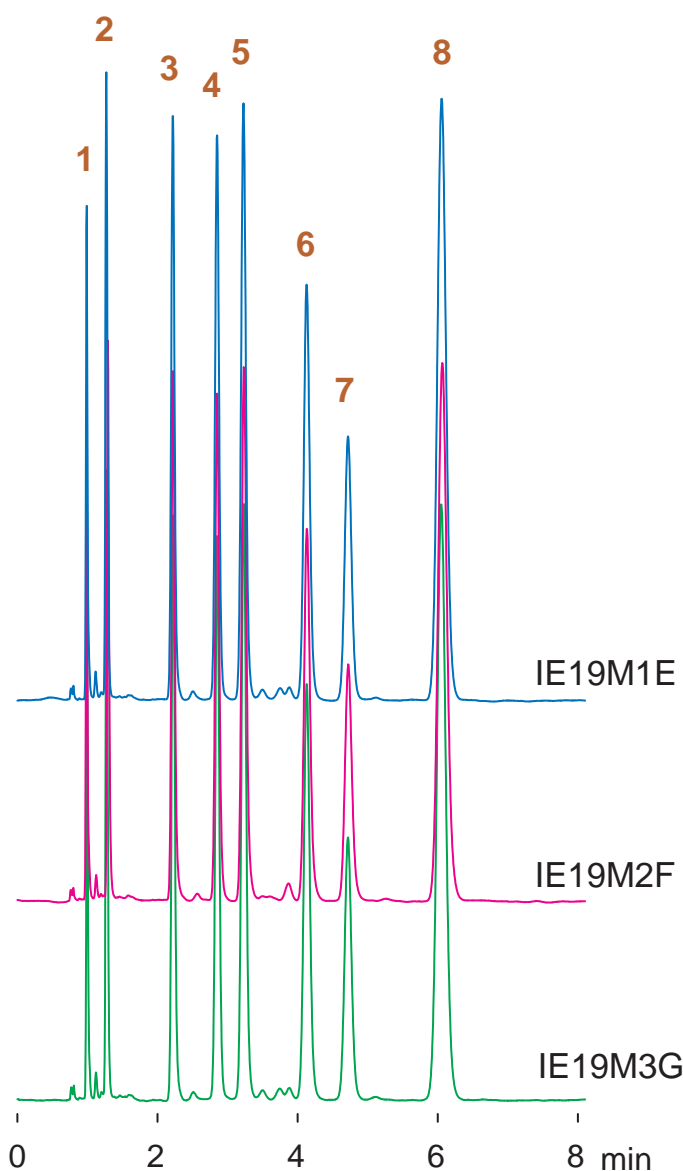

- 1 CC(=O)Nc1ccc(O)cc1 acetaminophen
- 2 Oc1nc2nncn2c1 1-hydroxy-7-azabenzotriazole
- 3 CC12CC[C@H]3[C@@H]([C@@H]1CC[C@@H]2O)C(=O)O prednisolone
- 4 CC(=O)Nc1ccsc1 methyl 3-amino-2-thiophenecarboxylate
- 5 CC12CC[C@H]3[C@@H]([C@@H]1CC[C@@H]2O)C(=O)O 6-alpha-methylprednisolone
- 6 CC12CC[C@H]3[C@@H]([C@@H]1CC[C@@H]2O)C(=O)O corticosterone
- 7 NC1=CC=C(C=C1)C(=O)c2ccccc2 4-aminobenzophenone
- 8 CCCOC(=O)c1ccc(O)cc1 propylparaben

Scherzo SM-C18, 75 x 4.6 mm
 10 mM ammonium acetate / acetonitrile = 65 / 35
 1 mL/min, 37 deg.C, 254 nm