
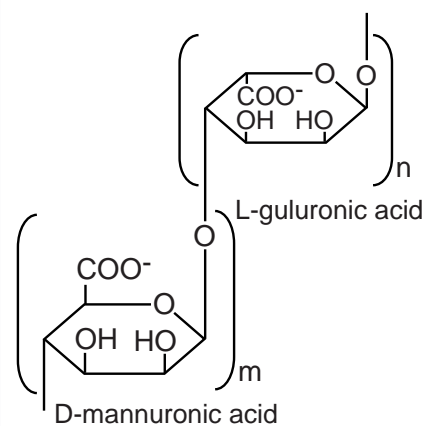


Presto FF-C18

150 x 4.6 mm

Application

Alginic acid
アルギン酸

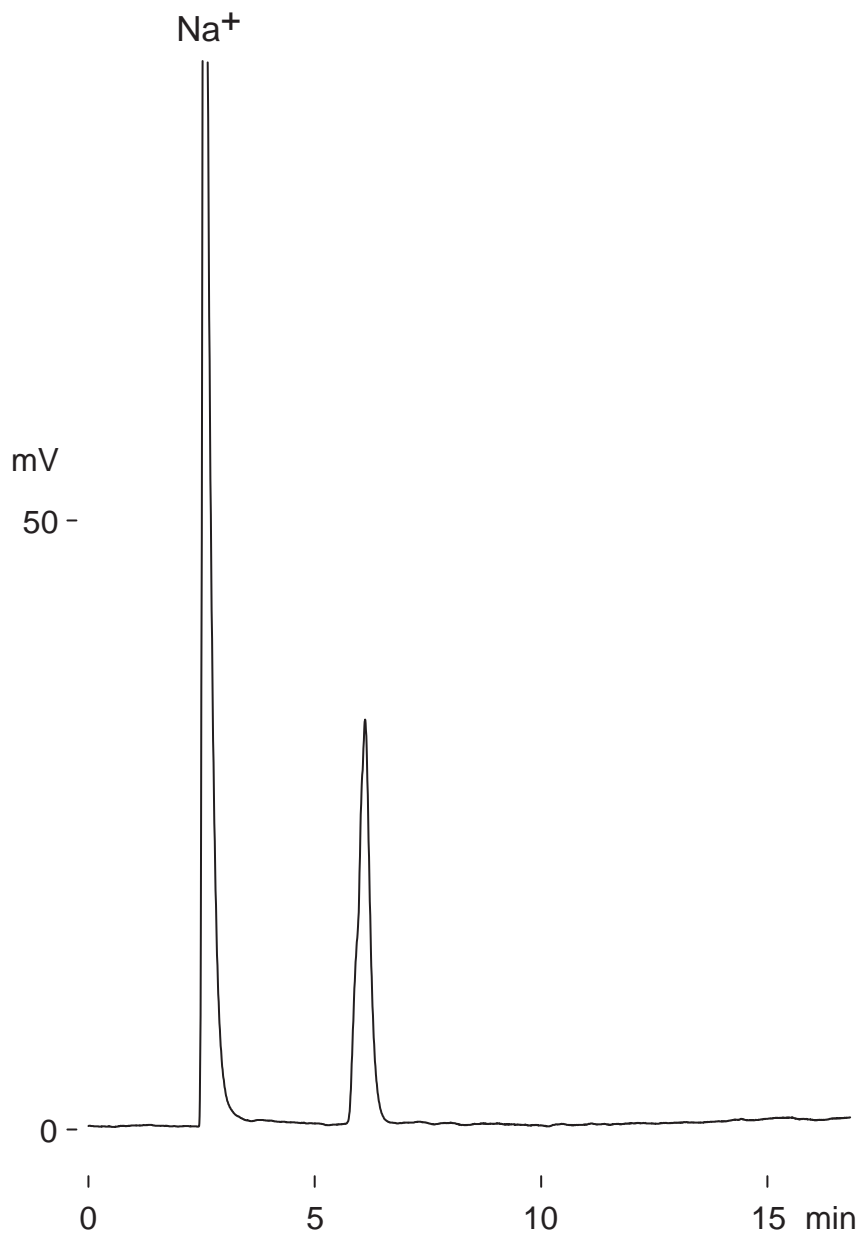

alginic acid sodium salt
(ca.200 kDa)

Presto FF-C18, 150 x 4.6mm

A: water / trifluoroacetic acid = 100 / 0.1

B: acetonitrile / trifluoroacetic acid = 100 / 0.1

0-90 %B (0-15min), 0.4 mL/min (23MPa), 37 deg.C

ELSD (spray chamber 50deg.C, drift tube 100deg.C)

5 uL (25ug)