


Presto FF-C18

150 x 4.6 mm


Application

### Hyaluronic acid

ヒアルロン酸(鶏冠由来)

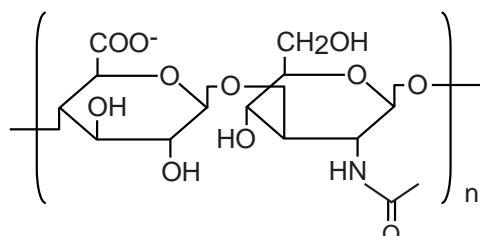

hyaluronic acid sodium salt  
(from rooster comb, ca.1 MDa)

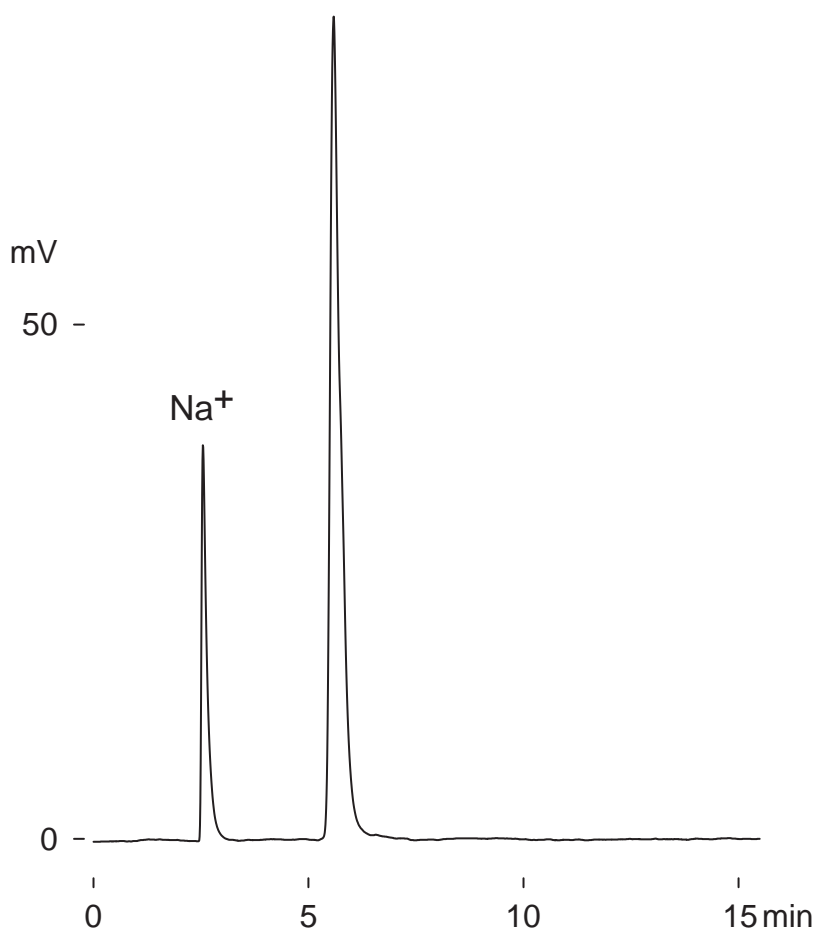

Presto FF-C18, 150 x 4.6mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-90 %B (0-15min), 0.5 mL/min (23MPa), 37 deg.C

ELSD (spray chamber 50deg.C, drift tube 100deg.C)

3 uL (7.5ug)