



Presto FF-C18

30 x 4.6 mm
150 x 4.6 mm

Application

Mannan

マンナン(出芽酵母由来)

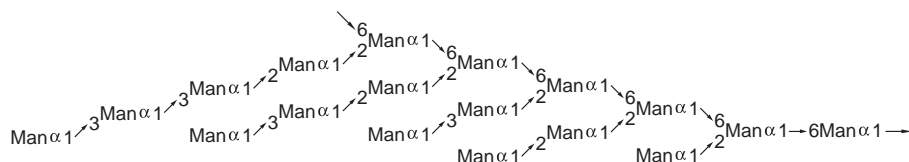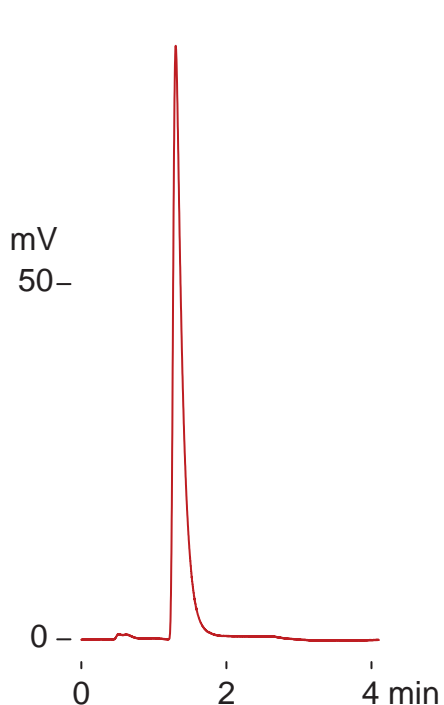

mannan
(from *Saccharomyces cerevisiae*, ca.130 kDa)

Presto FF-C18, 30 x 4.6 mm

A: water


B: acetonitrile

0-70 %B (0-2min)

0.5 mL/min (5 MPa), 37 deg.C,

ELSD (spray chamber 50 deg.C, drift tube 100 deg.C)

1 uL (5 ug)

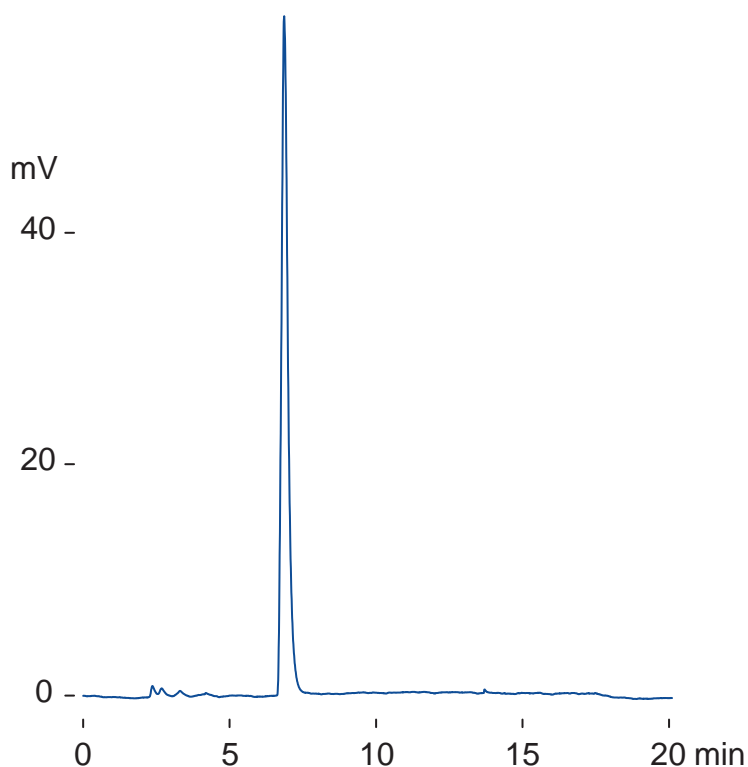

Presto FF-C18, 150 x 4.6 mm

A: water

B: acetonitrile

0-50 %B (0-15min)

0.4 mL/min (17 MPa), 37 deg.C

ELSD (spray chamber 50 deg.C
drift tube 100 deg.C)

2 uL (10 ug)