


Presto FF-C18

150 x 4.6 mm

Application


PEGylated protein (2)

PEG化タンパク質

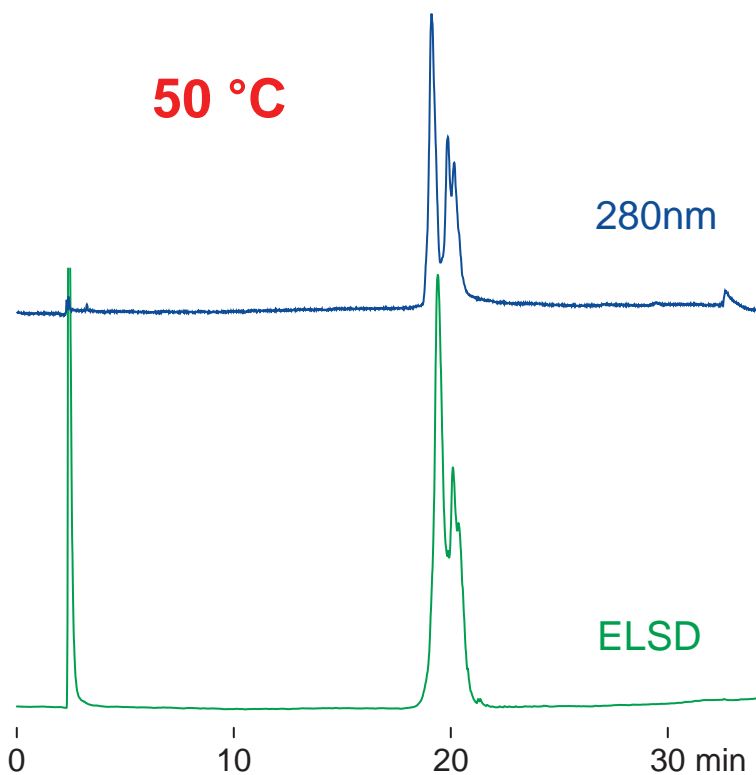

PEGylated protein, 60 kDa
(20 kDa protein + 40 kDa PEG)

Presto FF-C18, 150 x 4.6 mm
A: water /trifluoroacetic acid
= 100 / 0.1
B: acetonitrile /trifluoroacetic acid
= 100 / 0.1
20-70 %B (0-30min)
0.4 mL/min (15 MPa)
1 μ L

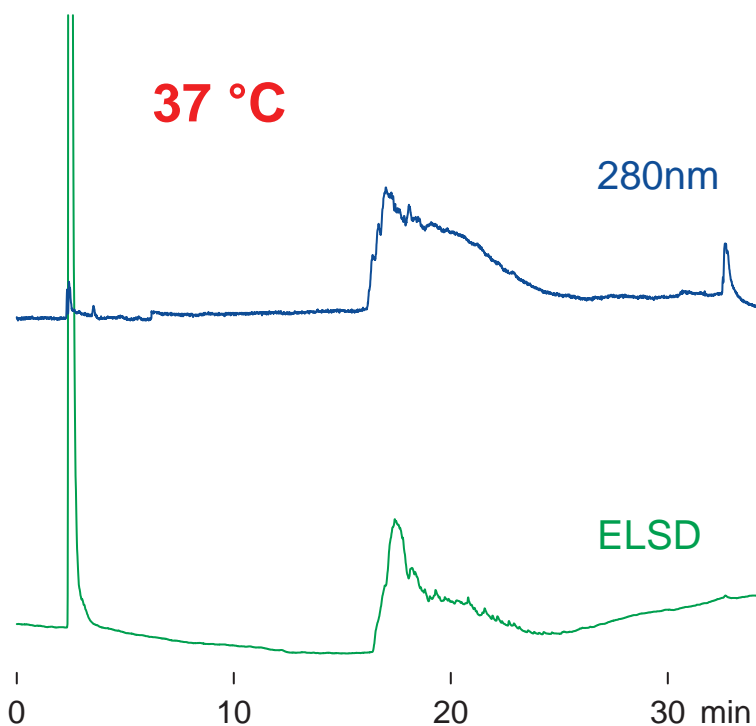

PEGylated protein, 60 kDa
(20 kDa protein + 40 kDa PEG)

Presto FF-C18, 150 x 4.6 mm
A: water /trifluoroacetic acid
= 100 / 0.1
B: acetonitrile /trifluoroacetic acid
= 100 / 0.1
20-70 %B (0-30min)
0.4 mL/min (19 MPa)
1 μ L