


Presto FF-C18

50 x 4.6 mm
250 x 4.6 mm

Application

The effect of gradient slope and column length for alpha-lactalbumin separation

α -ラクトアルブミン の分離に与えるグラジエント勾配とカラム長さの影響

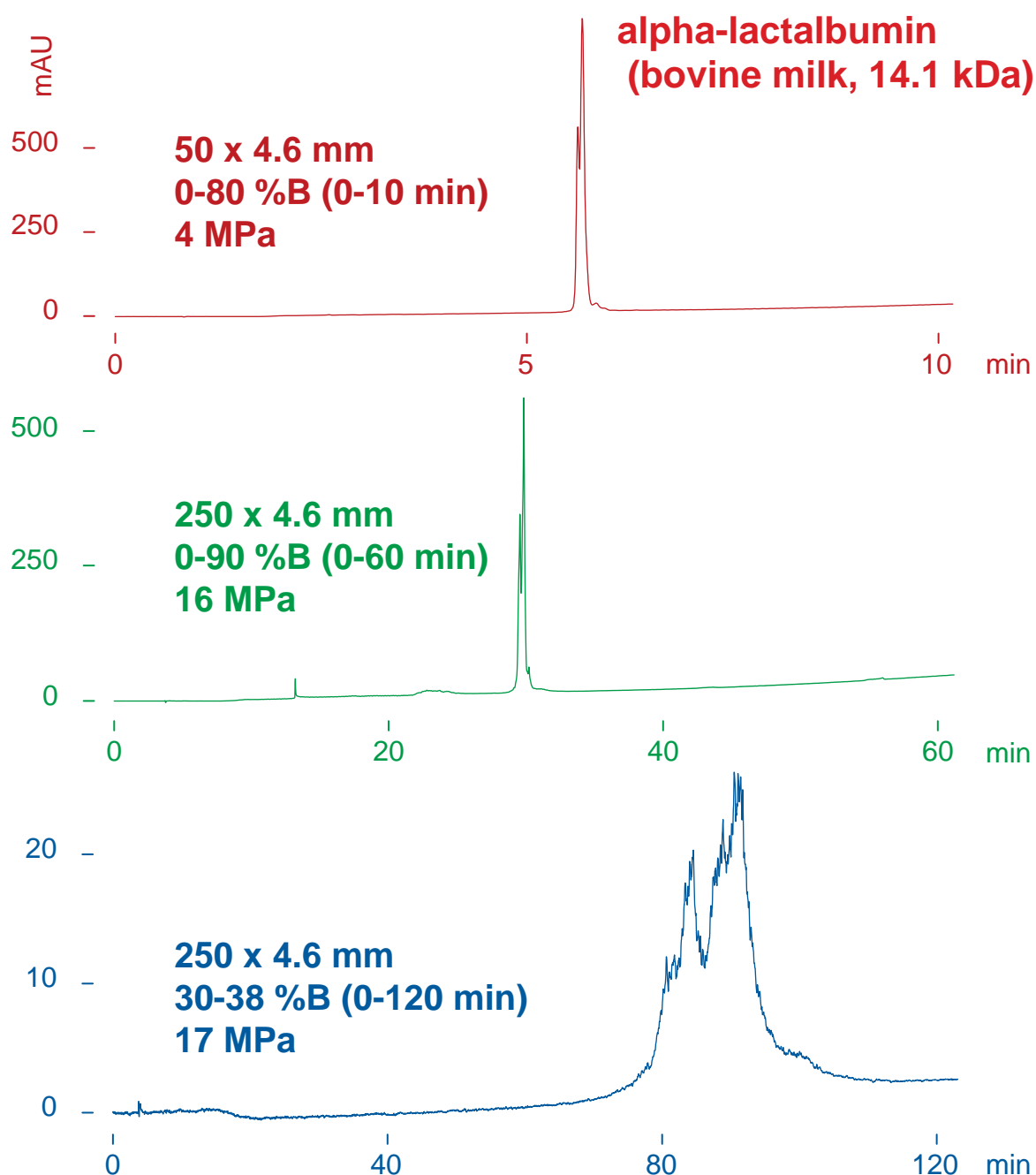

Presto FF-C18

A: water / trifluoroacetic acid = 100 / 0.1

B: acetonitrile / trifluoroacetic acid = 100 / 0.1

0.4 mL/min, 80 deg.C, 220 nm, 0.4-1.4 μ L (4-14 μ g)