


Scherzo SM-C18
 Scherzo SW-C18
 Unison UK-C18

250 x 3 mm

Application

Sialylglycopeptide in egg yolk extract

卵黄抽出物中のシアリルグリコペプチド

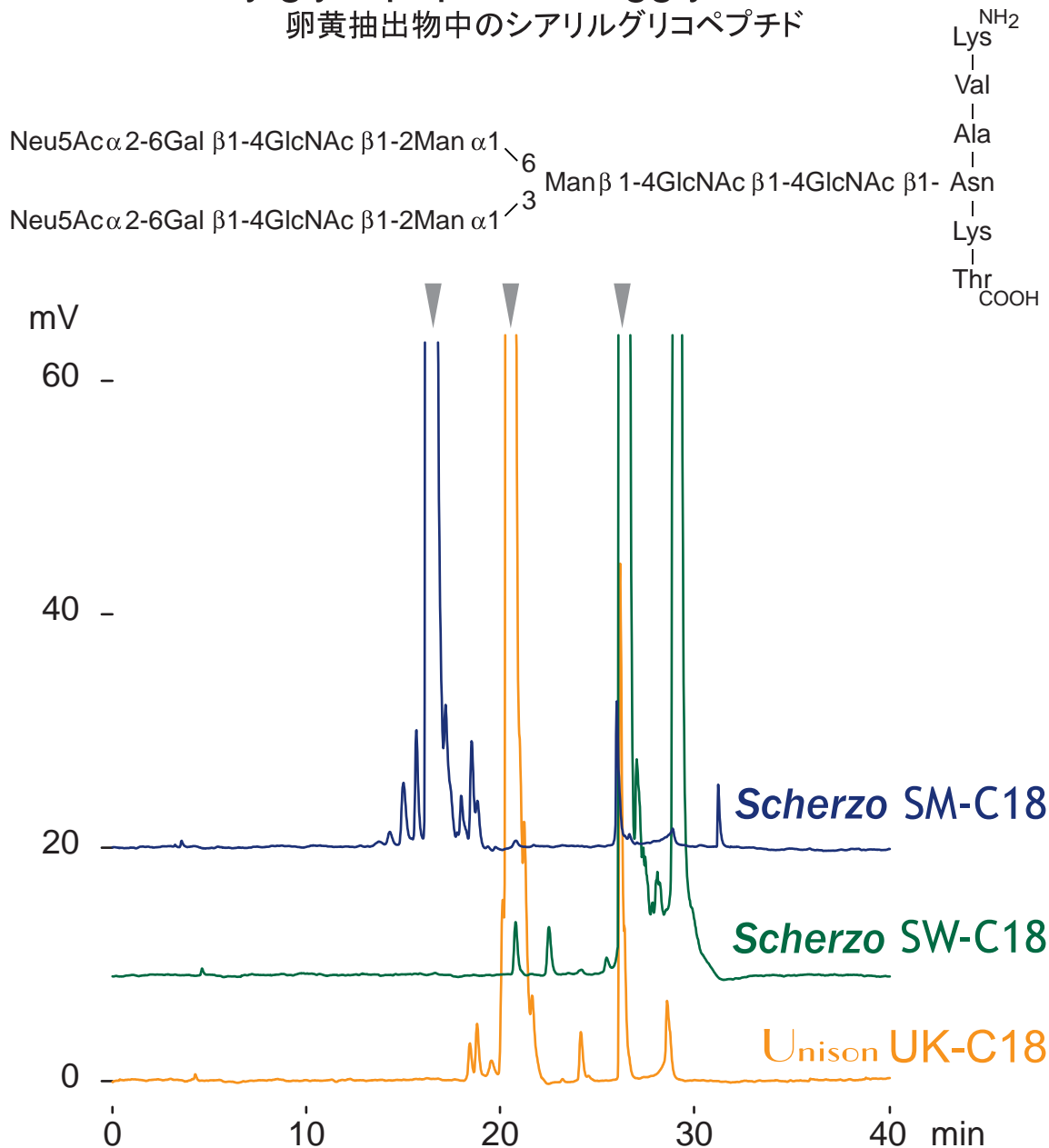

250 x 3 mm

A: water /HCOOH /TFA = 100 /0.5 /0.02

B: acetonitrile /HCOOH /TFA = 100 /0.5 /0.02

0-5 %B(0-20min), 5-90 %B(20-25min), 0 %B(25-40min)

0.3 mL/min (8-9MPa), 37 deg.C

ELSD (spray chamber 40 deg.C, drift tube 100 deg.C)

5 μ L

Courtesy of FUSHIMI Pharmaceutical Co., LTD.