

Cadenza CD-C18 HT



150 x 4.6 mm

Application


Comparison between 3um and sub-2um ODS columns

3um と sub-2um ODSカラムの比較

**3um Cadenza CD-C18 HT, 150 x 4.6 mm**

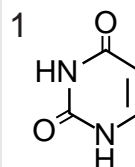
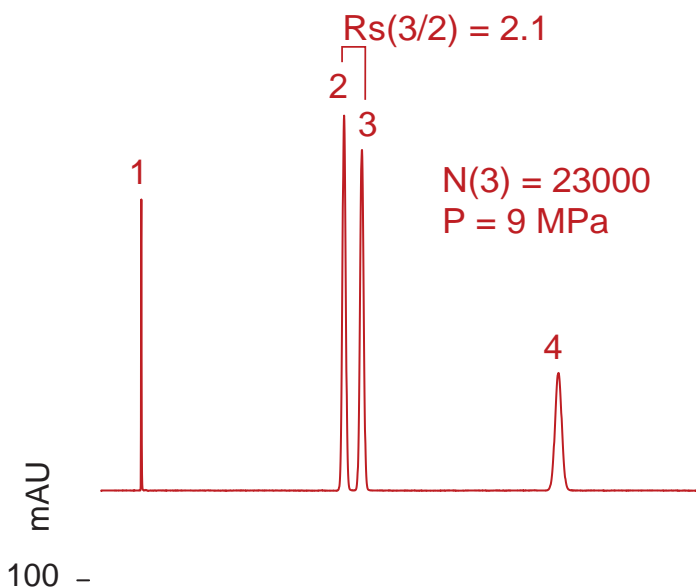


uracil

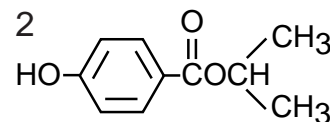


isopropylparaben

**1.7um ODS, 100 x 4.6 mm**

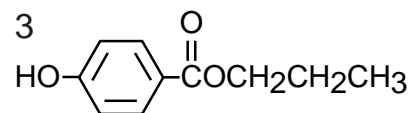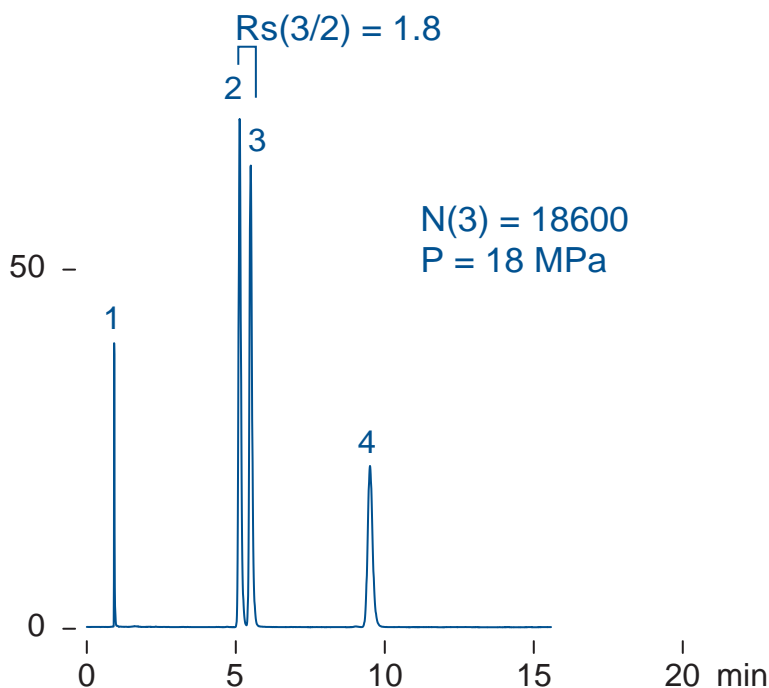

propylparaben

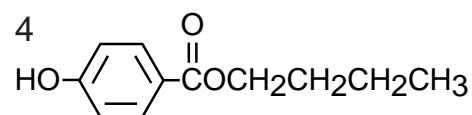

butylparaben

water / acetonitrile / acetic acid = 60 / 40 / 0.1  
 1 mL/min, 37 deg.C, 260 nm, 1 uL ( 0.02- 0.16ug)