


Scherzo SW-C18

150 x 3 mm

Application

Aliphatic amines

脂肪族アミン

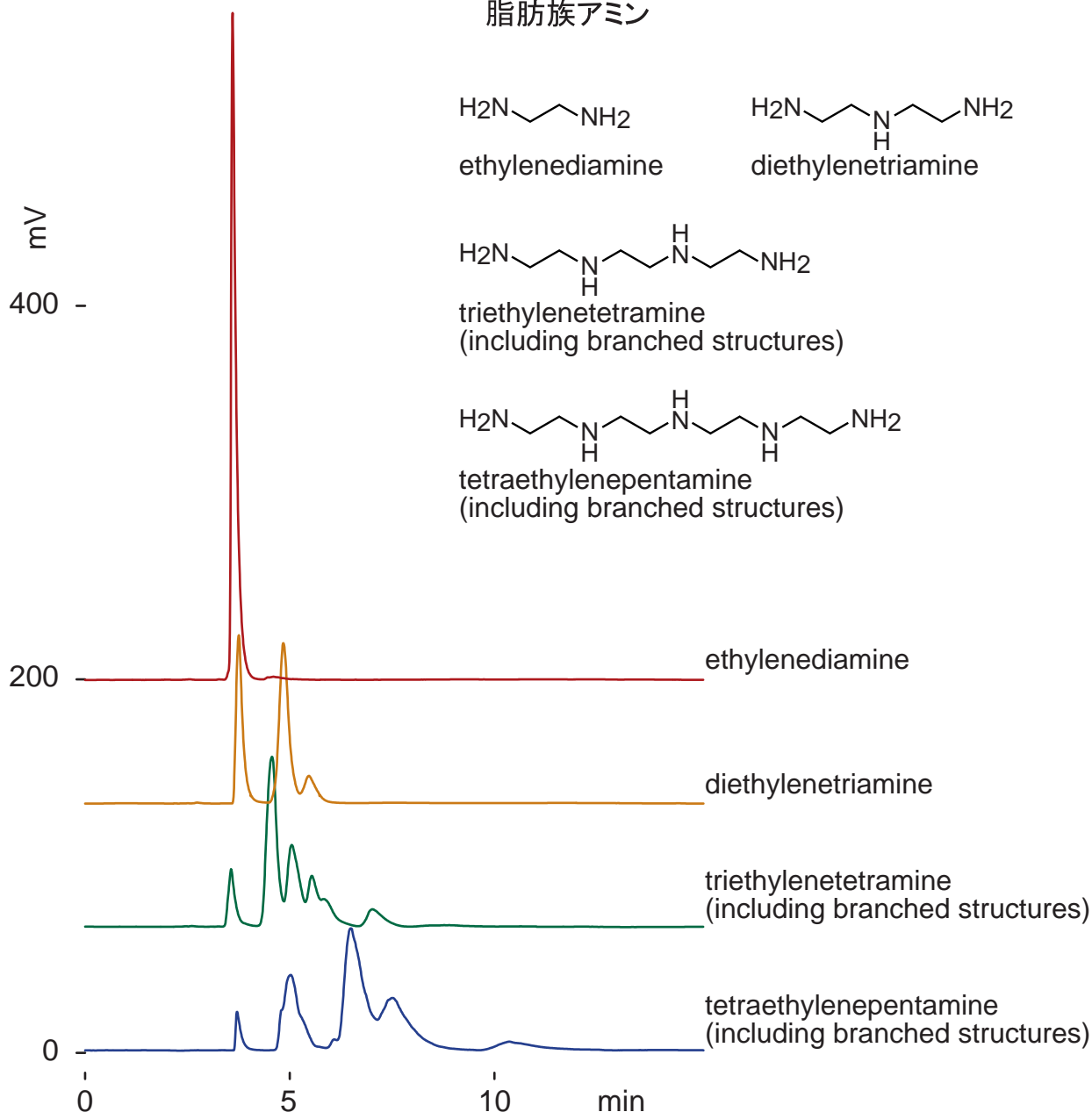

Scherzo SW-C18, 150 x 3 mm

A: 5mM acetic acid / tetrahydrofuran = 90 / 10

B: 100mM ammonium acetate / tetrahydrofuran = 90 / 10

0-80 %B (0-10min)

0.3 mL/min (8MPa), 37 deg.C

ELSD (spray chamber 10 deg.C, drift tube 40 deg.C)

1-2 uL (5-10nL, 50mM acetic acid)