


Scherzo SW-C18
Unison UK-C18

150 x 3 mm

Application

Oligopeptides (angiotensins)

オリゴペプチド（アンジオテンシン類）

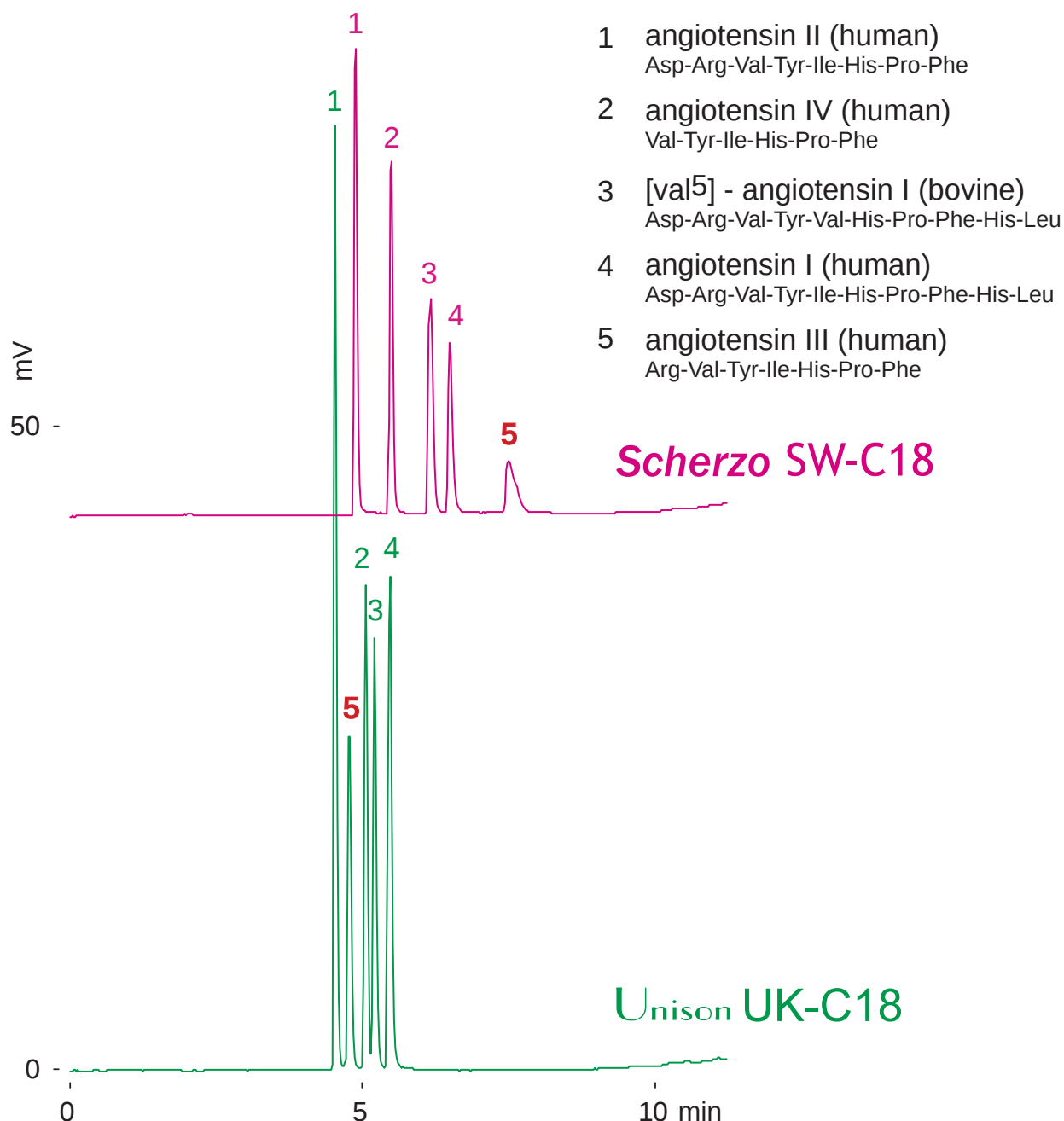

150 x 3 mm

A: 10 mM ammonium formate,

B: 50 mM ammonium formate / acetonitrile = 25 / 75

20 - 90 %B (0 - 10 min)

0.4 mL/min (11MPa), 37 deg.C, ELSD (spray chamber 50 deg.C, drift tube 100 deg.C)

5.0 uL (1 ug)