


Unison UK-Phenyl
Unison UK-C18

150 x 3 mm


Application

Alkylbenzenesulfonates

アルキルベンゼンスルホン酸ナトリウム

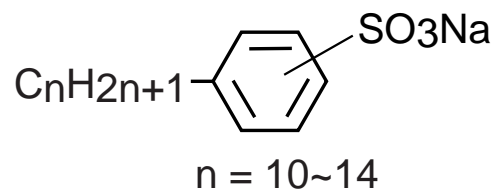

sodium alkylbenzenesulfonates
(including branched structure)

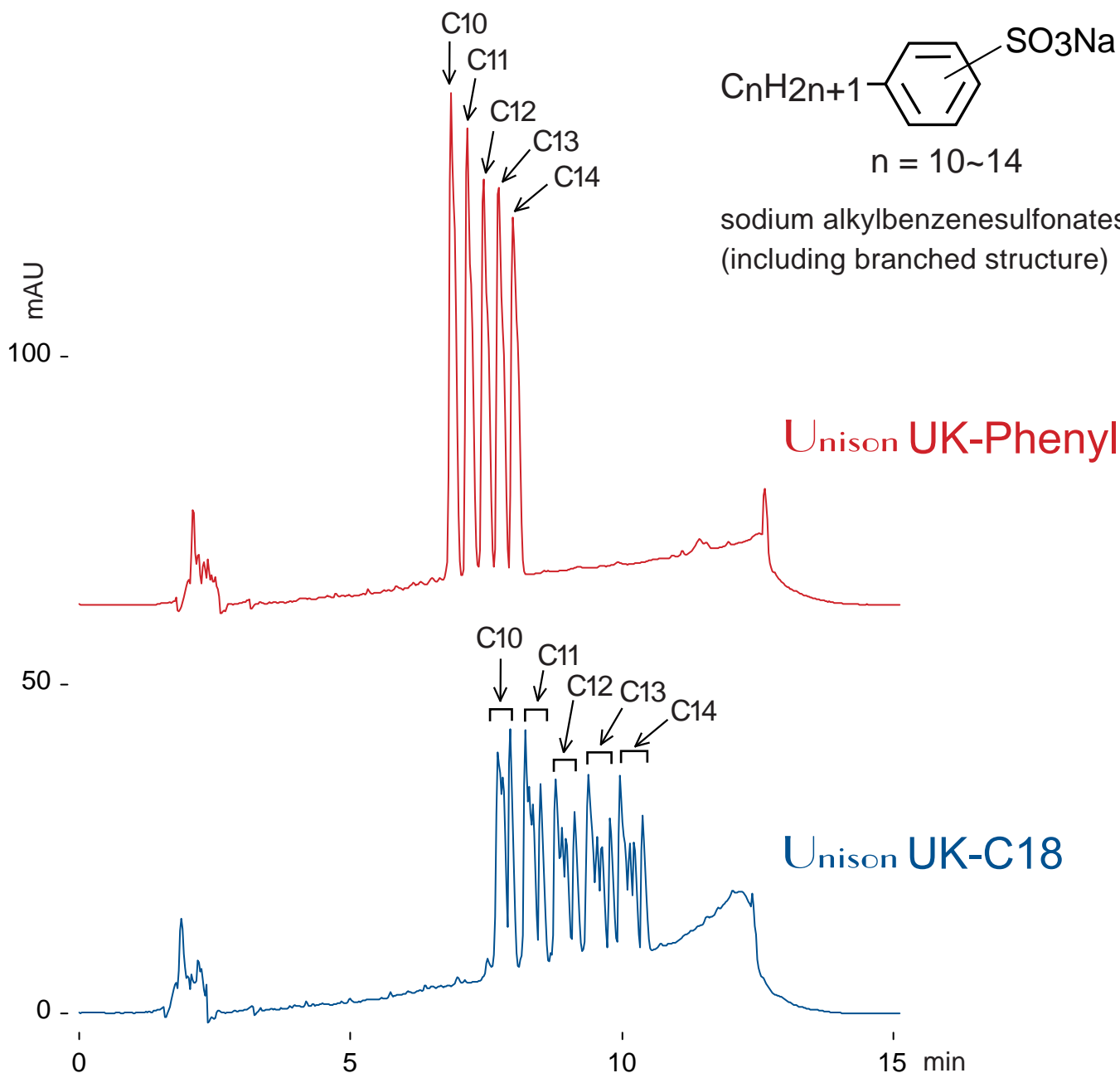

150 x 3 mm

A: 50 mM ammonium formate

B: acetonitrile / tetrahydrofuran = 80 / 20

25-95 %B (0-10 min)

0.4 mL/min (9-10 MPa), 37 deg.C, 260 nm, 15 uL (8 ug, 50% MeOH)