

Cadenza CW-C18

Unison UK-C8

Unison UK-Phenyl

Scherzo SM-C18








Unison UK-C18

50 x 3 mm

Application

LC-MS: Thyroid hormone

甲状腺ホルモン

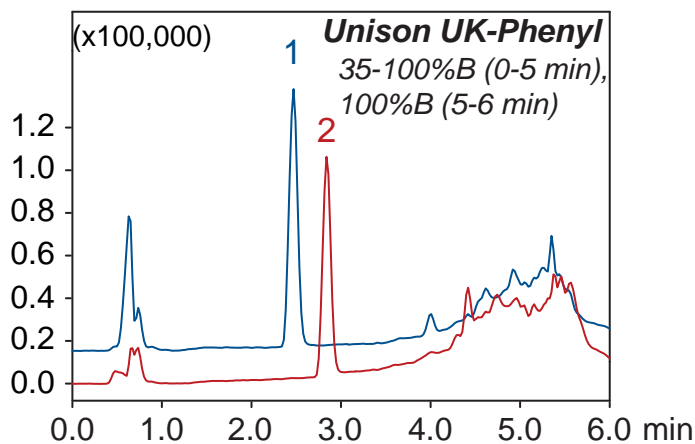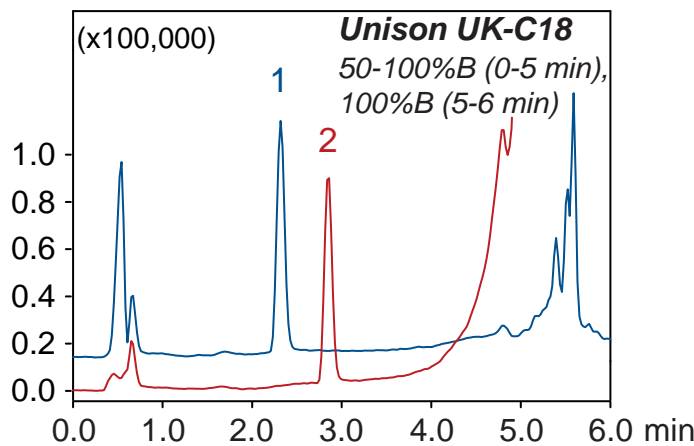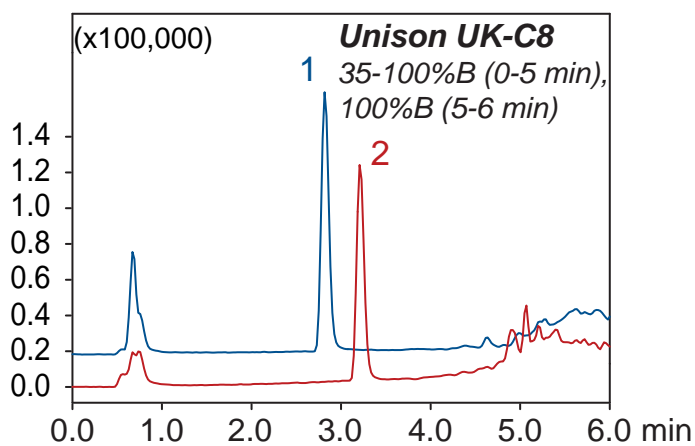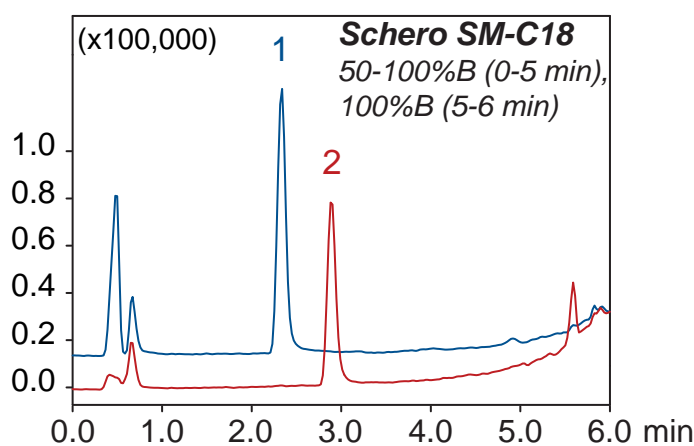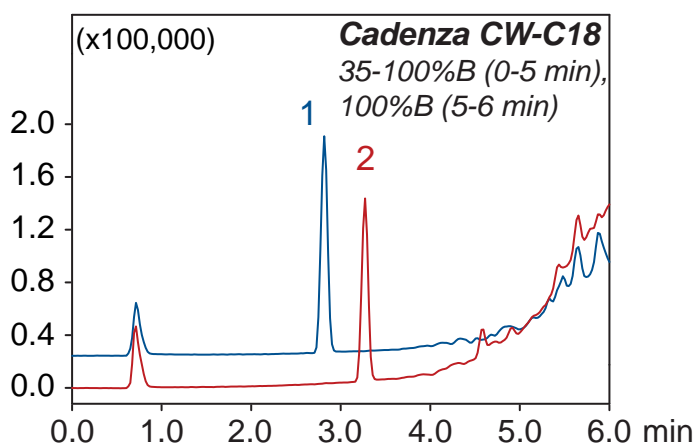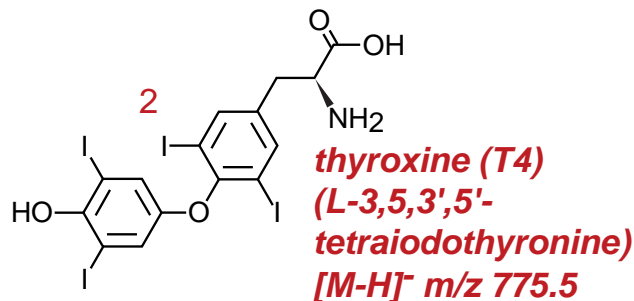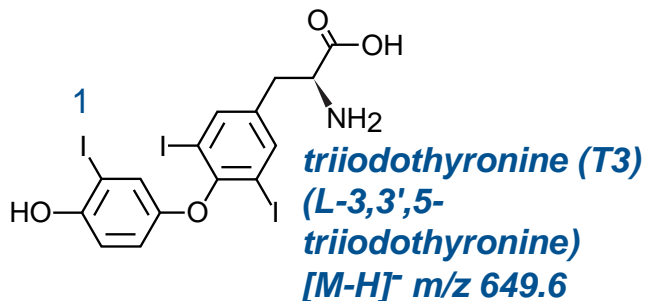

50 x 3 mm
A: water / AcOH = 100 / 0.1
B: MeOH / AcOH = 100 / 0.1
0.4 mL/min (6-7 MPa), 37°C
20 µL (50ng/mL)
MeOH / water / AcOH = 50 / 50 / 0.1
Single Quad. MS (ESI)