


Unison UK-Silica

250 x 4.6 mm

Application

DMF (dimethyl formamide), DMSO (dimethyl sulfoxide)

DMF (ジメチルホルムアミド), DMSO (ジメチルスルホキシド)

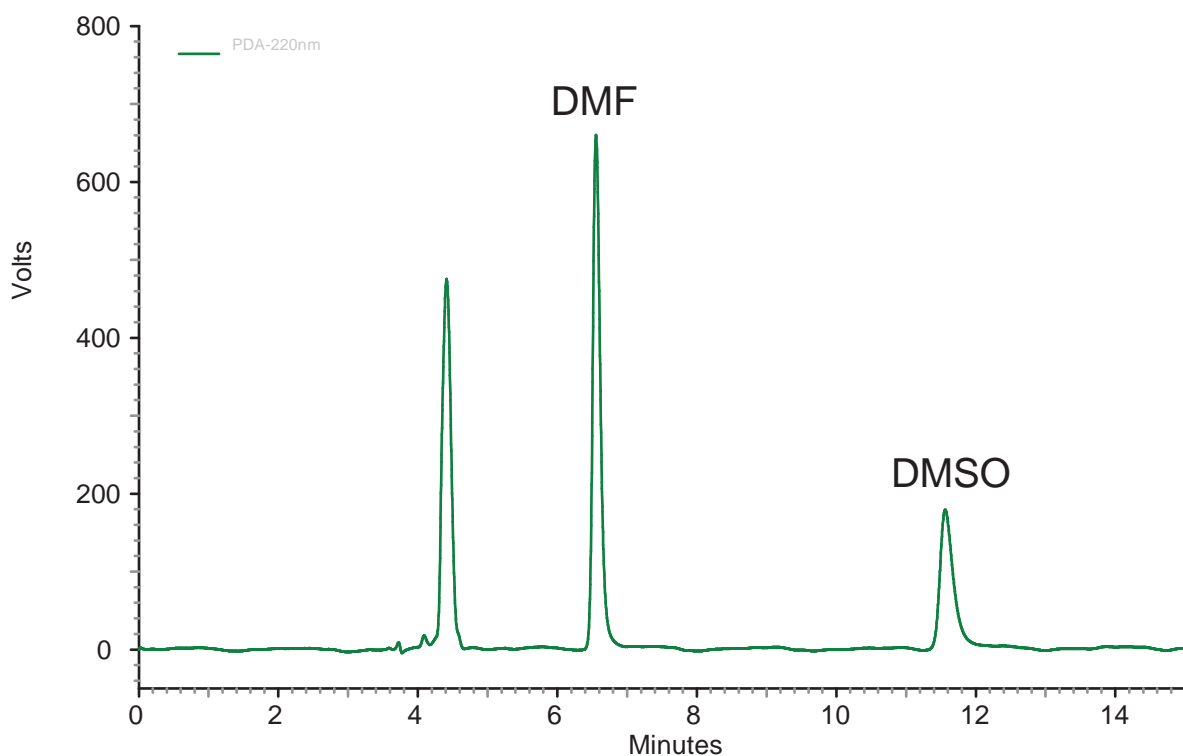

Unison UK-Silica, 250 x 4.6 mm  
hexane / ethanol = 50 / 50  
1.0 mL/min, 45°C  
220nm  
10  $\mu$ L (10ppm for DMF, 50ppm for DMSO)

Courtesy of PEAKMAN SP, KOREA