




Scherzo SM-C18

150 x 4.6 mm

Application

**Histamine**

ヒスタミン

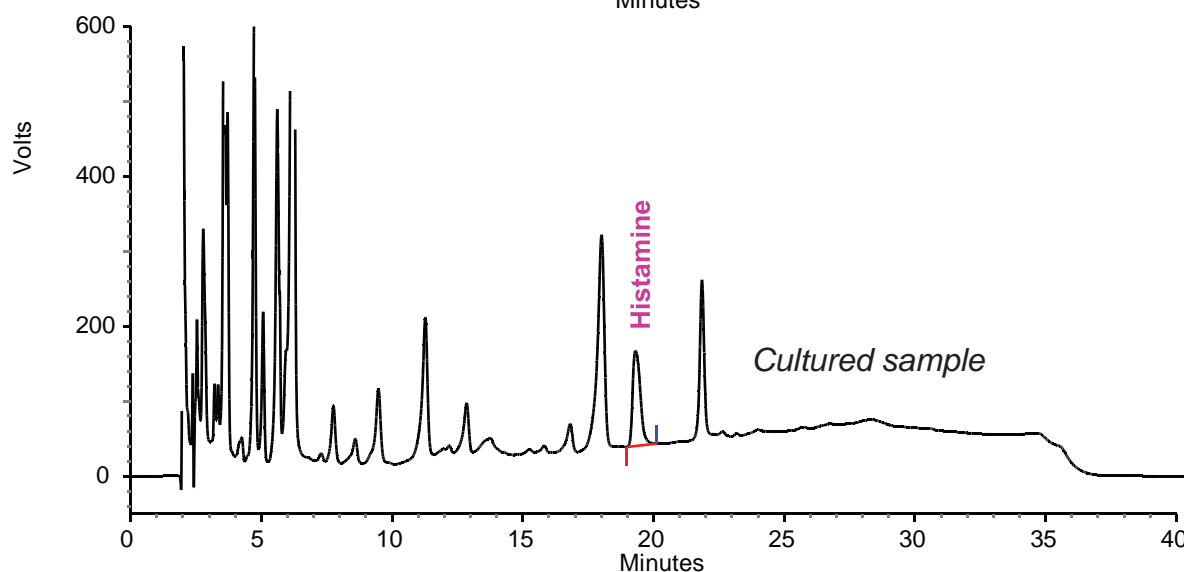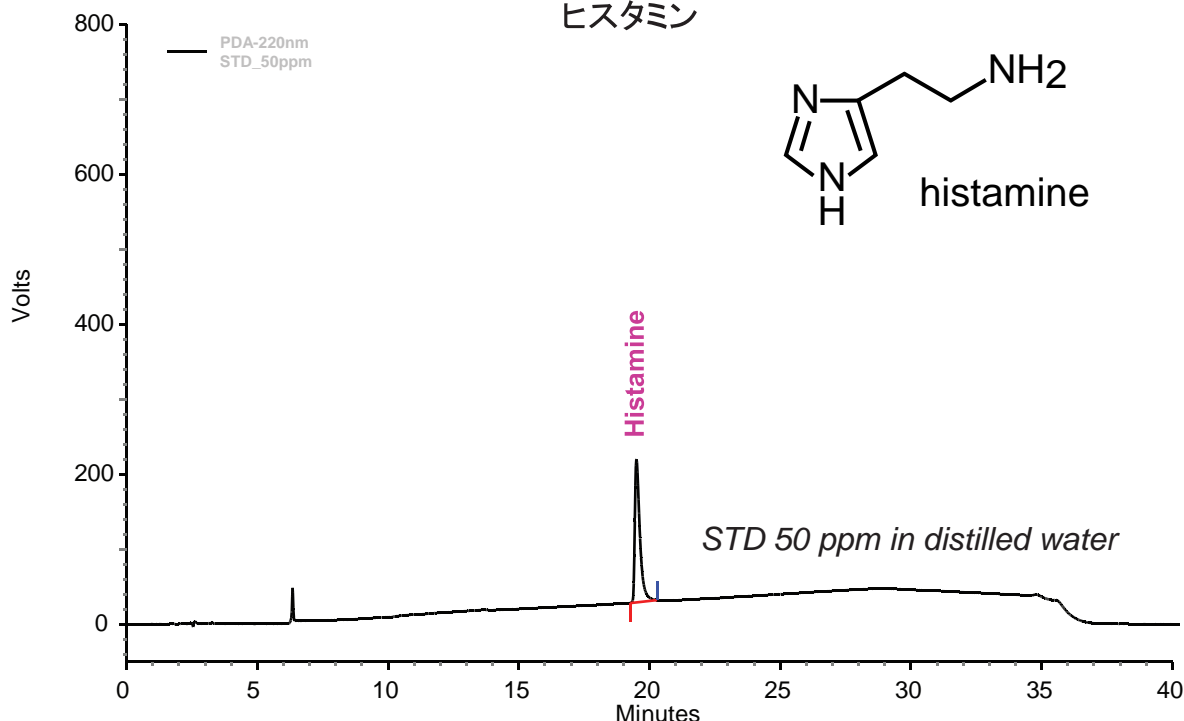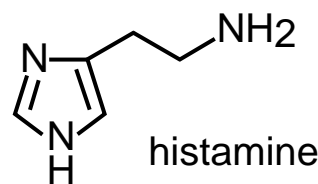

Scherzo SM-C18, 150 x 4.6 mm

A: 10mM KH<sub>2</sub>PO<sub>4</sub>

B: 50mM phosphate buffer (pH 7.4) / acetonitrile = 95 / 5

0-80%B (0-25min), 80%B (25-32min), 0%B (32-40min)

0.8 mL/min, 40°C

10 μL

220nm

Courtesy of PEAKMAN SP, KOREA