



Scherzo SM-C18

150 x 4.6 mm

Application

### S-allyl-L-cysteine in Garlic Extract

ニンニク抽出液中のアリルシステイン

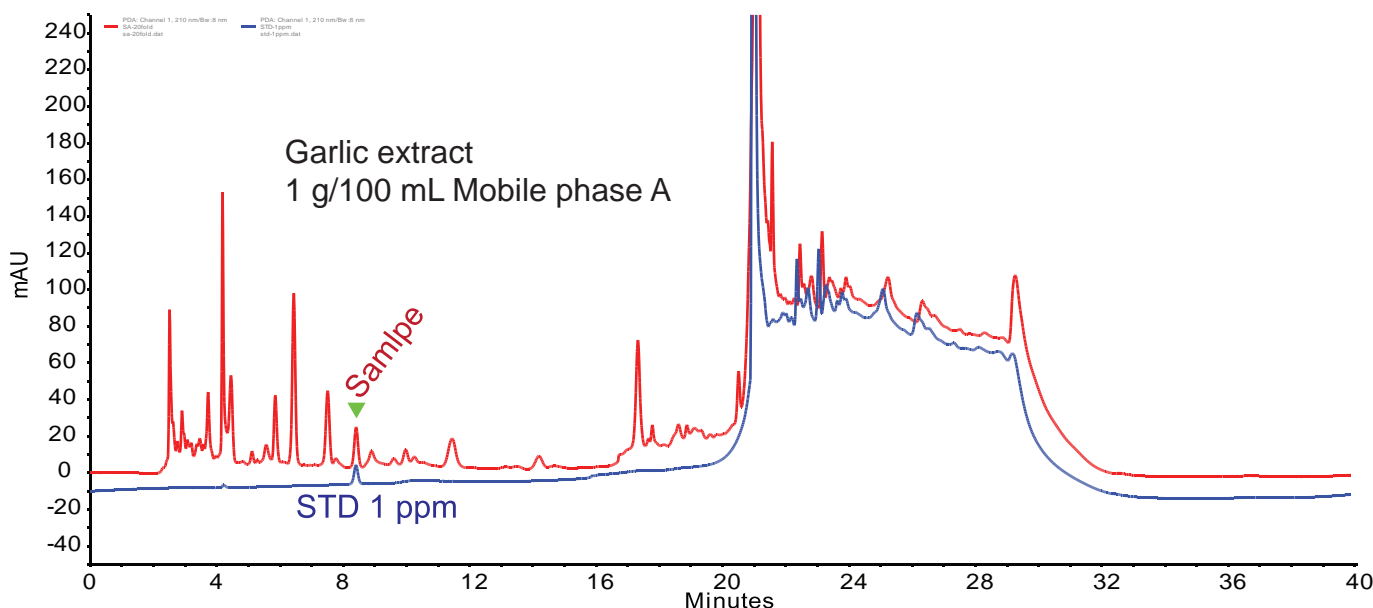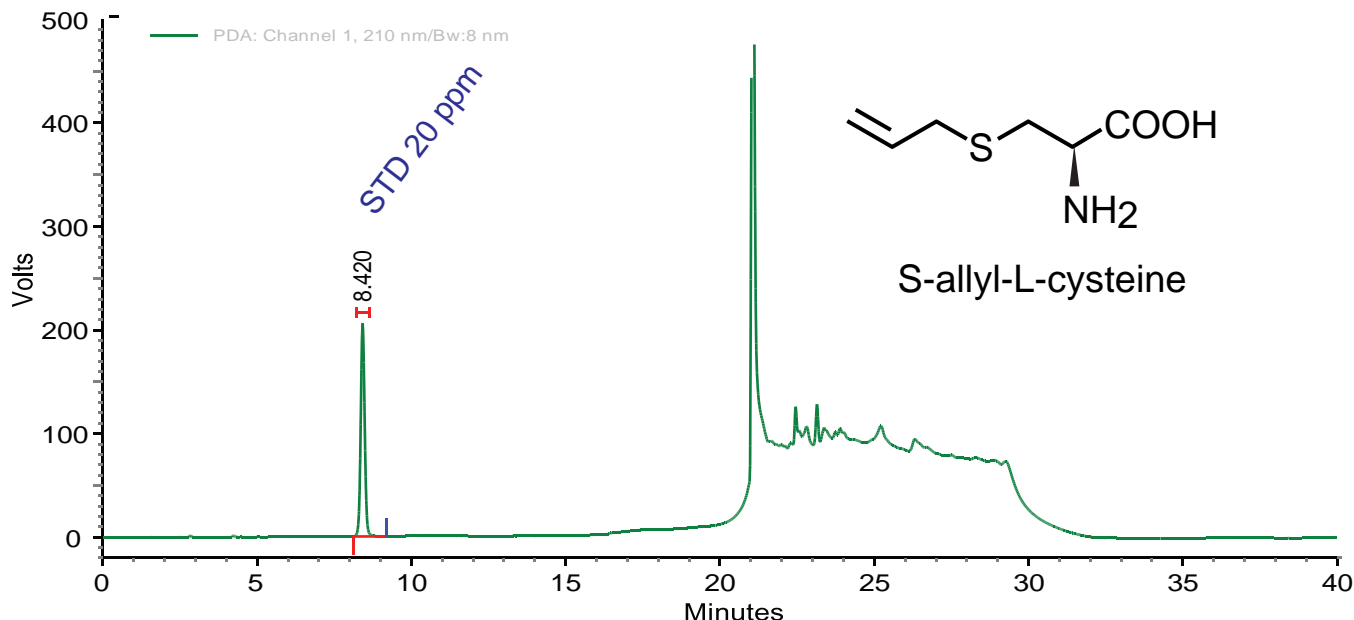

Scherzo SM-C18, 150 x 4.6 mm

A: 0.1% H<sub>3</sub>PO<sub>4</sub> in water

B: 0.1% H<sub>3</sub>PO<sub>4</sub> in 80% acetonitrile

0%B (0-7min), 0-10%B (7-15min), 10-90%B (15-17min), 90%B (17-25min),  
0%B (25-40min)

0.5 mL/min, 40°C

5 μL

210nm

Courtesy of PEAKMAN SP, KOREA