

Unison UK-C8


150 x 4.6 mm

Application

### Sunscreen compounds

紫外線吸収剤

- 1: Diethylamino hydroxybenzoyl hexyl benzoate (INCI) - UVA
- 2: Octocrylene (2-ethylhexyl 2-cyano-3,3-diphenyl-2-propenoate) - UVA
- 3: Ethylhexyl methoxycinnamate (INCI) - UVB
- 4: 2-Ethylhexyl salicylate - UVB
- 5: Bis-Ethylhexyloxyphenol Methoxyphenyl Triazine (INCI) - UVA

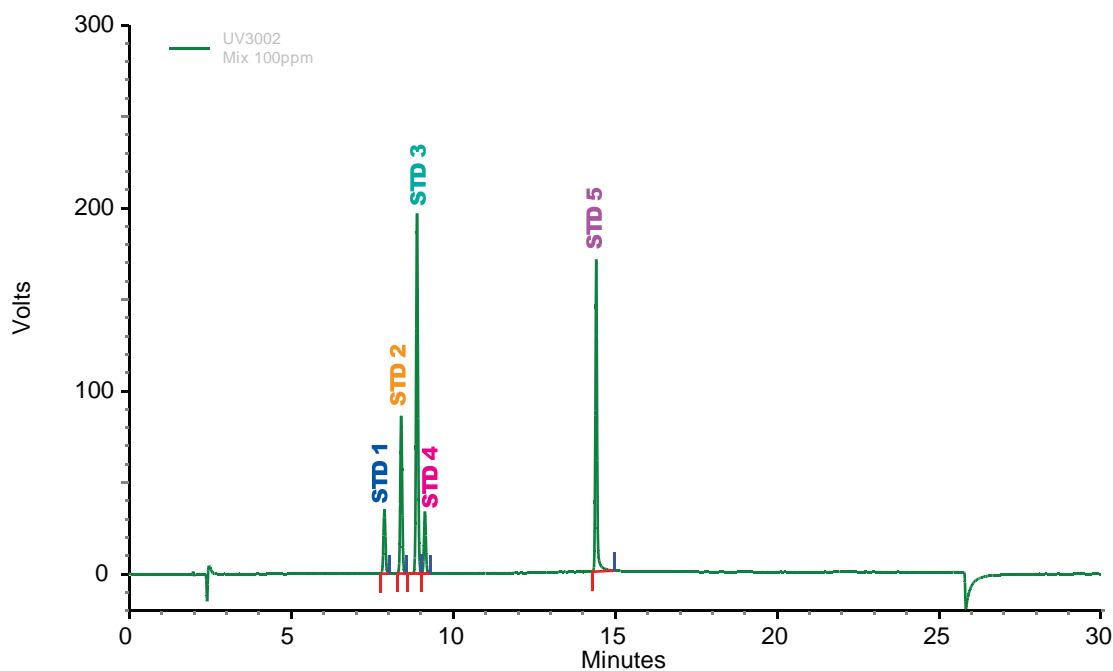

Unison UK-C8, 150 x 4.6 mm

A: 0.1% formic acid

B: methanol

70%B (0-12min), 70-100%B (12-23min), 70%B (23-30min)

1.0 mL/min, 35°C

2 µL (THF solution)

320nm

Courtesy of PEAKMAN SP, KOREA