

Intrada SEC

250 x 4.6 mm

Application

SEC-MS: Elastase, from Porcine pancreas

エラスターゼ(ブタ膵臓由来)のSEC-MS分析





SEC-MS: Elastase

Intrada SEC

Elastase, from Porcine pancreas

20170807_PLUS_IntradaSEC_100mMAF_2017... 8/7/2017 9:25:00 AM

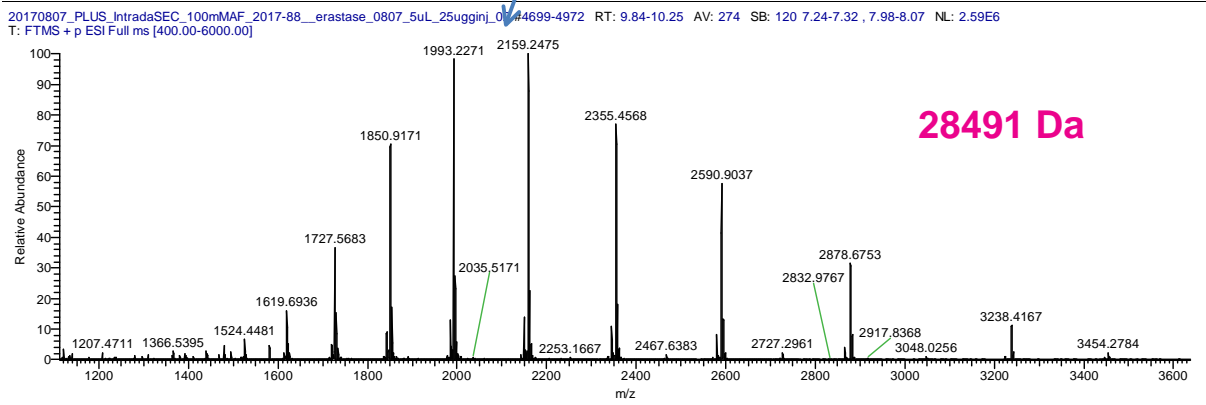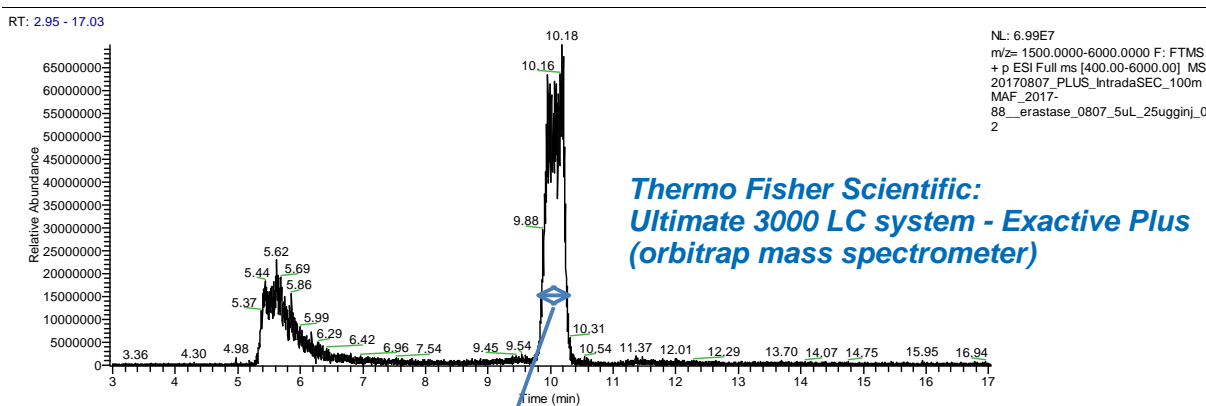

Intrada SEC, 250x4.6mm
100mM ammonium formate
0.3 mL/min
37°C, 5µL (25µg)

Courtesy of Dr. Keiko KUWATA, Nagoya University, JAPAN