

Unison UK-C18 MF

150 x 2 mm


Application

Metal Free Column: Nucleoside triphosphate (NTP), Trinucleotide


メタルフリーカラム: ヌクレオシド三リン酸, トリヌクレオチド

Metal Free Column

Unison UK-C18 MF

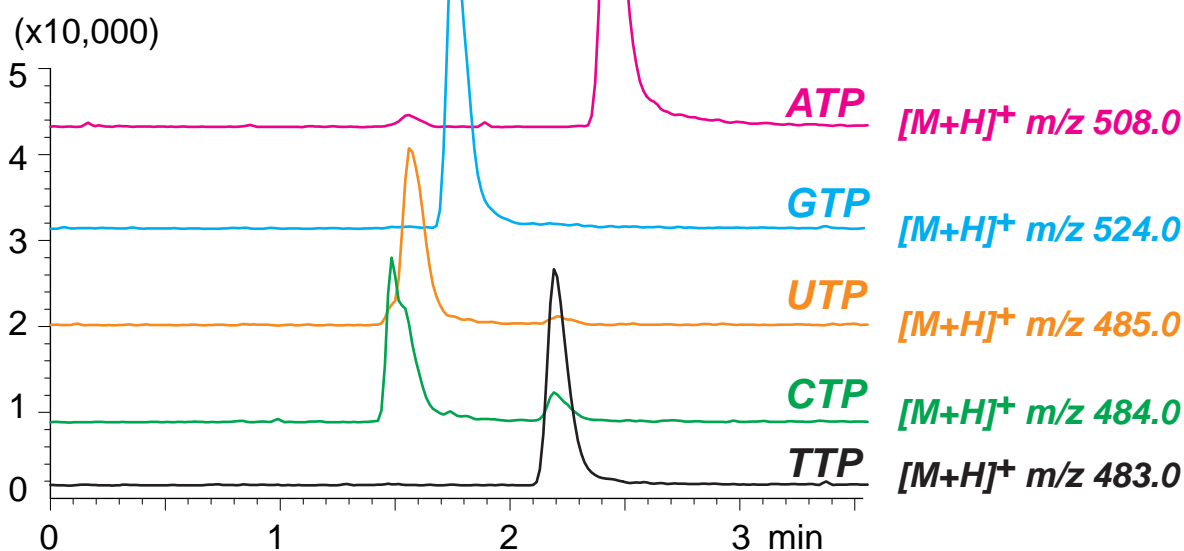

Metal Column

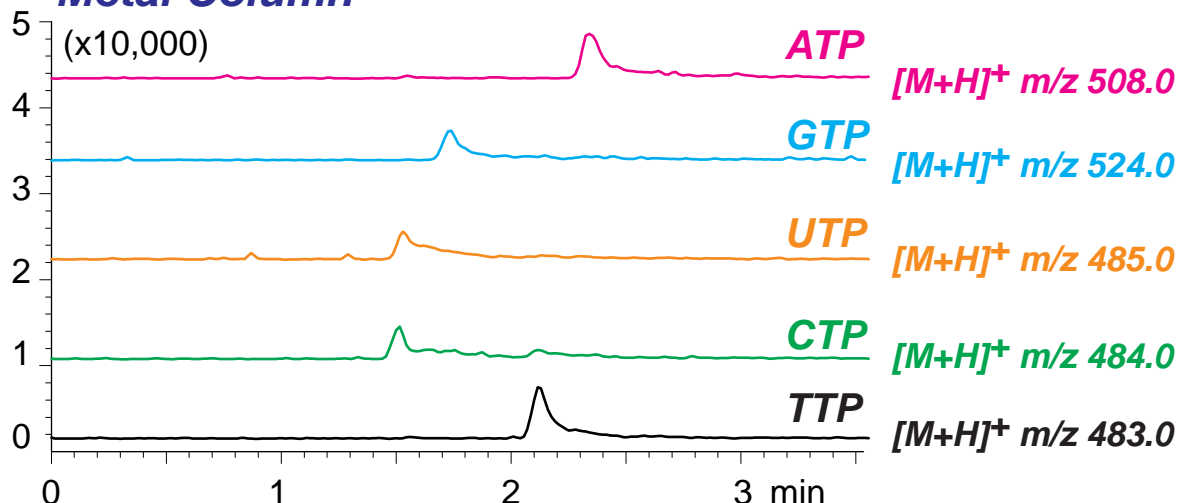

Unison UK-C18 MF, 150 x 2 mm
 5mM ammonium formate
 0.2 mL/min (8MPa), 37°C
 2 μ L (10 μ g/mL in 5mM ammonium formate)
 Single Quad. MS (ESI)