



**Scherzo SM-C18 MF**

100 x 2 mm

Application

**Metal-Free Column: phosphoserine**

メタルフリーカラム: ホスホセリン

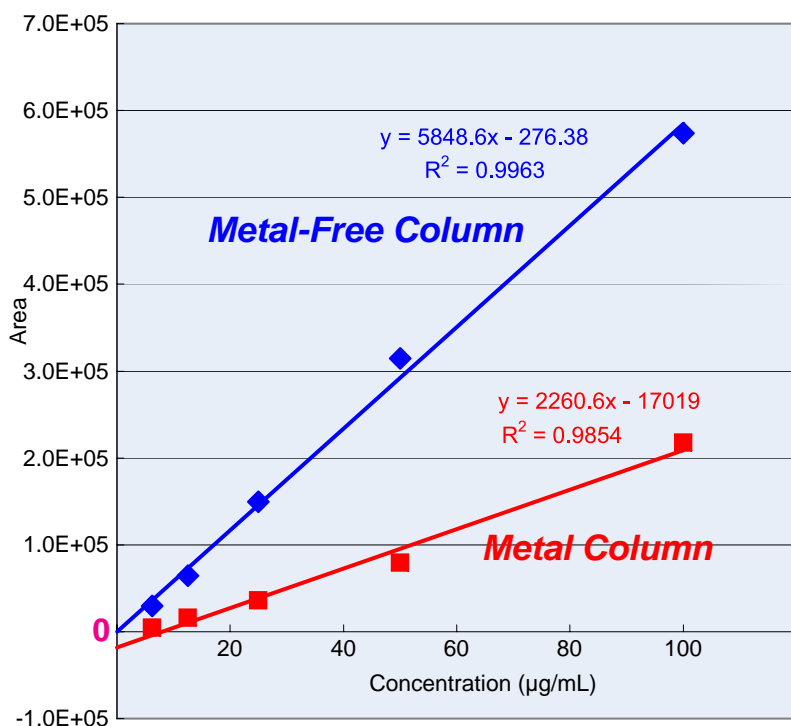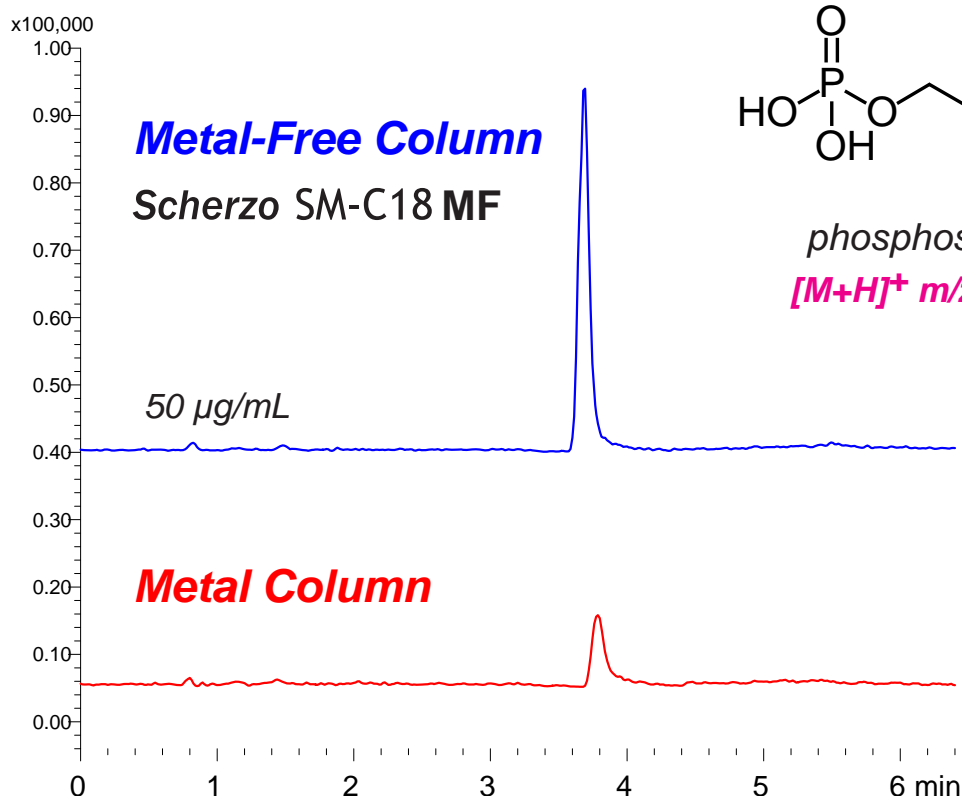

Scherzo SM-C18 MF, 100 x 2 mm

A: water / HCOOH = 100 / 0.1  
 B: ACN / 100mM HCOONH4 = 20 / 80

0%B (0-2min)  
 0-100%B (2-4 min)  
 100%B (4-7min)

0.2 mL/min (6MPa), 37°C  
 5 µL (1% HCOOH)

Single Quad. MS (ESI, positive)