



Cadenza HS-C18

250 x 4.6 mm

Application

Nitrovin (animal medicine)

ニトロビン (動物用医薬品)

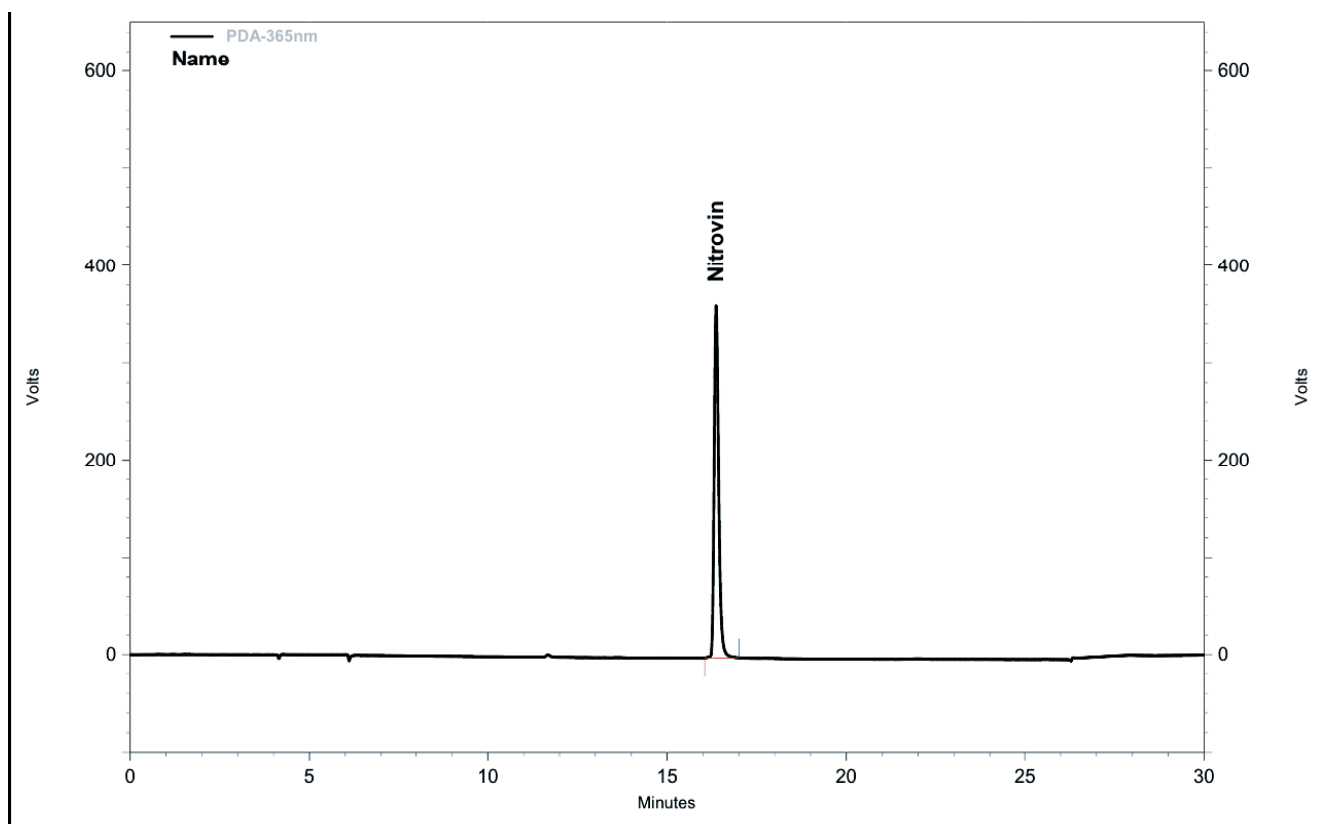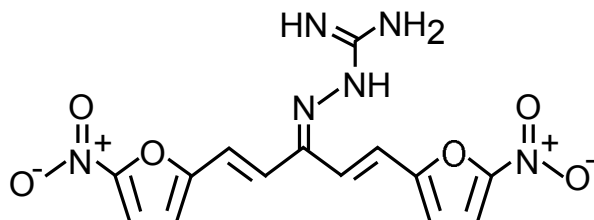

Cadenza HS-C18, 250 x 4.6 mm
 A: 0.1% formic acid in water
 B: 0.1% formic acid in acetonitrile
 0-80 % B(0~23 min), 80-0% B(23~23.1 min), 0% B(23.1~30 min)
 0.8 mL/min, 35°C
 5 µL
 365nm

Courtesy of PEAKMAN SP, KOREA