


Unison UK-C1  
 Unison UK-C18  
 Cadenza CD-C18

50 x 3 mm

Application

**Tizanidine (muscle-relaxant) impurity**

チザニジン(筋緊張緩和剤)不純物

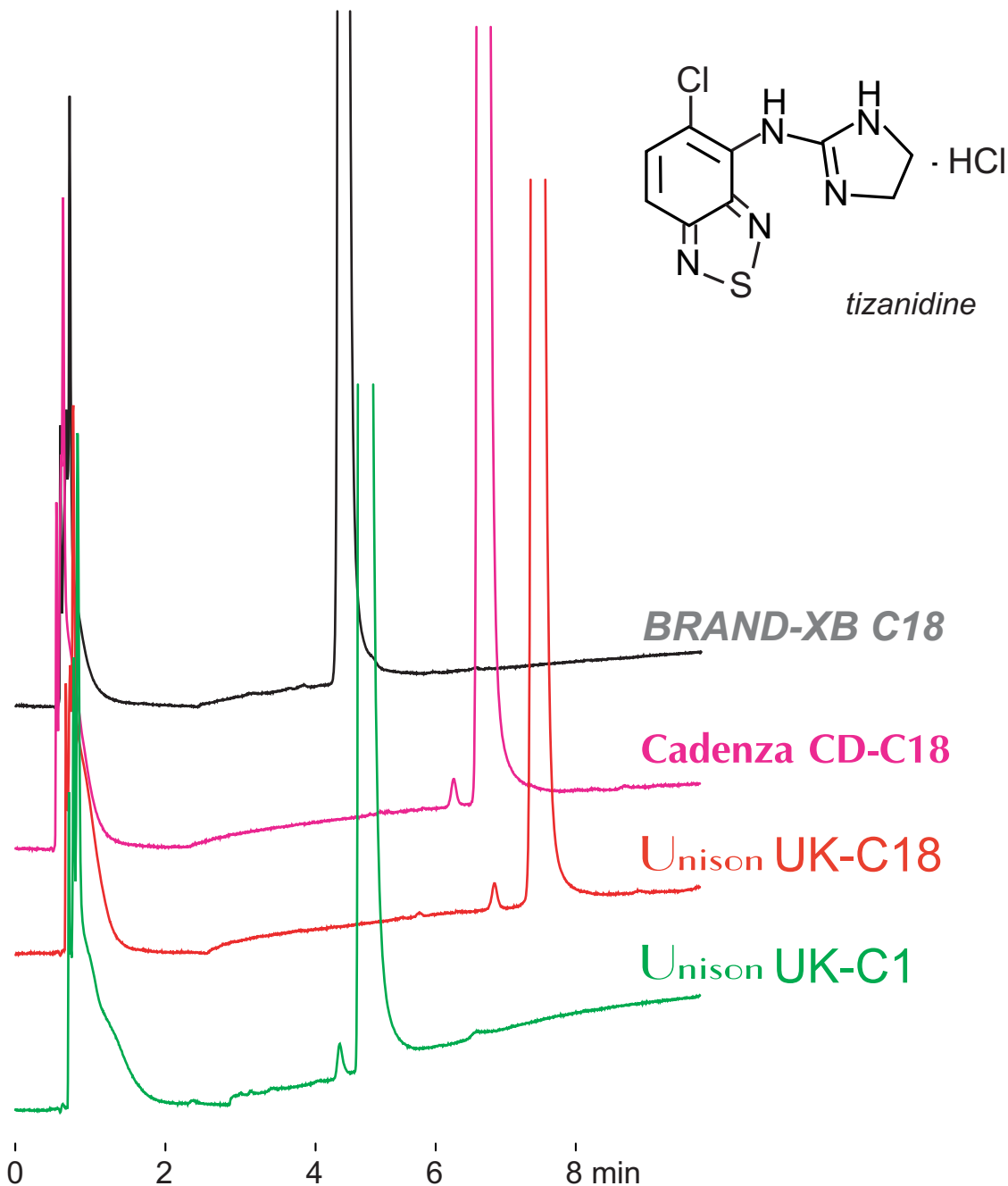

50 x 3 mm  
 A: 10mM HCOONH<sub>4</sub>  
 B: 50mM HCOONH<sub>4</sub> / acetonitrile = 80 / 20  
 0 - 10%B (0 - 10min)  
 0.4mL/min (5MPa), 37°C  
 2 μL (0.1mg/mL, DMSO), 320nm