


DACAPO DX-C18

250 x 3 mm

Application

LC-MS: Omega-3 and Omega-6 fatty acids

オメガ-3およびオメガ-6脂肪酸のLC-MS分析

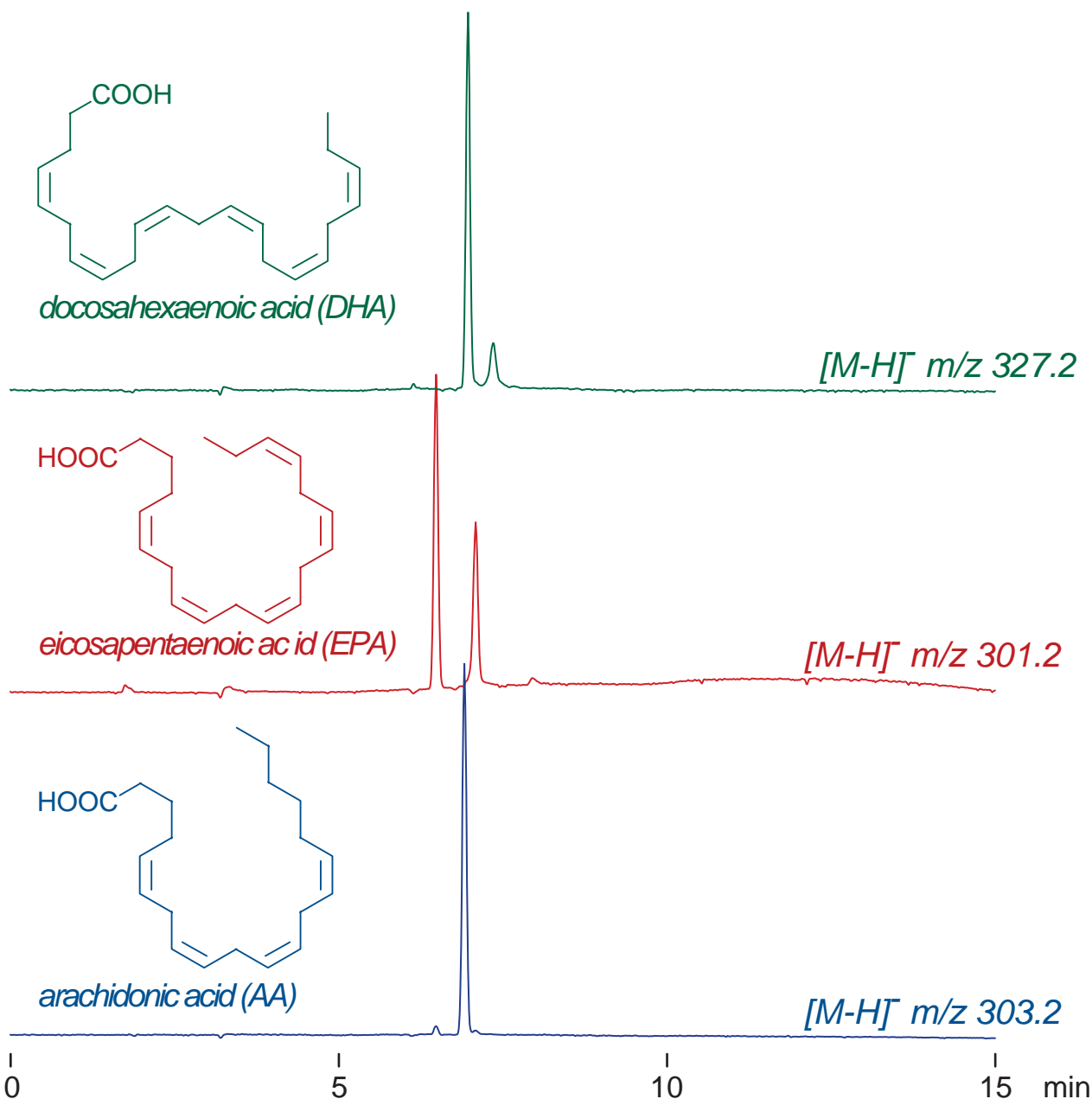

Dacapo DX-C18, 250 x 3 mm

A: 10mM ammonium hydroxide

B: acetonitrile

10-100%B (0-15min)

0.4mL/min (21MPa), 37°C, MS (ESI, SIM, Negative)

10µL (1µg/mL, ACN/50mM NH₄OH = 50/50)