


Cadenza CX-C18  
 Unison UK-C18  
 Cadenza CD-C18

75 x 4.6 mm


Application

**Levofloxacin (new quinolone, antibiotics)**

レボフロキサシン (ニューキノロン系抗生物質)

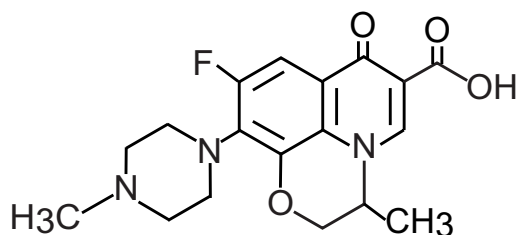

levofloxacin

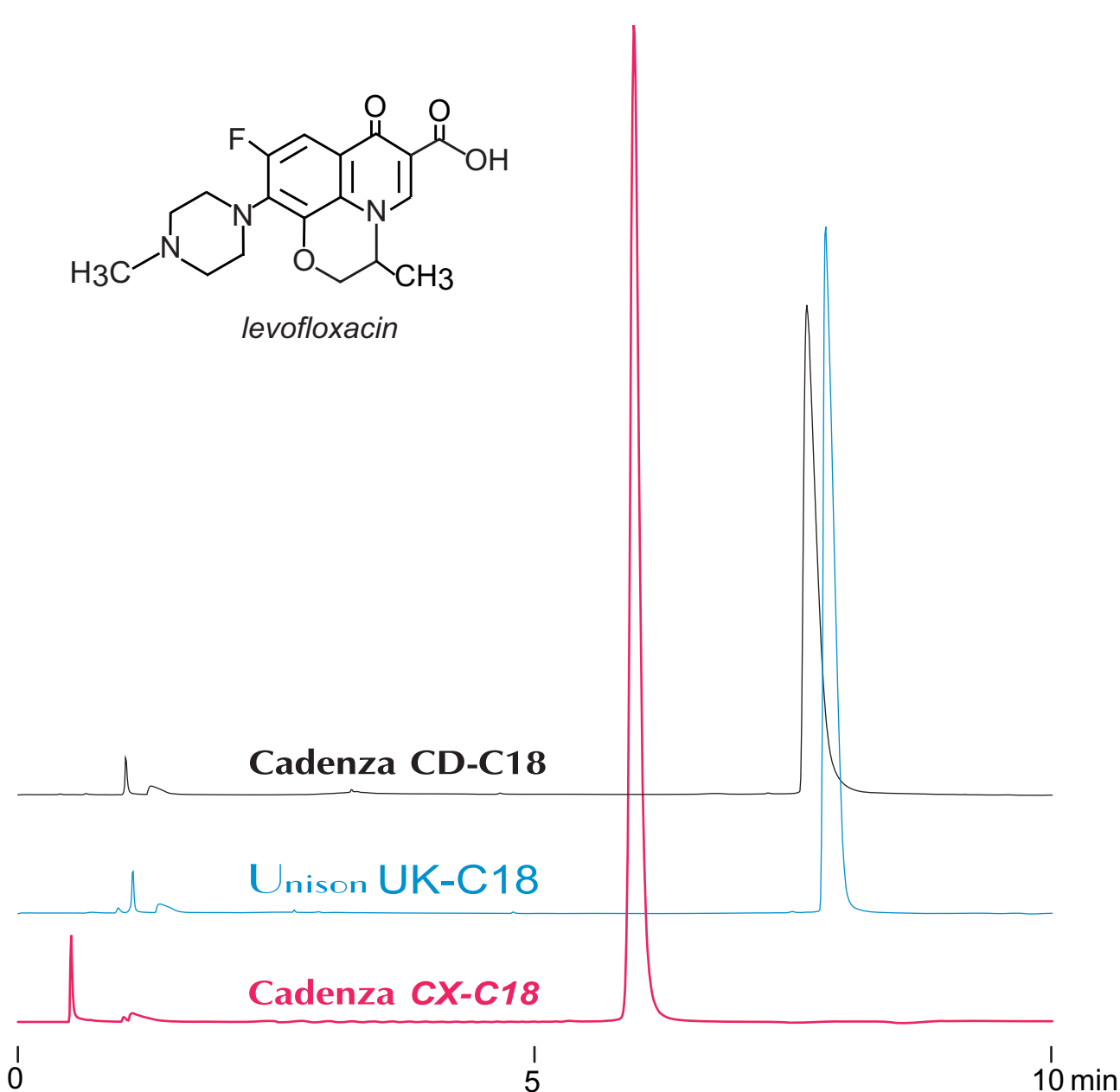

75 x 4.6 mm

A: water / formic acid= 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-50%B (0-10min)

1.0 mL/min (7MPa), 37°C

2 µL (0.1mg/mL, 1%HCOOH), 260nm