



Cadenza CX-C18
 Unison UK-C18
 Cadenza CD-C18

75 x 4.6 mm

Application

Propiverine (Anticholinergic, Ca-antagonist, Antimuscarinic)

プロピペリン (抗コリン薬, Ca拮抗薬, 抗ムスカリン作用薬)

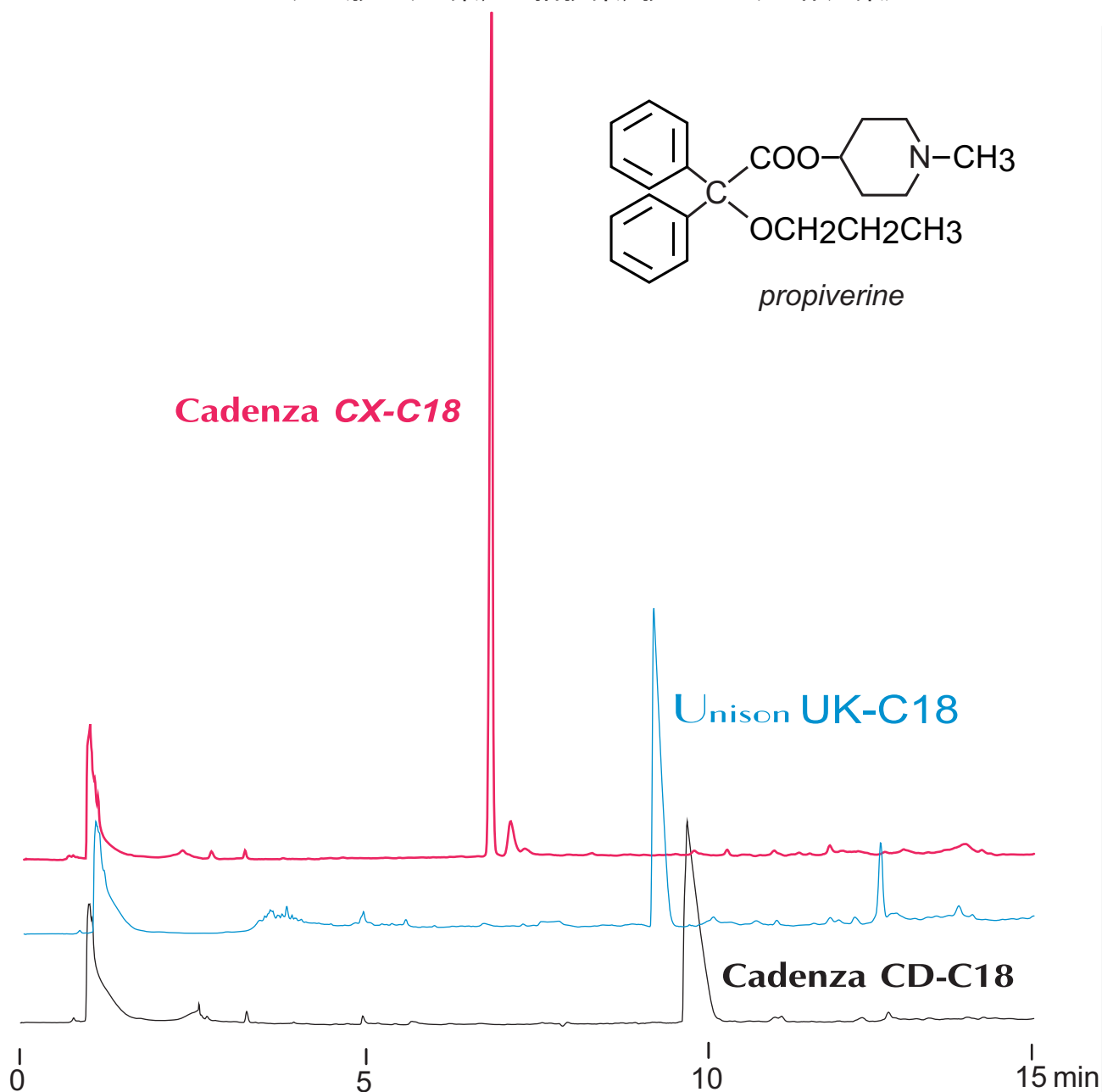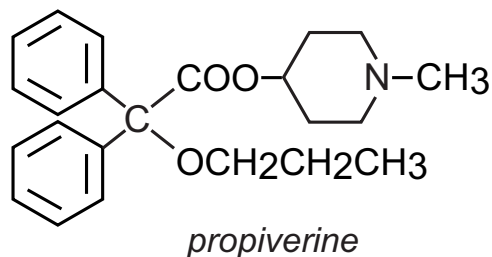

75 x 4.6 mm

A: water / formic acid= 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-75%B (0-15min)

1.0 mL/min (7MPa), 37°C

10 µL (0.5 mg/mL, 1% HCOOH), 260nm