



Unison UK-C18
Cadenza CD-C18
Cadenza CX-C18

75 x 4.6 mm

Application

Etodolac (nonsteroidal anti-inflammatory drug, NSAIDs)

エトドラク (非ステロイド性鎮痛薬・抗炎症薬)

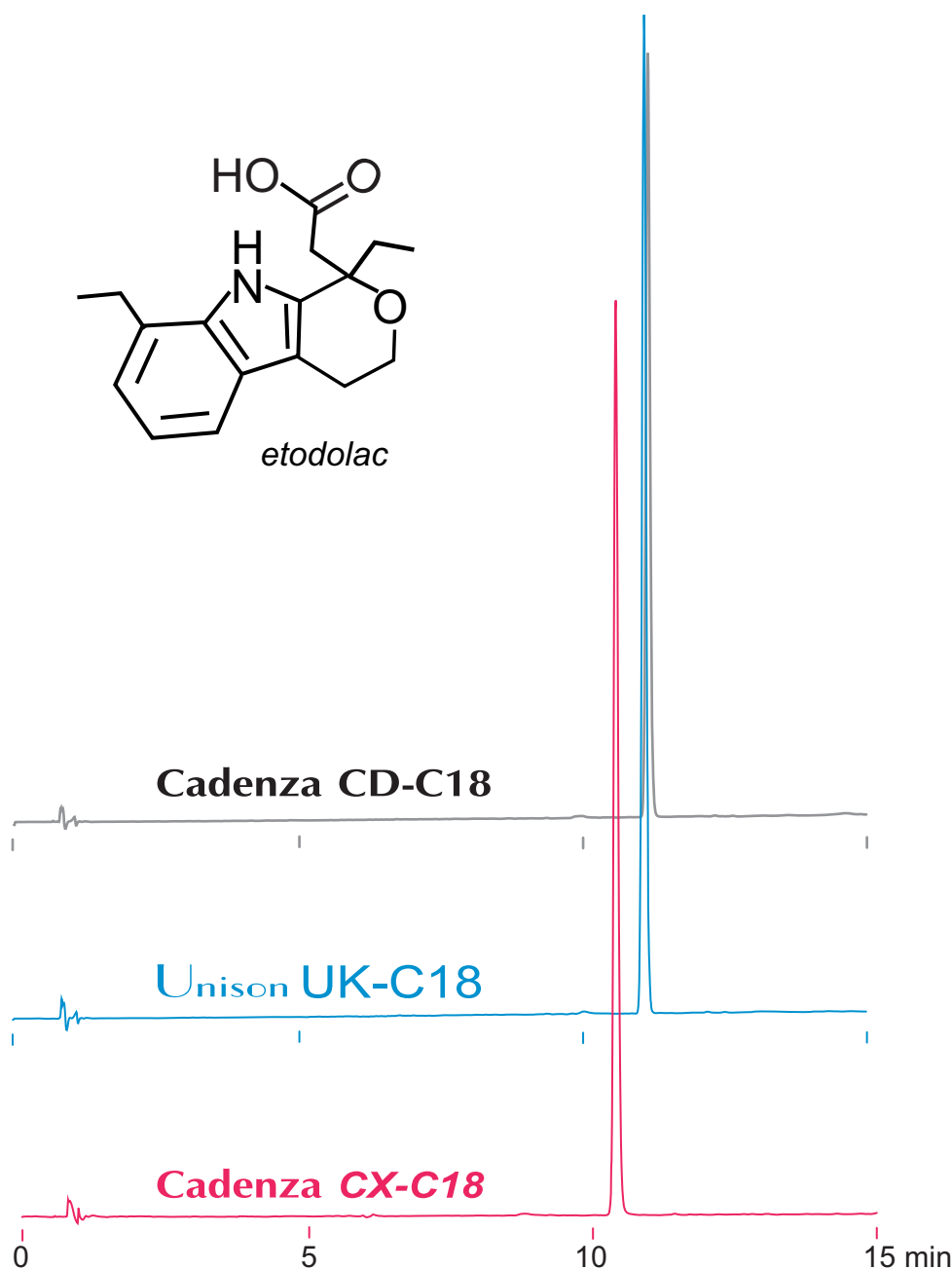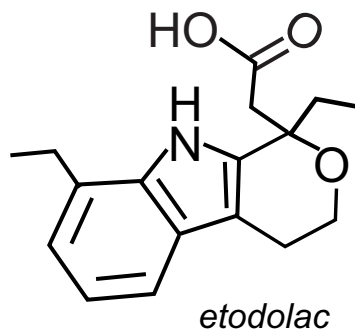

75 x 4.6 mm

A: water / formic acid= 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-75%B (0-15min)

1.0 mL/min (7MPa), 37°C, 260nm

2 µL (0.5 mg/mL, 0.1% HCOOH in MeOH)