


Unison UK-C18  
 Unison UK-C8  
 Unison UK-C1

75 x 4.6 mm

Application

Tipepidine hibenzate (Antitussive, Expectorant)

チペピジン ヒベンズ酸塩 (鎮咳去痰薬)

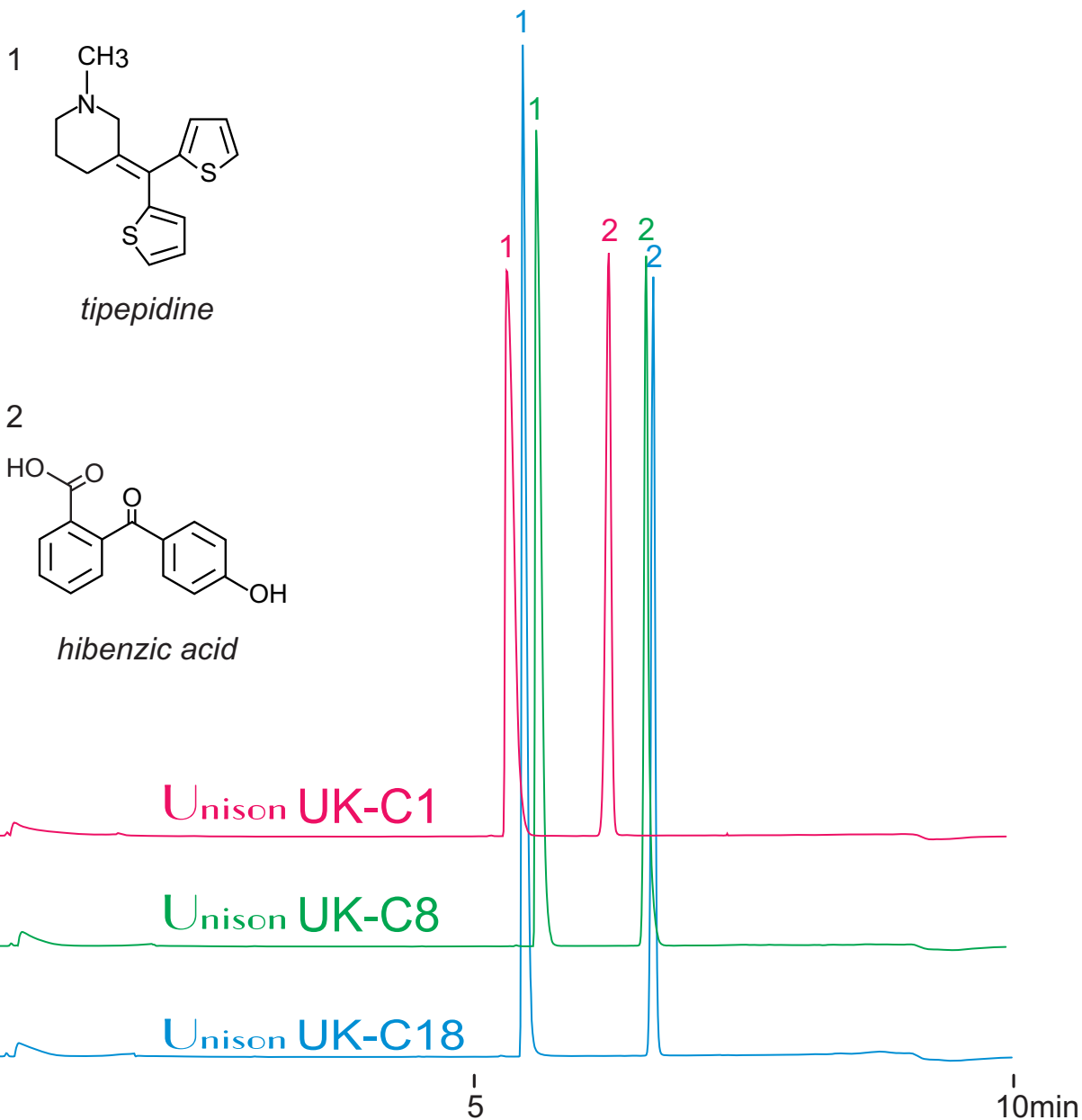

75 x 4.6 mm

A: water / formic acid= 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-50%B (0-10min)

1.0 mL/min (7MPa), 37°C, 260nm

2 μL (1 mg/mL, DMSO)