



Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18

50 x 3 mm

Application

Antibiotic (vancomycin)

抗生物質 (バンコマイシン)

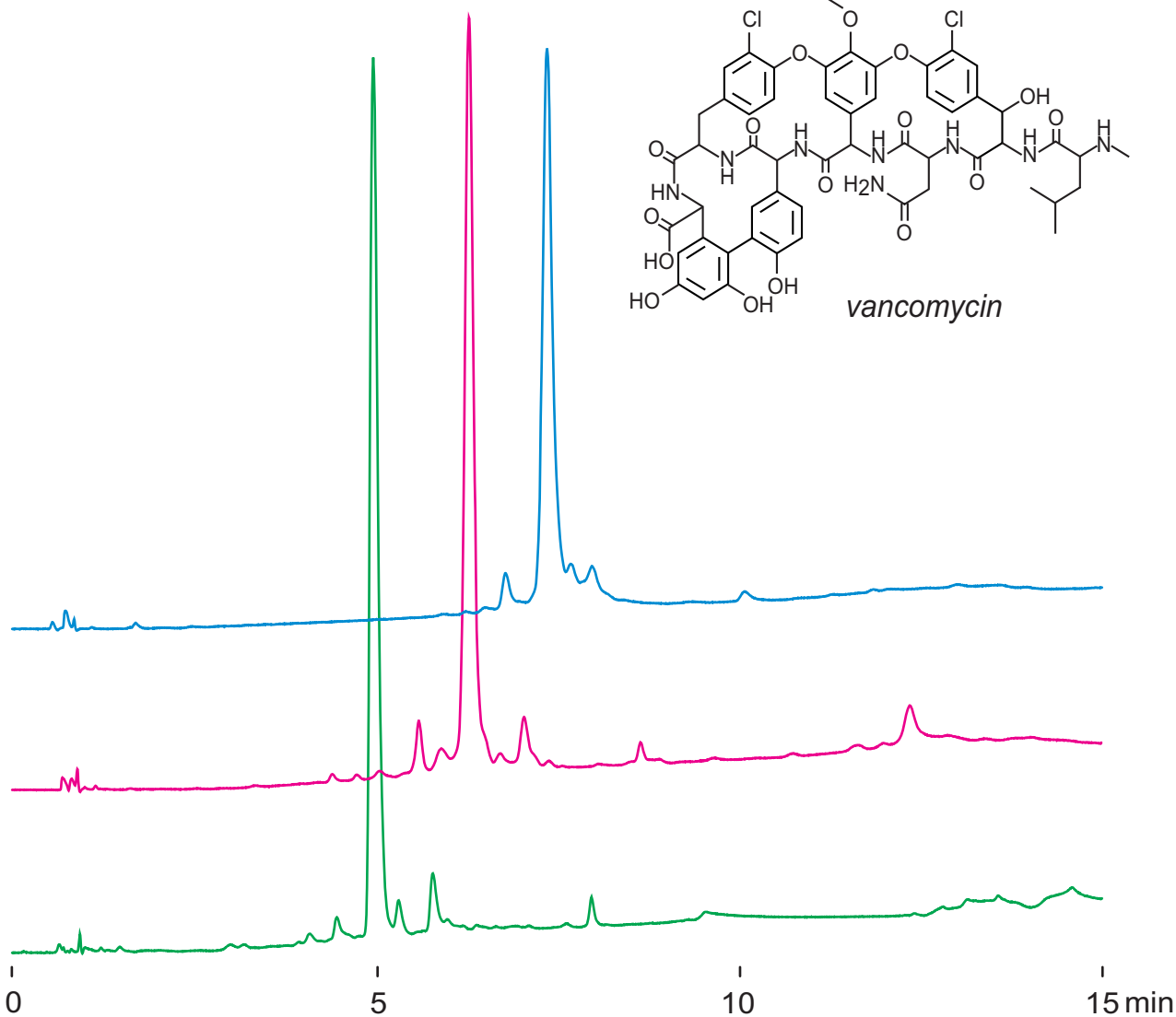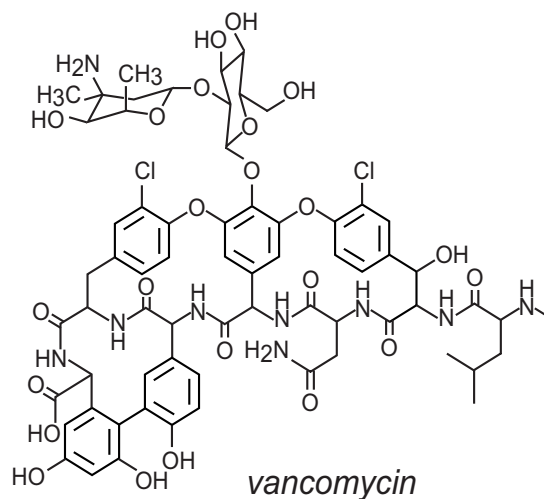

50 x 3 mm

A: 10mM ammonium acetate / acetonitrile = 90 / 10

B: 100mM ammonium acetate / acetonitrile = 50 / 50

0-100%B (0-15 min)

0.4mL/min (5MPa) , 37°C, 260nm

2μL (1mg/mL, 0.1% HCOOH)